

Sample name: **M01**  
 Assay name: **UV-metric pKa**  
 Assay ID: **171-16001**  
 Filename: **C:\Sirius\_T3\Mehtap\20170915\_exp03\_uv\_M01-M14\171-16001\_M01\_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 12:18:58 AM**  
 Analyst: **Dorothy Leverse**  
 Instrument ID: **T311053**

## Results

pKa 1 **9.53**  
 RMSD **0.002 0.002**  
 Chi squared **0.0063**  
 PCA calculated number of pKas **1**  
 Average ionic strength **0.158 M**  
 Average temperature **24.9°C**  
 Analyte concentration range **115.2 µM to 104.2 µM**

Number of pKas source **Predicted**  
 Wavelength clipping **230.0 nm to 450.0 nm**  
 pH clipping **1.267 to 12.716**


## Warnings and errors

Errors None  
 Warnings None

## Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
<b>Assay Medium</b>				
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			









## Graphs

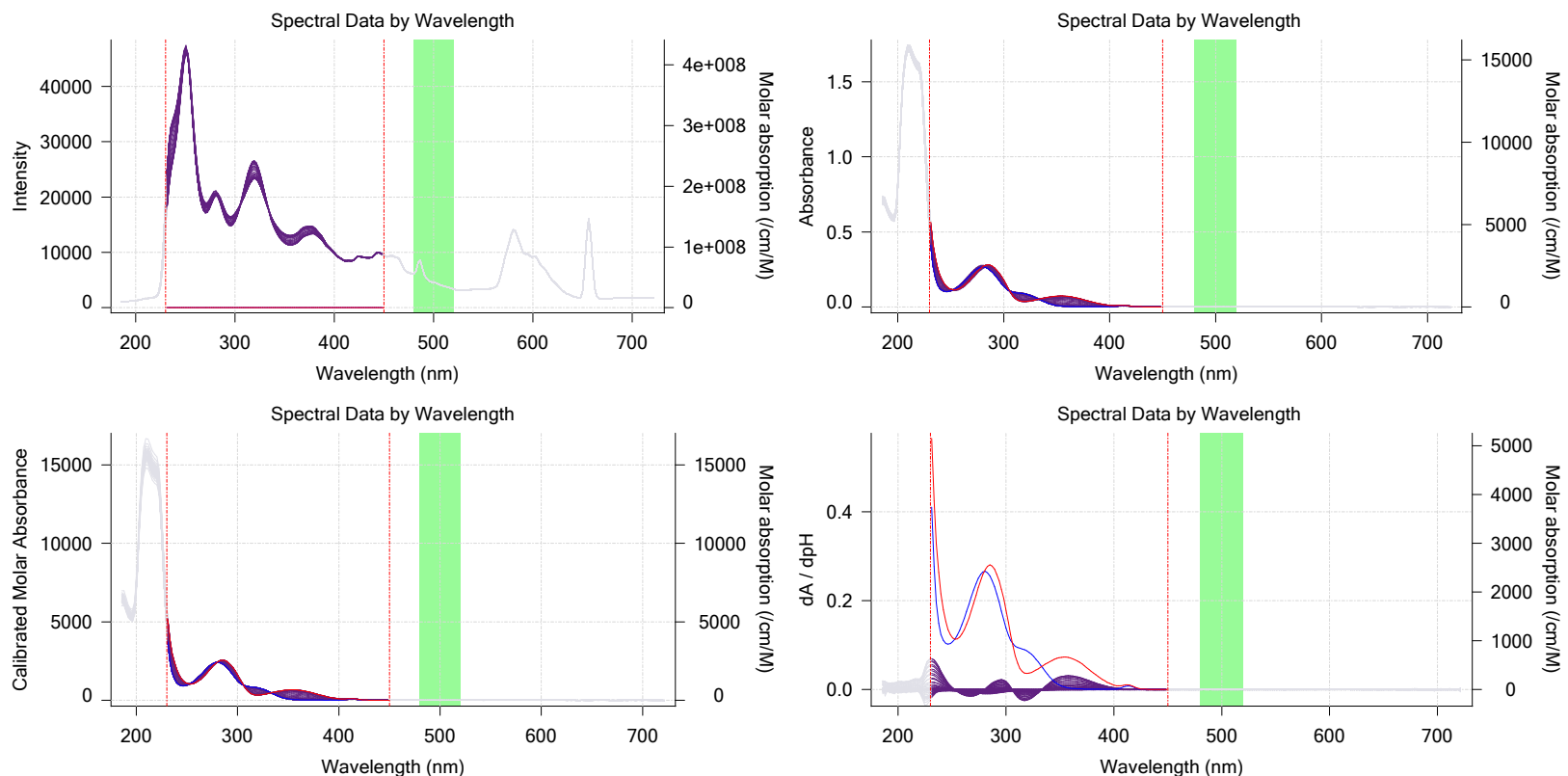

Sample name: **M01**  
 Assay name: **UV-metric pKa**  
 Assay ID: **171-16001**  
 Filename: **C:\Sirius\_T3\Mehtap\20170915\_exp03\_uv\_M01-M14\171-16001\_M01\_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 12:18:58 AM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**


## Graphs (continued)

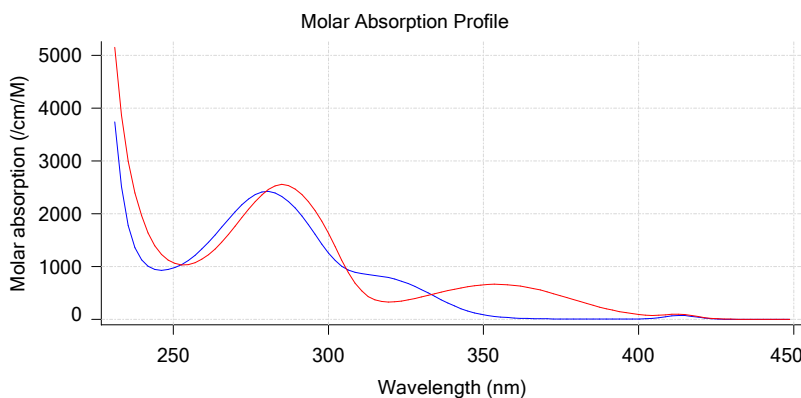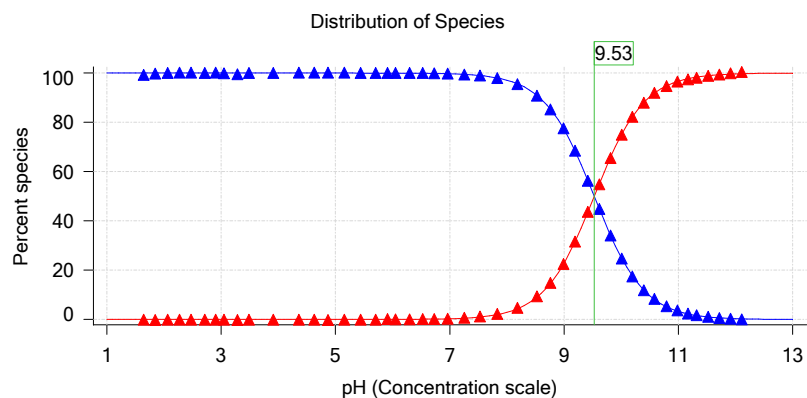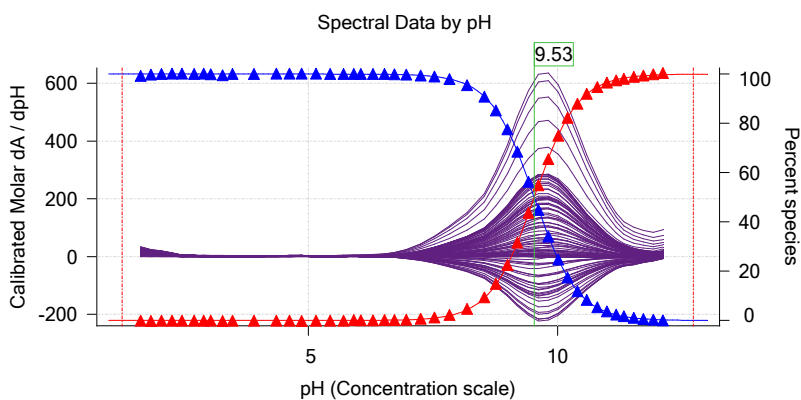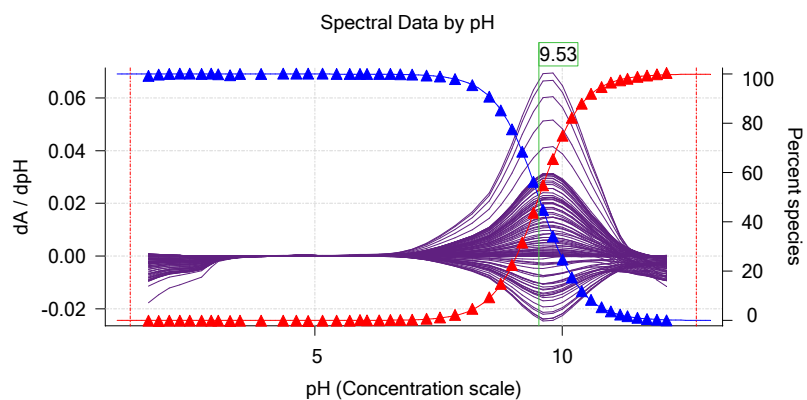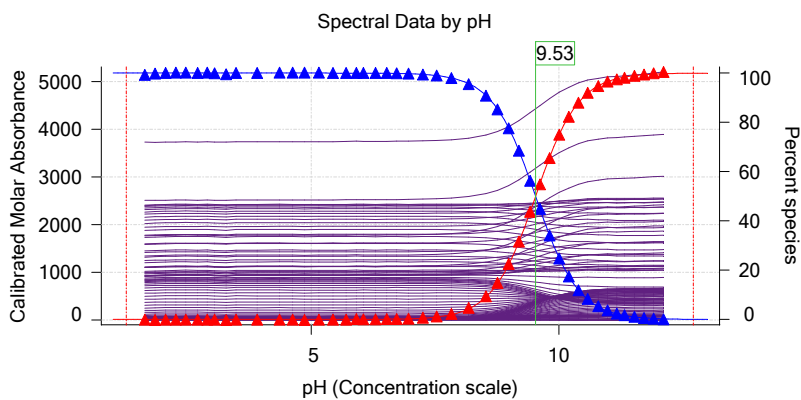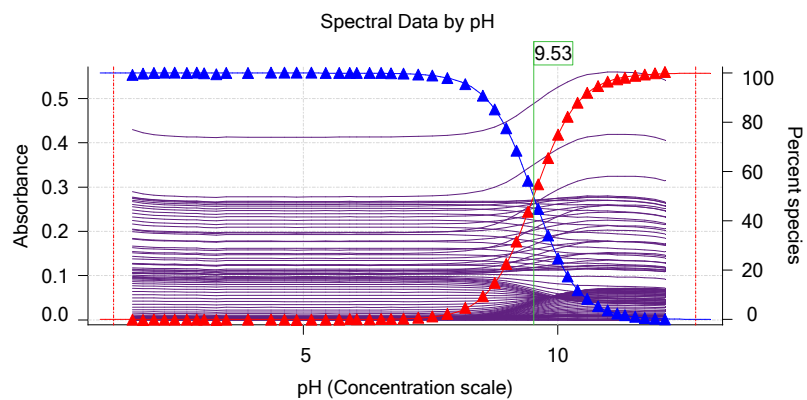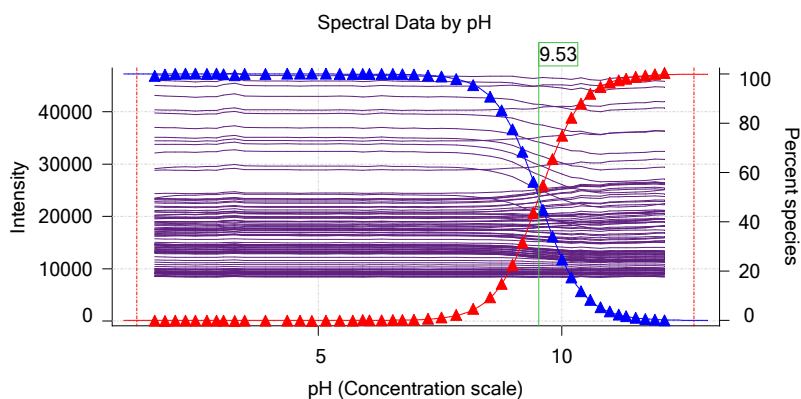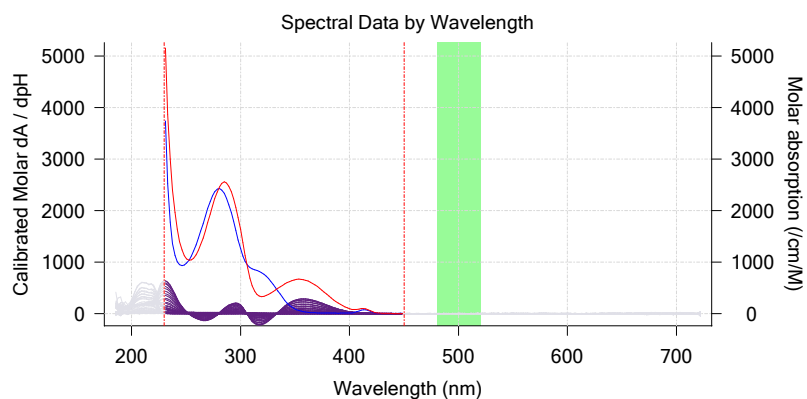

Sample name: **M01**  
 Assay name: **UV-metric pKa**  
 Assay ID: **171-16001**  
 Filename: **C:\Sirius\_T3\Mehtap\20170915\_exp03\_uv\_M01-M14\171-16001\_M01\_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 12:18:58 AM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Graphs (continued)

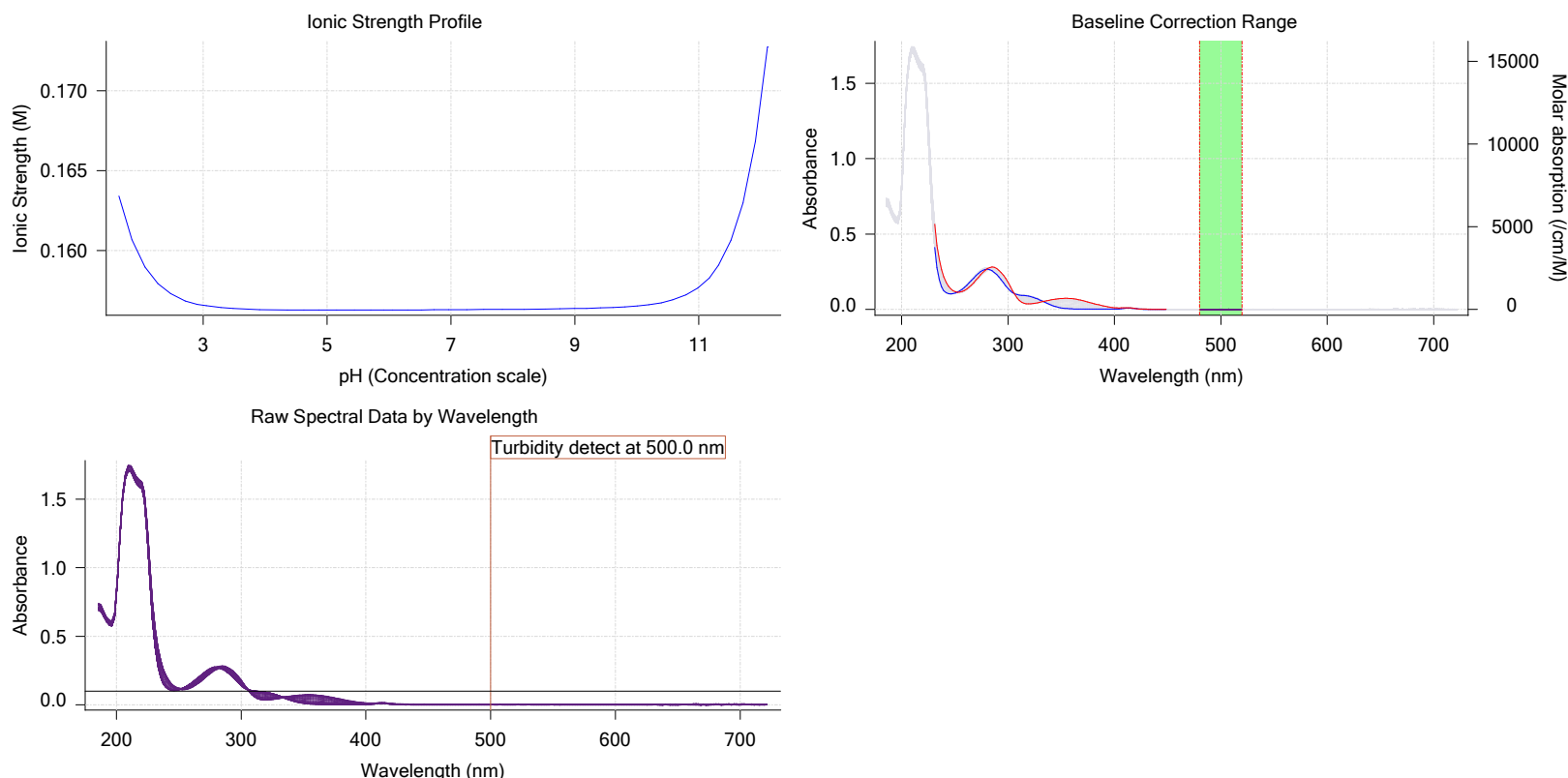

## Events

Time	Event	Water	Acid	Base	Buffer	pH	dpH/dt	pH R-squared	pH SD
3:14.8	Dark spectrum								
3:16.3	Reference spectrum								
3:43.9	Volume reset due to vial change								
5:14.4	Initial pH = 7.58								
6:27.5	Data point 4	1.50000 mL	0.07173 mL	0.00000 mL	0.02500 mL	1.767	-0.00977	0.78966	0.0005
6:56.4	Data point 5	1.50000 mL	0.07173 mL	0.02585 mL	0.02500 mL	1.968	0.00319	0.12145	0.0004
7:13.3	Data point 6	1.50000 mL	0.07173 mL	0.04292 mL	0.02500 mL	2.178	0.00256	0.02371	0.0008
7:30.1	Data point 7	1.50000 mL	0.07173 mL	0.05343 mL	0.02500 mL	2.387	0.00134	0.04947	0.0003
7:47.0	Data point 8	1.50000 mL	0.07173 mL	0.05988 mL	0.02500 mL	2.587	0.00808	0.74621	0.0004
8:03.8	Data point 9	1.50000 mL	0.07173 mL	0.06395 mL	0.02500 mL	2.824	0.00421	0.55048	0.0002
8:35.8	Data point 10	1.50000 mL	0.07173 mL	0.06649 mL	0.02500 mL	3.018	0.00305	0.27409	0.0002
8:52.3	Data point 11	1.50000 mL	0.07173 mL	0.06799 mL	0.02500 mL	3.163	0.00985	0.83321	0.0005
9:24.4	Data point 12	1.50000 mL	0.07173 mL	0.06928 mL	0.02500 mL	3.398	-0.00071	0.02190	0.0002
9:56.2	Data point 13	1.50000 mL	0.07173 mL	0.06992 mL	0.02500 mL	3.599	0.01372	0.92930	0.0007
10:12.7	Data point 14	1.50000 mL	0.07173 mL	0.07032 mL	0.02500 mL	4.025	0.03398	0.95611	0.0017
10:34.4	Data point 15	1.50000 mL	0.07173 mL	0.07060 mL	0.02500 mL	4.475	0.05093	0.94773	0.0025
11:00.9	Data point 16	1.50000 mL	0.07173 mL	0.07074 mL	0.02500 mL	4.740	0.09335	0.93864	0.0047
11:23.0	Data point 17	1.50000 mL	0.07173 mL	0.07081 mL	0.02500 mL	4.981	0.09598	0.97415	0.0048
11:47.2	Data point 18	1.50000 mL	0.07173 mL	0.07086 mL	0.02500 mL	5.268	0.09687	0.96095	0.0049
12:11.9	Data point 19	1.50000 mL	0.07173 mL	0.07091 mL	0.02500 mL	5.554	0.06861	0.89376	0.0035
12:33.6	Data point 20	1.50000 mL	0.07173 mL	0.07095 mL	0.02500 mL	5.820	0.04899	0.77276	0.0027
12:55.3	Data point 21	1.50000 mL	0.07173 mL	0.07100 mL	0.02500 mL	6.027	0.04894	0.51513	0.0033
13:11.8	Data point 22	1.50000 mL	0.07173 mL	0.07105 mL	0.02500 mL	6.167	0.07505	0.73890	0.0043
13:38.7	Data point 23	1.50000 mL	0.07173 mL	0.07112 mL	0.02500 mL	6.414	0.03445	0.45553	0.0025
14:05.3	Data point 24	1.50000 mL	0.07173 mL	0.07121 mL	0.02500 mL	6.643	0.02580	0.56317	0.0017
14:31.9	Data point 25	1.50000 mL	0.07173 mL	0.07131 mL	0.02500 mL	6.842	0.03556	0.62571	0.0022

Sample name: **M01**  
 Assay name: **UV-metric pKa**  
 Assay ID: **17I-16001**  
 Filename: **C:\Sirius\_T3\Mehtap\20170915\_exp03\_uv\_M01-M14\17I-16001\_M01\_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 12:18:58 AM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Events (continued)

Time	Event	Water	Acid	Base	Buffer	pH	dpH/dt	pH R-squared	pH SD	dpH/dt time
15:03.8	Data point 26	1.50000 mL	0.07173 mL	0.07140 mL	0.02500 mL	7.082	0.08064	0.86740	0.00427	10.0 s
15:35.7	Data point 27	1.50000 mL	0.07173 mL	0.07150 mL	0.02500 mL	7.372	0.08220	0.81083	0.00451	11.0 s
16:03.6	Data point 28	1.50000 mL	0.07173 mL	0.07157 mL	0.02500 mL	7.648	0.09094	0.85201	0.00486	12.5 s
16:32.9	Data point 29	1.50000 mL	0.07173 mL	0.07164 mL	0.02500 mL	7.948	0.08823	0.79089	0.00489	13.0 s
17:02.5	Data point 30	1.50000 mL	0.07173 mL	0.07171 mL	0.02500 mL	8.299	0.07729	0.74716	0.00441	13.0 s
17:32.3	Data point 31	1.50000 mL	0.07173 mL	0.07180 mL	0.02500 mL	8.643	0.08364	0.71607	0.00488	10.5 s
17:54.4	Data point 32	1.50000 mL	0.07173 mL	0.07190 mL	0.02500 mL	8.879	0.07257	0.74836	0.00418	10.5 s
18:21.9	Data point 33	1.50000 mL	0.07173 mL	0.07199 mL	0.02500 mL	9.106	0.04595	0.80980	0.00252	10.0 s
18:48.6	Data point 34	1.50000 mL	0.07173 mL	0.07211 mL	0.02500 mL	9.308	0.02912	0.71643	0.00172	10.0 s
19:15.3	Data point 35	1.50000 mL	0.07173 mL	0.07227 mL	0.02500 mL	9.534	0.01072	0.49674	0.00075	10.0 s
19:42.2	Data point 36	1.50000 mL	0.07173 mL	0.07248 mL	0.02500 mL	9.734	-0.00153	0.04936	0.00034	10.0 s
20:09.1	Data point 37	1.50000 mL	0.07173 mL	0.07274 mL	0.02500 mL	9.929	-0.00778	0.61653	0.00049	10.0 s
20:40.9	Data point 38	1.50000 mL	0.07173 mL	0.07312 mL	0.02500 mL	10.126	-0.01119	0.79580	0.00062	10.0 s
21:12.6	Data point 39	1.50000 mL	0.07173 mL	0.07361 mL	0.02500 mL	10.316	-0.01593	0.93993	0.00082	10.0 s
21:44.4	Data point 40	1.50000 mL	0.07173 mL	0.07437 mL	0.02500 mL	10.512	-0.01977	0.94622	0.00100	10.0 s
22:16.5	Data point 41	1.50000 mL	0.07173 mL	0.07549 mL	0.02500 mL	10.701	-0.01701	0.94144	0.00087	10.0 s
22:43.7	Data point 42	1.50000 mL	0.07173 mL	0.07747 mL	0.02500 mL	10.912	-0.01150	0.87274	0.00061	10.0 s
22:10.6	Data point 43	1.50000 mL	0.07173 mL	0.07963 mL	0.02500 mL	11.102	-0.01194	0.90435	0.00062	10.0 s
23:27.3	Data point 44	1.50000 mL	0.07173 mL	0.08286 mL	0.02500 mL	11.288	-0.01112	0.90441	0.00058	10.0 s
23:44.0	Data point 45	1.50000 mL	0.07173 mL	0.08779 mL	0.02500 mL	11.440	-0.01149	0.89534	0.00060	10.0 s
24:16.2	Data point 46	1.50000 mL	0.07173 mL	0.09675 mL	0.02500 mL	11.637	-0.01180	0.89407	0.00062	10.0 s
24:43.6	Data point 47	1.50000 mL	0.07173 mL	0.10985 mL	0.02500 mL	11.830	-0.00514	0.70833	0.00031	10.0 s
25:16.4	Data point 48	1.50000 mL	0.07173 mL	0.13307 mL	0.02500 mL	12.026	-0.00754	0.76863	0.00042	10.0 s
25:44.2	Data point 49	1.50000 mL	0.07173 mL	0.16905 mL	0.02500 mL	12.216	-0.00132	0.13261	0.00018	10.0 s
27:44.6	Assay volumes	1.75000 mL	0.25035 mL	0.16905 mL	0.02500 mL					

## Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
<b>General Settings</b>				
Analyst name	Dorothy Levorse			
Separate reference vial	Yes			
<b>Standard Experiment Settings</b>				
Number of titrations	1			
Minimum pH	1.800			
Maximum pH	12.200			
pH step between points of	0.200			
Minimum titrant addition	0.00002 mL			
Maximum titrant addition	0.10000 mL			
Argon flow rate	100%			
Start titration using	Cautious pH adjust			
<b>Advanced General Settings</b>				
Detect turbidity using	Spectrometer			
Monitor at a wavelength of	500.0 nm			
Absorbance threshold of	0.100			
Collect turbidity sensor data	No			
Stir after titrant addition for	5 seconds			
For titrant addition, stir at	15%			
<b>Titration Pre-Dose</b>				
Titration pre-dose	None			
<b>Assay Medium</b>				
Cosolvent in use	No			
ISA water volume	1.50 mL			
Water added	Automatic			
After water addition, stir for	5 seconds			
At a speed of	15%			
Buffer in use	Yes			
Buffer type	Phosphate Buffer			

Sample name: **M01**  
 Assay name: **UV-metric pKa**  
 Assay ID: **171-16001**  
 Filename: **C:\Sirius\_T3\Mehtap\20170915\_exp03\_uv\_M01-M14\171-16001\_M01\_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 12:18:58 AM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			
After medium addition, stir for	5 seconds			
<b>Sample Sonication</b>				
Sonicate	No			
<b>Sample Dissolution</b>				
Perform a dissolution stage	No			
<b>Carbonate purge</b>				
Perform a carbonate purge	No			
<b>Temperature Control</b>				
Wait for temperature	Yes			
Required start temperature	25.0°C			
Acceptable deviation	0.5°C			
Time to wait	60 seconds			
Stir speed of	15%			
<b>Titration 1</b>				
Titrate from	Low to high pH			
Adjust to start pH	Yes			
After pH adjust stir for	10 seconds			
<b>Data Point Stability</b>				
Stir during data point collection	Yes			
For point collection, stir at	15%			
Delay before data point collection	0 seconds			
Number of points to average	20 points			
Time interval between points	0.50 seconds			
Required maximum standard deviation	0.00500 dpH/dt			
Stability timeout after	60 seconds			
<b>Experiment cleanup</b>				
Adjust pH to cleanup	To start pH			
And then stir for	60 seconds			
For cleaning, stir at	20%			
Then add water volume	0.25 mL			
And then stir for	30 seconds			

## Calibration Settings

Setting	Value	Date/Time changed	Imported from
Four-Plus alpha	0.112	9/16/2017 12:18:58 AM	C:\Sirius_T3\HCl17115.t3r
Four-Plus S	1.0006	9/16/2017 12:18:58 AM	C:\Sirius_T3\HCl17115.t3r
Four-Plus jH	0.7	9/16/2017 12:18:58 AM	C:\Sirius_T3\HCl17115.t3r
Four-Plus jOH	-0.6	9/16/2017 12:18:58 AM	C:\Sirius_T3\HCl17115.t3r
Base concentration factor	1.015	9/16/2017 12:18:58 AM	C:\Sirius_T3\KOH17111.t3r
Acid concentration factor	1.003	9/16/2017 12:18:58 AM	C:\Sirius_T3\HCl17115.t3r

## Instrument Settings

Setting	Value	Batch Id	Install date
Instrument owner	Merck		
Instrument ID	T311053		
Instrument type	T3 Simulator		
Software version	1.1.3.0		
Dispenser module		T3DM1100253	3/31/2009 6:24:52 AM
Dispenser 0	Water		3/31/2009 6:25:05 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Water (0.15 M KCl)	8-18-17	9/8/2017 9:22:43 AM
Dispenser 2	Acid		3/31/2009 6:25:11 AM
Syringe volume	0.5 mL		

Sample name: **M01**  
 Assay name: **UV-metric pKa**  
 Assay ID: **171-16001**  
 Filename: **C:\Sirius\_T3\Mehtap\20170915\_exp03\_uv\_M01-M14\171-16001\_M01\_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 12:18:58 AM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Instrument Settings (continued)

Setting	Value	Batch Id	Install date
Firmware version	1.2.1(r2)		
Titrant	Acid (0.5 M HCl)	166940	9/8/2017 9:21:27 AM
Dispenser 1	Base		3/31/2009 6:25:21 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Base (0.5 M KOH)	01/06/17	9/8/2017 9:20:03 AM
Dispenser 5	Cosolvent		3/31/2009 6:26:24 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Distribution valve 5	Distribution Valve		3/31/2009 6:28:19 AM
Firmware version	1.1.3		
Port A	Methanol (80%, 0.15 M KCl)	8-15-17	9/13/2017 12:23:11 PM
Dispenser 3	Buffer		8/3/2010 6:05:16 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Phosphate Buffer		9/12/2017 12:32:29 PM
Dispenser 6	Octanol		10/22/2010 11:52:43 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Octanol	9-14-17	9/14/2017 10:30:38 AM
Titrator		T3TM1100153	3/31/2009 6:24:17 AM
Horizontal axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Probe I/O firmware version	1.1.1		
Electrode	T3 Electrode	T3E0769	8/15/2017 10:21:54 AM
E0 calibration	-8.54 mV		9/16/2017 12:19:22 AM
Filling solution	3M KCl	KCL095	9/13/2017 9:16:19 AM
Liquids			
Wash 1	50% IPA:50% Water		9/15/2017 9:38:18 AM
Wash 2	0.5% Triton X-100 in H2O		9/15/2017 9:38:22 AM
Buffer position 1	pH7 Wash		9/15/2017 9:38:24 AM
Buffer position 2	pH 7		9/15/2017 9:38:27 AM
Storage position			9/15/2017 9:38:55 AM
Wash water	3.7e+003 mL	9-11-17	9/11/2017 4:28:43 PM
Waste	6.4e+003 mL		9/11/2017 4:28:49 PM
Temperature controller			8/5/2010 7:35:13 AM
Turbidity detector			3/31/2009 6:24:45 AM
Spectrometer		072390	11/23/2010 12:22:28 PM
Dip probe		11086	
Wavelength coefficient A0	185.563		
Wavelength coefficient A1	2.17439		
Wavelength coefficient A2	-0.000285622		
Total lamp lit time	114:03:31		11/23/2010 12:22:28 PM
Calibrated on	9/6/2017 9:33:02 AM		
Integration time	11		
Scans averaged	10		
Autoloader		T3AL1100237	11/10/2015 10:34:13 AM
Left-right axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Front-back axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Configuration			
Alternate titration position	Titration position		
Alternate reference position	Reference position		
Maximum standard vial volume	3.50 mL		
Maximum alternate vial volume	25.00 mL		
Automatic action idle period	5 minute(s)		





## Assay Settings

Sample name: **M01** Experiment start time: **9/16/2017 12:18:58 AM**  
Assay name: **UV-metric pKa** Analyst: **Dorothy Levorse**  
Assay ID: **17I-16001** Instrument ID: **T311053**  
Filename: **C:\Sirius\_T3\Mehtap\20170915\_exp03\_uv\_M01-M14\17I-16001\_M01\_UV-metric pKa.t3r**

### Instrument Settings (continued)

Setting	Value	Batch Id	Install date
Titrant tube volume	1.3 mL		
Syringe flush count	3.50		
Flowing wash pump volume	20.0 mL		
Flowing wash stir duration	5 s		
Flowing wash stir speed	30%		
Solvent wash stir duration	5 s		
Solvent wash stir speed	30%		
Surfactant wash stir duration	5 s		
Surfactant wash stir speed	30%		
E0 calibration minimum number of points	10		
E0 calibration maximum standard deviation	0.01500		
E0 calibration timeout period	60 s		
E0 calibration stir duration	5 s		
E0 calibration preparation stir speed	30%		
E0 calibration buffer wash stir duration	5 s		
E0 calibration buffer wash stir speed	30%		
E0 calibration reading stir speed	0%		
Spectrometer calibration stir duration	5 s		
Spectrometer calibration stir speed	30%		
Spectrometer calibration wash pump volume	20.0 mL		
Spectrometer calibration wash stir duration	5 s		
Spectrometer calibration wash stir speed	30%		
Overhead dispense height	10000		

### Refinement Settings

Setting	Value	Default value
Turbidity detection method	Spectrometer	Spectrometer
Turbidity wavelength to assess	500.0 nm	500.0 nm
Turbidity maximum absorbance	0.100	0.100
Turbidity probe threshold	50.00	50.00
Exclude turbid points	Yes	Yes
Low intensity warning threshold	100	100
Minimum absorbance change threshold	0.100	0.100
Eigenvector autocorrelation threshold	0.80	0.80
Maximum RMSD severe warning	0.250	0.250
Maximum RMSD warning	0.050	0.050