

Sample name: **M19**  
 Assay name: **UV-metric pKa**  
 Assay ID: **17J-02024**  
 Filename: **C:\Sirius\_T3\Mehtap\20171002\_exp11\_pKa\17J-02024\_M19\_UV-metric pKa.t3r**

Experiment start time: **10/2/2017 7:55:10 PM**  
 Analyst: **Dorothy Leverse**  
 Instrument ID: **T311053**

## Results

pKa 1 **3.15**  
 RMSD **0.002 0.002**  
 Chi squared **0.0085**  
 PCA calculated number of pKas **2**  
 Average ionic strength **0.158 M**  
 Average temperature **24.9°C**  
 Analyte concentration range **80.9 µM to 73.2 µM**

Number of pKas source **Predicted**  
 Wavelength clipping **230.0 nm to 450.0 nm**  
 pH clipping **1.276 to 12.726**


## Warnings and errors

Errors None  
 Warnings PCA calculation disagrees with predicted number of pKas

## Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
<b>Assay Medium</b>				
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			


## Graphs

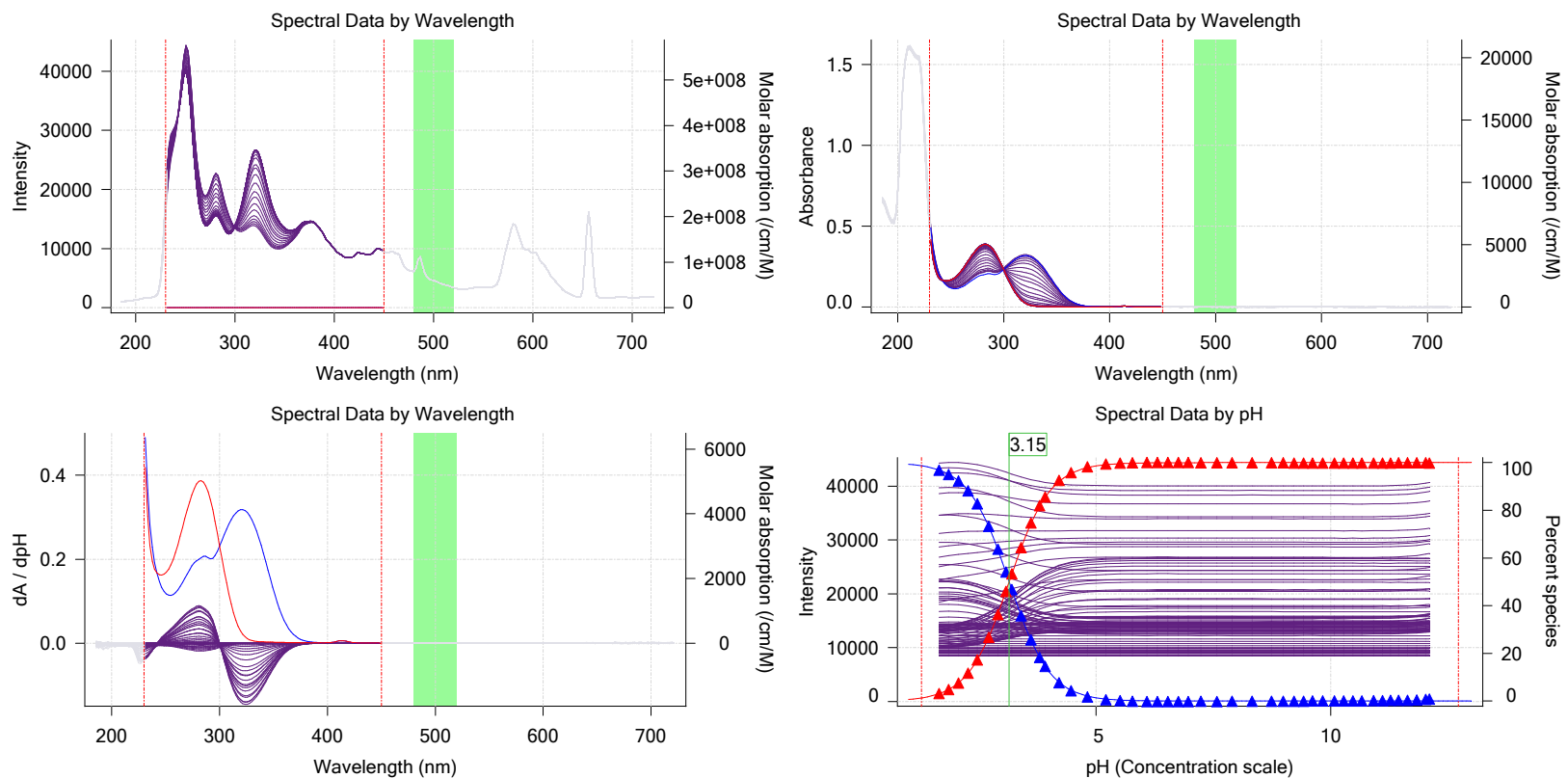

Sample name: **M19**  
 Assay name: **UV-metric pKa**  
 Assay ID: **17J-02024**  
 Filename: **C:\Sirius\_T3\Mehtap\20171002\_exp11\_pKa\17J-02024\_M19\_UV-metric pKa.t3r**

Experiment start time: **10/2/2017 7:55:10 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Graphs (continued)

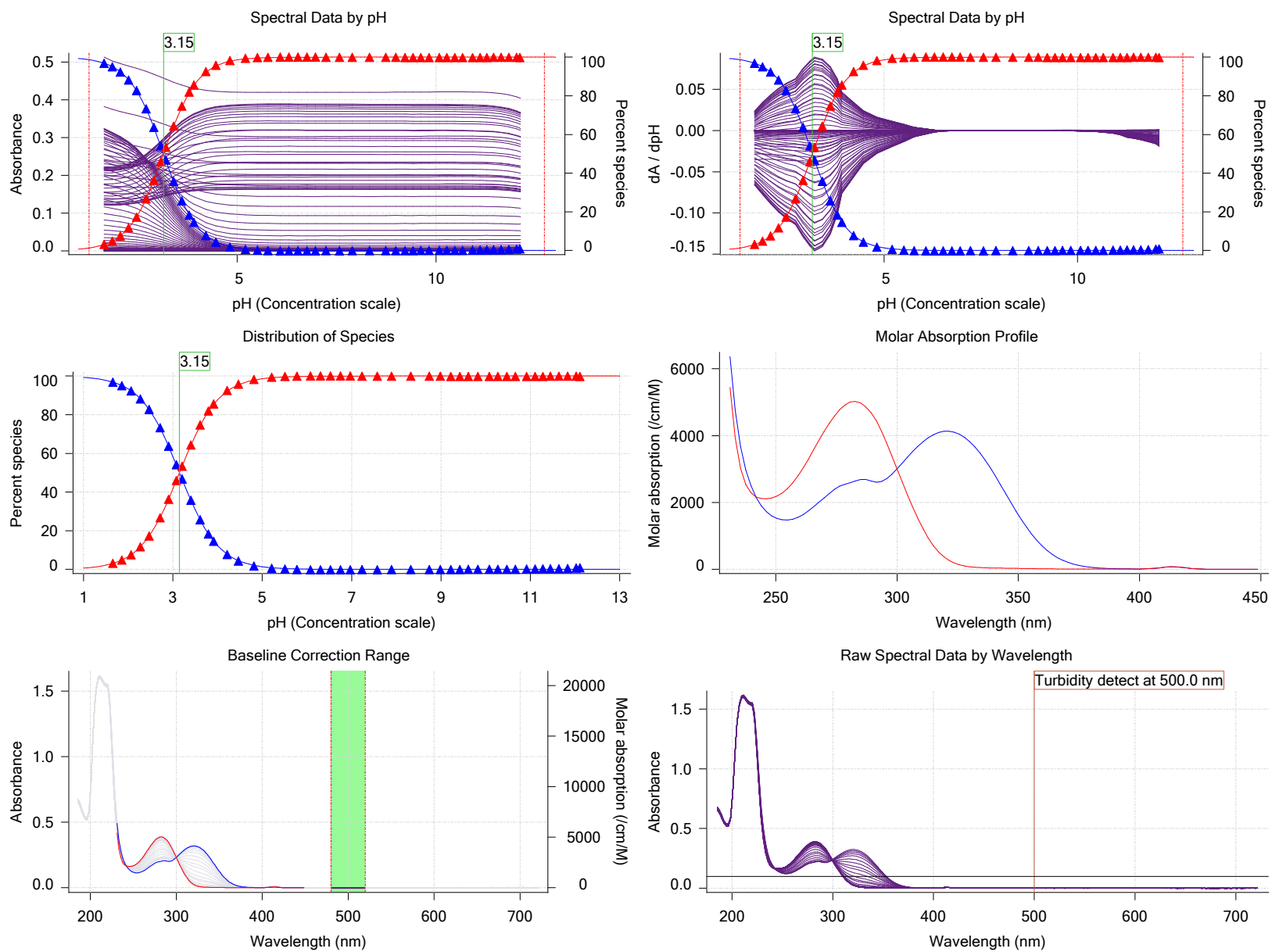

## Assay Model

### Settings

Settings	Value	Date/Time changed	Imported from
Sample name	M19	9/29/2017 6:36:04 PM	User entered value
Sample by	Volume		Default value
Sample volume	0.0020 mL	9/29/2017 6:36:04 PM	User entered value
Solvent	DMSO		Default value
Sample concentration	0.064600 M	10/2/2017 3:11:22 PM	User entered value
Solubility	Unknown		Default value
Molecular weight	269.32	9/29/2017 6:36:11 PM	User entered value
Individual pKa ionic environments	No		Default value
Number of pKas	1	9/29/2017 6:36:04 PM	User entered value
Sample is a	Base	9/29/2017 6:36:04 PM	User entered value
pKa 1	0.99	9/29/2017 6:36:04 PM	User entered value
logp (XH +)	-10.00		Default value
logP (neutral X)	-10.00	9/29/2017 6:36:04 PM	User entered value

Sample name: **M19** Experiment start time: **10/2/2017 7:55:10 PM**  
 Assay name: **UV-metric pKa** Analyst: **Dorothy Levorse**  
 Assay ID: **17J-02024** Instrument ID: **T311053**  
 Filename: **C:\Sirius\_T3\Mehtap\20171002\_exp11\_pKa\17J-02024\_M19\_UV-metric pKa.t3r**

## Events

Time	Event	Water	Acid	Base	Buffer	pH	dpH/dt	pH R-squared	pH SD
3:15.3	Dark spectrum								
3:16.7	Reference spectrum								
3:44.3	Volume reset due to vial change								
5:14.6	Initial pH = 7.70								
6:27.6	Data point 4	1.50000 mL	0.07048 mL	0.00000 mL	0.02500 mL	1.776	-0.00971	0.69270	0.000
6:56.4	Data point 5	1.50000 mL	0.07048 mL	0.02573 mL	0.02500 mL	1.978	0.00934	0.58344	0.000
7:13.3	Data point 6	1.50000 mL	0.07048 mL	0.04236 mL	0.02500 mL	2.182	0.01395	0.84699	0.000
7:30.1	Data point 7	1.50000 mL	0.07048 mL	0.05270 mL	0.02500 mL	2.389	0.00863	0.79531	0.000
7:46.7	Data point 8	1.50000 mL	0.07048 mL	0.05915 mL	0.02500 mL	2.584	0.01169	0.79530	0.000
8:03.4	Data point 9	1.50000 mL	0.07048 mL	0.06326 mL	0.02500 mL	2.828	0.00804	0.78537	0.000
8:30.3	Data point 10	1.50000 mL	0.07048 mL	0.06562 mL	0.02500 mL	3.022	0.00361	0.47395	0.000
8:46.8	Data point 11	1.50000 mL	0.07048 mL	0.06714 mL	0.02500 mL	3.195	0.00908	0.73478	0.000
9:03.2	Data point 12	1.50000 mL	0.07048 mL	0.06816 mL	0.02500 mL	3.322	0.00802	0.68399	0.000
9:24.9	Data point 13	1.50000 mL	0.07048 mL	0.06893 mL	0.02500 mL	3.519	0.01407	0.85553	0.000
9:41.4	Data point 14	1.50000 mL	0.07048 mL	0.06940 mL	0.02500 mL	3.728	0.01634	0.90995	0.000
9:57.9	Data point 15	1.50000 mL	0.07048 mL	0.06971 mL	0.02500 mL	3.912	0.02914	0.91510	0.001
10:14.4	Data point 16	1.50000 mL	0.07048 mL	0.06990 mL	0.02500 mL	4.029	0.03610	0.94295	0.001
10:36.1	Data point 17	1.50000 mL	0.07048 mL	0.07011 mL	0.02500 mL	4.331	0.06468	0.92988	0.003
10:57.7	Data point 18	1.50000 mL	0.07048 mL	0.07023 mL	0.02500 mL	4.585	0.09671	0.97562	0.004
11:21.4	Data point 19	1.50000 mL	0.07048 mL	0.07030 mL	0.02500 mL	4.934	0.09755	0.98235	0.004
11:51.1	Data point 20	1.50000 mL	0.07048 mL	0.07034 mL	0.02500 mL	5.328	0.09972	0.98673	0.004
12:23.6	Data point 21	1.50000 mL	0.07048 mL	0.07039 mL	0.02500 mL	5.630	0.09042	0.97228	0.004
12:46.0	Data point 22	1.50000 mL	0.07048 mL	0.07044 mL	0.02500 mL	5.897	0.06182	0.92519	0.003
13:12.6	Data point 23	1.50000 mL	0.07048 mL	0.07051 mL	0.02500 mL	6.202	0.06550	0.86922	0.003
13:39.2	Data point 24	1.50000 mL	0.07048 mL	0.07058 mL	0.02500 mL	6.425	0.07203	0.89919	0.003
14:05.7	Data point 25	1.50000 mL	0.07048 mL	0.07067 mL	0.02500 mL	6.651	0.03581	0.81983	0.001
14:32.5	Data point 26	1.50000 mL	0.07048 mL	0.07077 mL	0.02500 mL	6.875	0.06171	0.85162	0.003
14:59.0	Data point 27	1.50000 mL	0.07048 mL	0.07084 mL	0.02500 mL	7.089	0.08363	0.88832	0.004
15:25.6	Data point 28	1.50000 mL	0.07048 mL	0.07091 mL	0.02500 mL	7.356	0.08563	0.88949	0.004
15:53.3	Data point 29	1.50000 mL	0.07048 mL	0.07098 mL	0.02500 mL	7.695	0.08671	0.87838	0.004
16:17.5	Data point 30	1.50000 mL	0.07048 mL	0.07103 mL	0.02500 mL	8.018	0.09217	0.87206	0.004
16:45.0	Data point 31	1.50000 mL	0.07048 mL	0.07107 mL	0.02500 mL	8.450	0.08031	0.75053	0.004
17:12.1	Data point 32	1.50000 mL	0.07048 mL	0.07112 mL	0.02500 mL	8.863	0.09451	0.89525	0.004
17:37.1	Data point 33	1.50000 mL	0.07048 mL	0.07117 mL	0.02500 mL	9.111	0.09101	0.87835	0.004
18:00.5	Data point 34	1.50000 mL	0.07048 mL	0.07124 mL	0.02500 mL	9.344	0.07685	0.77851	0.004
18:22.2	Data point 35	1.50000 mL	0.07048 mL	0.07131 mL	0.02500 mL	9.547	0.04349	0.82739	0.002
18:43.7	Data point 36	1.50000 mL	0.07048 mL	0.07140 mL	0.02500 mL	9.741	0.01204	0.48248	0.000
19:10.6	Data point 37	1.50000 mL	0.07048 mL	0.07157 mL	0.02500 mL	9.954	0.00007	0.00006	0.000
19:27.2	Data point 38	1.50000 mL	0.07048 mL	0.07183 mL	0.02500 mL	10.209	-0.00882	0.64883	0.000
19:58.8	Data point 39	1.50000 mL	0.07048 mL	0.07234 mL	0.02500 mL	10.401	-0.01399	0.88215	0.000
20:30.7	Data point 40	1.50000 mL	0.07048 mL	0.07312 mL	0.02500 mL	10.601	-0.01542	0.92853	0.000
21:02.6	Data point 41	1.50000 mL	0.07048 mL	0.07437 mL	0.02500 mL	10.794	-0.01496	0.85278	0.000
21:34.5	Data point 42	1.50000 mL	0.07048 mL	0.07728 mL	0.02500 mL	11.044	-0.00961	0.78726	0.000
22:01.3	Data point 43	1.50000 mL	0.07048 mL	0.08010 mL	0.02500 mL	11.235	-0.01346	0.83601	0.000
22:17.9	Data point 44	1.50000 mL	0.07048 mL	0.08438 mL	0.02500 mL	11.406	-0.01117	0.86267	0.000
22:34.6	Data point 45	1.50000 mL	0.07048 mL	0.09076 mL	0.02500 mL	11.574	-0.01112	0.85817	0.000
23:01.7	Data point 46	1.50000 mL	0.07048 mL	0.10207 mL	0.02500 mL	11.764	-0.01189	0.88000	0.000
23:28.9	Data point 47	1.50000 mL	0.07048 mL	0.12053 mL	0.02500 mL	11.956	-0.00770	0.76880	0.000
23:56.4	Data point 48	1.50000 mL	0.07048 mL	0.15028 mL	0.02500 mL	12.147	-0.00576	0.61762	0.000
24:18.7	Data point 49	1.50000 mL	0.07048 mL	0.16773 mL	0.02500 mL	12.226	-0.00321	0.30780	0.000
26:18.9	Assay volumes	1.75000 mL	0.24621 mL	0.16773 mL	0.02500 mL				

## Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
<b>General Settings</b>				
Analyst name	Dorothy Levorse			

Sample name: **M19**  
 Assay name: **UV-metric pKa**  
 Assay ID: **17J-02024**  
 Filename: **C:\Sirius\_T3\Mehtap\20171002\_exp11\_pKa\17J-02024\_M19\_UV-metric pKa.t3r**

Experiment start time: **10/2/2017 7:55:10 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
Separate reference vial	Yes			
<b>Standard Experiment Settings</b>				
Number of titrations	1			
Minimum pH	1.800			
Maximum pH	12.200			
pH step between points of	0.200			
Minimum titrant addition	0.00002 mL			
Maximum titrant addition	0.10000 mL			
Argon flow rate	100%			
Start titration using	Cautious pH adjust			
<b>Advanced General Settings</b>				
Detect turbidity using	Spectrometer			
Monitor at a wavelength of	500.0 nm			
Absorbance threshold of	0.100			
Collect turbidity sensor data	No			
Stir after titrant addition for	5 seconds			
For titrant addition, stir at	15%			
<b>Titration Pre-Dose</b>				
Titration pre-dose	None			
<b>Assay Medium</b>				
Cosolvent in use	No			
ISA water volume	1.50 mL			
Water added	Automatic			
After water addition, stir for	5 seconds			
At a speed of	15%			
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			
After medium addition, stir for	5 seconds			
<b>Sample Sonication</b>				
Sonicate	No			
<b>Sample Dissolution</b>				
Perform a dissolution stage	No			
<b>Carbonate purge</b>				
Perform a carbonate purge	No			
<b>Temperature Control</b>				
Wait for temperature	Yes			
Required start temperature	25.0°C			
Acceptable deviation	0.5°C			
Time to wait	60 seconds			
Stir speed of	15%			
<b>Titration 1</b>				
Titrate from	Low to high pH			
Adjust to start pH	Yes			
After pH adjust stir for	10 seconds			
<b>Data Point Stability</b>				
Stir during data point collection	Yes			
For point collection, stir at	15%			
Delay before data point collection	0 seconds			
Number of points to average	20 points			
Time interval between points	0.50 seconds			
Required maximum standard deviation	0.00500 dpH/dt			
Stability timeout after	60 seconds			
<b>Experiment cleanup</b>				
Adjust pH to cleanup	To start pH			
And then stir for	60 seconds			
For cleaning, stir at	20%			



## Assay Settings

Sample name: **M19** Experiment start time: **10/2/2017 7:55:10 PM**  
Assay name: **UV-metric pKa** Analyst: **Dorothy Levorse**  
Assay ID: **17J-02024** Instrument ID: **T311053**  
Filename: **C:\Sirius\_T3\Mehtap\20171002\_exp11\_pKa\17J-02024\_M19\_UV-metric pKa.t3r**

## Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
Then add water volume	0.25 mL			
And then stir for	30 seconds			

## Calibration Settings

Setting	Value	Date/Time changed	Imported from
Four-Plus alpha	0.114	10/2/2017 7:55:10 PM	C:\Sirius_T3\17J-02006_Blank standardisation.t3r
Four-Plus S	1.0012	10/2/2017 7:55:10 PM	C:\Sirius_T3\17J-02006_Blank standardisation.t3r
Four-Plus jH	0.4	10/2/2017 7:55:10 PM	C:\Sirius_T3\17J-02006_Blank standardisation.t3r
Four-Plus jOH	-0.5	10/2/2017 7:55:10 PM	C:\Sirius_T3\17J-02006_Blank standardisation.t3r
Base concentration factor	1.011	10/2/2017 7:55:10 PM	C:\Sirius_T3\KOH17I22.t3r
Acid concentration factor	1.018	10/2/2017 7:55:10 PM	C:\Sirius_T3\17J-02006_Blank standardisation.t3r

## Instrument Settings

Setting	Value	Batch Id	Install date
Instrument owner	Merck		
Instrument ID	T311053		
Instrument type	T3 Simulator		
Software version	1.1.3.0		
Dispenser module		T3DM1100253	3/31/2009 6:24:52 AM
Dispenser 0	Water		3/31/2009 6:25:05 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Water (0.15 M KCl)	8-18-17	9/26/2017 9:05:04 AM
Dispenser 2	Acid		3/31/2009 6:25:11 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Acid (0.5 M HCl)	166940	9/8/2017 9:21:27 AM
Dispenser 1	Base		3/31/2009 6:25:21 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Base (0.5 M KOH)	9-22-17	9/22/2017 4:02:42 PM
Dispenser 5	Cosolvent		3/31/2009 6:26:24 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Distribution valve 5	Distribution Valve		3/31/2009 6:28:19 AM
Firmware version	1.1.3		
Port A	Methanol (80%, 0.15 M KCl)	9-26-17	9/29/2017 9:58:40 AM
Port B	Cyclohexane		9/19/2017 2:15:02 PM
Port C	MeCN (50%, 0.15 M KCl)	10-2-17	10/2/2017 11:28:55 AM
Dispenser 3	Buffer		8/3/2010 6:05:16 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Phosphate Buffer		9/12/2017 12:32:29 PM
Dispenser 6	Octanol		10/22/2010 11:52:43 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Octanol	9-14-17	9/14/2017 10:30:38 AM
Titrator		T3TM1100153	3/31/2009 6:24:17 AM
Horizontal axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Probe I/O firmware version	1.1.1		
Electrode	T3 Electrode	T3E0769	8/15/2017 10:21:54 AM
E0 calibration	-8.84 mV		10/2/2017 7:55:34 PM
Filling solution	3M KCl	KCL095	10/2/2017 9:26:59 AM
Liquids			

Sample name: **M19** Experiment start time: **10/2/2017 7:55:10 PM**  
 Assay name: **UV-metric pKa** Analyst: **Dorothy Levorse**  
 Assay ID: **17J-02024** Instrument ID: **T311053**  
 Filename: **C:\Sirius\_T3\Mehtap\20171002\_exp11\_pKa\17J-02024\_M19\_UV-metric pKa.t3r**

## Instrument Settings (continued)

Setting	Value	Batch Id	Install date
Wash 1	50% IPA:50% Water		10/2/2017 9:38:49 AM
Wash 2	0.5% Triton X-100 in H2O		10/2/2017 9:38:52 AM
Buffer position 1	pH7 Wash		10/2/2017 9:38:54 AM
Buffer position 2	pH 7		10/2/2017 9:38:57 AM
Storage position			10/2/2017 9:36:04 AM
Wash water	3e+003 mL	9-27-17	9/27/2017 4:24:06 PM
Waste	7e+003 mL		9/27/2017 4:24:14 PM
Temperature controller			8/5/2010 7:35:13 AM
Turbidity detector			3/31/2009 6:24:45 AM
Spectrometer		072390	11/23/2010 12:22:28 PM
Dip probe		11086	
Wavelength coefficient A0	185.563		
Wavelength coefficient A1	2.17439		
Wavelength coefficient A2	-0.000285622		
Total lamp lit time	313:32:06		11/23/2010 12:22:28 PM
Calibrated on	9/26/2017 9:22:07 AM		
Integration time	11		
Scans averaged	10		
Autoloader		T3AL1100237	11/10/2015 10:34:13 AM
Left-right axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Front-back axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Configuration			
Alternate titration position	Titration position		
Alternate reference position	Reference position		
Maximum standard vial volume	3.50 mL		
Maximum alternate vial volume	25.00 mL		
Automatic action idle period	5 minute(s)		
Titrant tube volume	1.3 mL		
Syringe flush count	3.50		
Flowing wash pump volume	20.0 mL		
Flowing wash stir duration	5 s		
Flowing wash stir speed	30%		
Solvent wash stir duration	5 s		
Solvent wash stir speed	30%		
Surfactant wash stir duration	5 s		
Surfactant wash stir speed	30%		
E0 calibration minimum number of points	10		
E0 calibration maximum standard deviation	0.01500		
E0 calibration timeout period	60 s		
E0 calibration stir duration	5 s		
E0 calibration preparation stir speed	30%		
E0 calibration buffer wash stir duration	5 s		
E0 calibration buffer wash stir speed	30%		
E0 calibration reading stir speed	0%		
Spectrometer calibration stir duration	5 s		
Spectrometer calibration stir speed	30%		
Spectrometer calibration wash pump volume	20.0 mL		
Spectrometer calibration wash stir duration	5 s		
Spectrometer calibration wash stir speed	30%		
Overhead dispense height	10000		

## Refinement Settings

Setting	Value	Default value
Turbidity detection method	Spectrometer	Spectrometer
Turbidity wavelength to assess	500.0 nm	500.0 nm





Sample name: **M19**

Assay name: **UV-metric pKa**

Assay ID: **17J-02024**

Filename: **C:\Sirius\_T3\Mehtap\20171002\_exp11\_pKa\17J-02024\_M19\_UV-metric pKa.t3r**

Experiment start time: **10/2/2017 7:55:10 PM**

Analyst: **Dorothy Levorse**

Instrument ID: **T311053**

Refinement Settings (continued)

Setting	Value	Default value
Turbidity maximum absorbance	0.100	0.100
Turbidity probe threshold	50.00	50.00
Exclude turbid points	Yes	Yes
Low intensity warning threshold	100	100
Minimum absorbance change threshold	0.100	0.100
Eigenvector autocorrelation threshold	0.80	0.80
Maximum RMSD severe warning	0.250	0.250
Maximum RMSD warning	0.050	0.050

Tray Information

Title  
Location H1