

Sample name: **D01**  
 Assay name: **UV-metric psKa**  
 Assay ID: **171-28003**  
 Filename: **C:\Sirius\_T3\Mehtap\20170928\_exp09\_uv\_pKa\171-28003\_D01\_UV-metric psKa.t3r**

Experiment start time: **9/28/2017 4:59:03 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**




## Yasuda-Shedlovsky result

Extrapolation type	pKa 0%	SD	Intercept	Slope	R <sup>2</sup>	Ionic strength	Temperature
Yasuda-Shedlovsky	2.19	±0.09	5.62	-132.3950	0.9829	0.167 M	25.0°C
Yasuda-Shedlovsky	9.63	±0.08	11.63	-19.9043	0.6090	0.167 M	25.0°C
Yasuda-Shedlovsky	11.11	±0.12	12.76	7.3790	0.0904	0.167 M	25.0°C

## Component assay results

Titration	Methanol weight%	Direction	Result type	Dielectric constant	[H2O]	Ionic strength	Temperature	psKa 1	psKa 2
171-28003 Points 4 to 32	59.26 %	Up	UV-metric pKa	52.0	19.5 M	0.158 M	25.0°C	✓	1.80 ✓
171-28003 Points 34 to 72	49.27 %	Up	UV-metric pKa	56.7	24.8 M	0.168 M	25.0°C	✓	1.86 ✓
171-28003 Points 74 to 112	39.51 %	Up	UV-metric pKa	61.2	30.3 M	0.176 M	25.0°C	✓	1.99 ✓

## Graphs

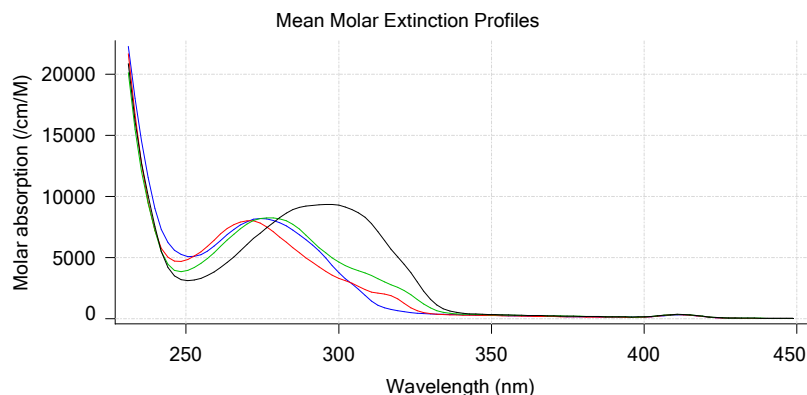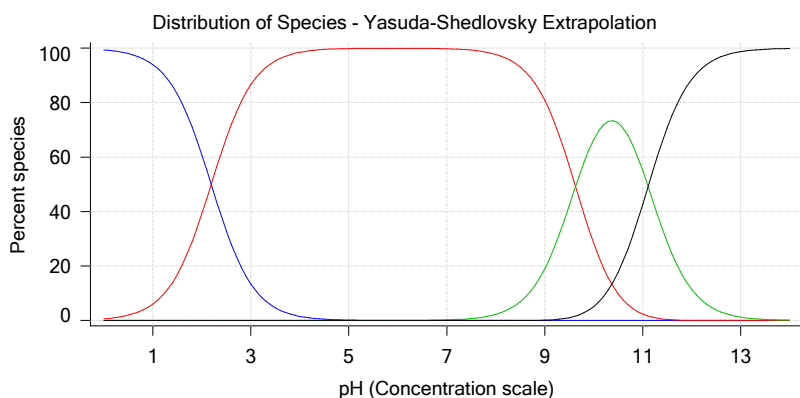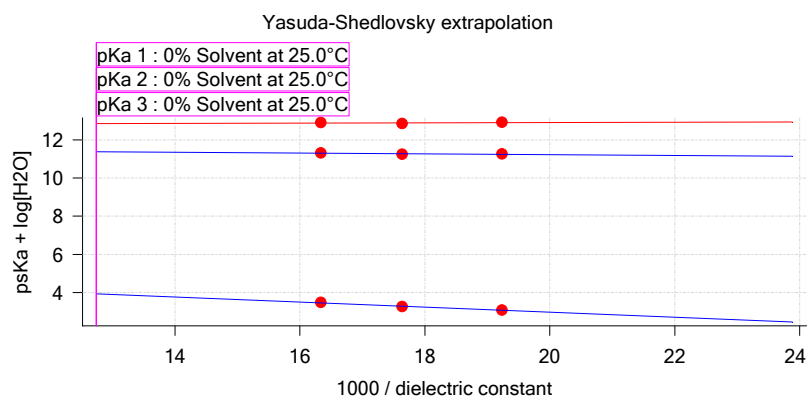

## UV-metric psKa Titration 1 of 3 171-28003 Points 4 to 32

## Results

pKa 1 **1.80**  
 pKa 2 **9.97**  
 pKa 3 **11.63**  
 RMSD **0.002 0.005 0.008 0.001**  
 Chi squared **0.0407**  
 PCA calculated number of pKas **3**  
 Average ionic strength **0.158 M**  
 Average temperature **25.0°C**  
 Analyte concentration range **39.9 µM to 37.6 µM**  
 Methanol weight % **59.3 %**  
 Dielectric constant **52.0**

Sample name: **D01**  
 Assay name: **UV-metric psKa**  
 Assay ID: **171-28003**  
 Filename: **C:\Sirius\_T3\Mehtap\20170928\_exp09\_uv\_pKa\171-28003\_D01\_UV-metric psKa.t3r**

Experiment start time: **9/28/2017 4:59:03 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Results (continued)

Water concentration **19.5 M**

Number of pKas source **Manual (3)**  
 Wavelength clipping **230.0 nm to 450.0 nm**  
 pH clipping **1.469 to 12.553**


## Warnings and errors

Errors **None**  
 Warnings **PCA calculation disagrees with predicted number of pKas**

## Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
<b>Assay Medium</b>				
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			







## Graphs

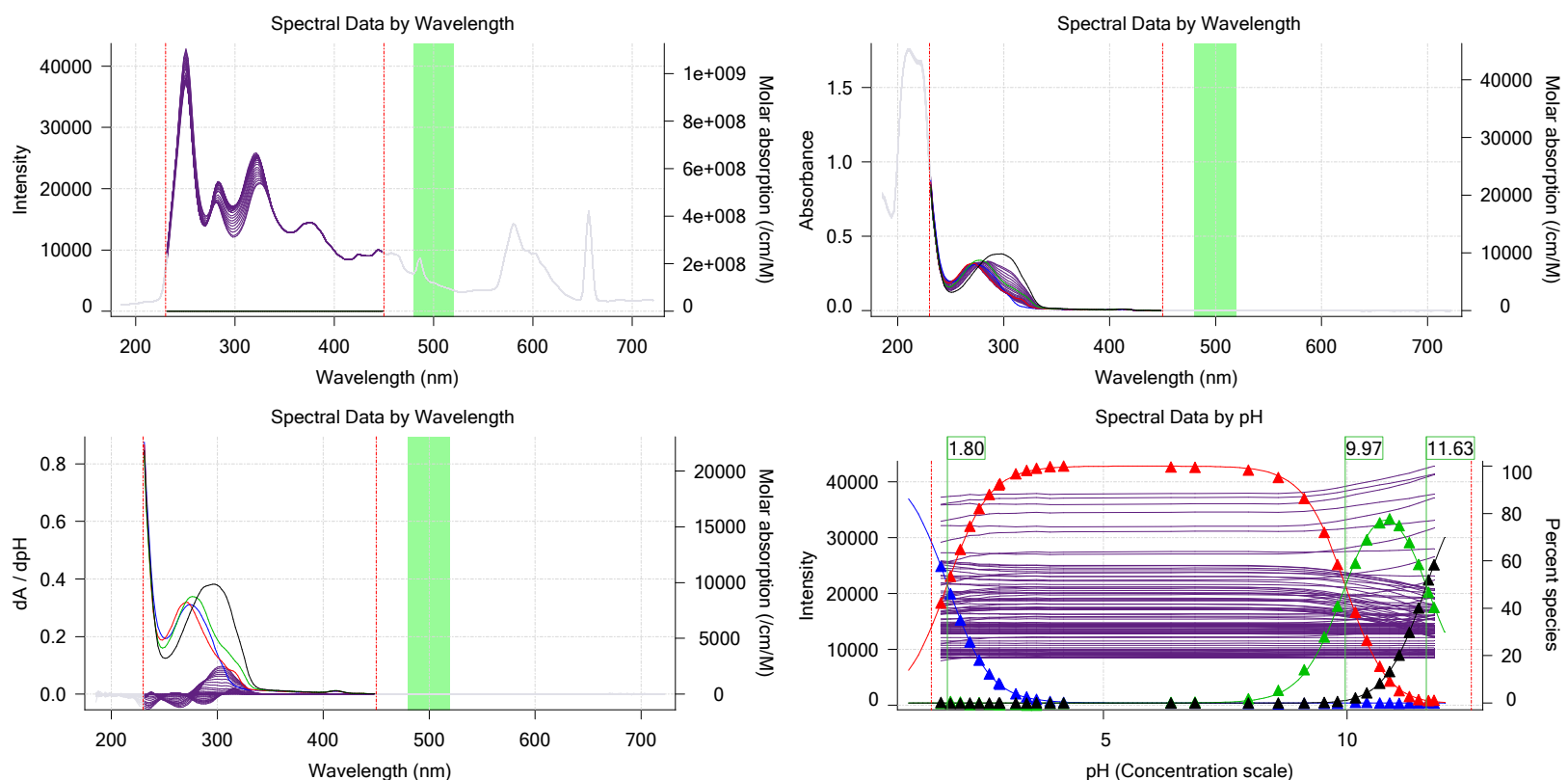

Sample name: **D01**  
 Assay name: **UV-metric psKa**  
 Assay ID: **171-28003**  
 Filename: **C:\Sirius\_T3\Mehtap\20170928\_exp09\_uv\_pKa\171-28003\_D01\_UV-metric psKa.t3r**

Experiment start time: **9/28/2017 4:59:03 PM**  
 Analyst: **Dorothy Levorso**  
 Instrument ID: **T311053**

## Graphs (continued)

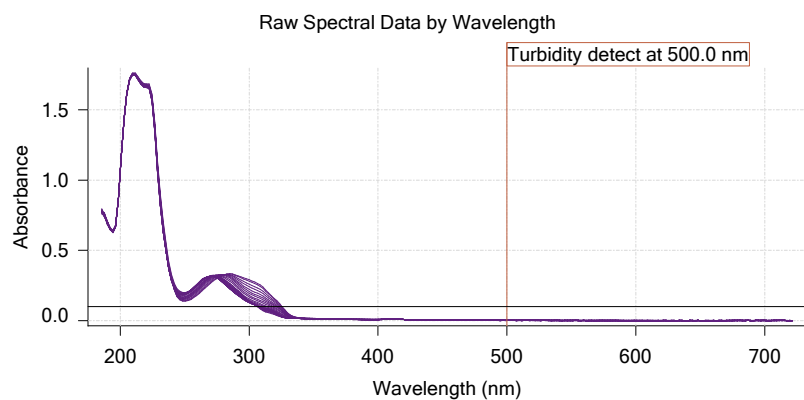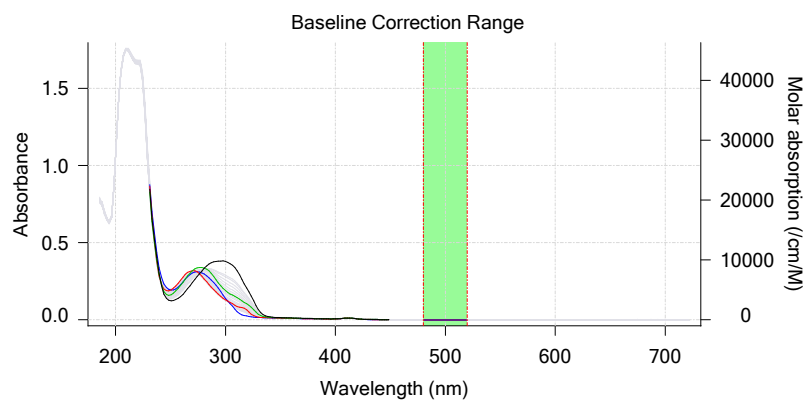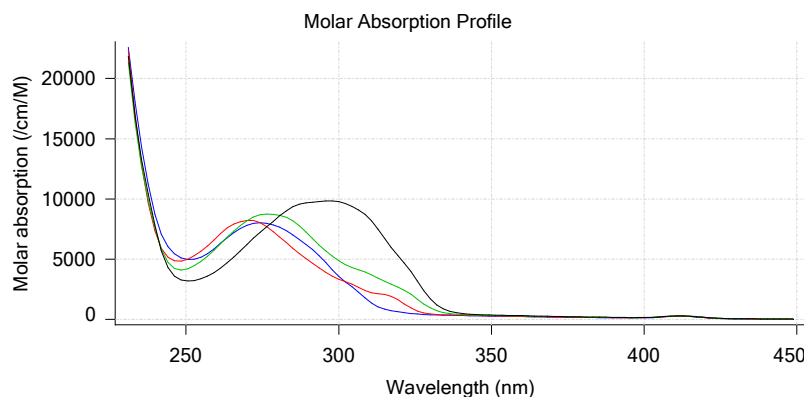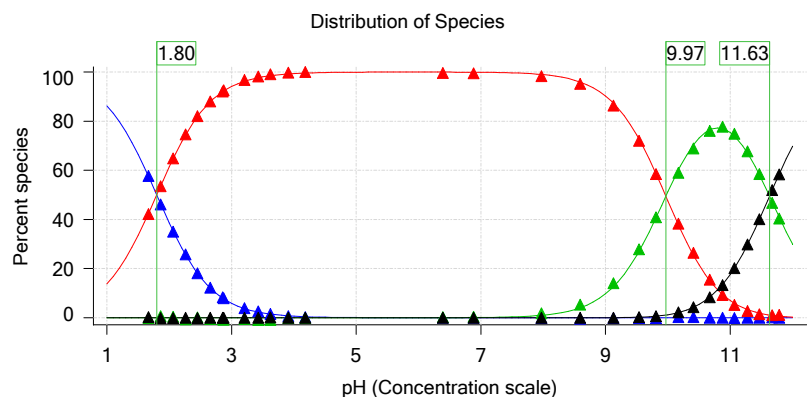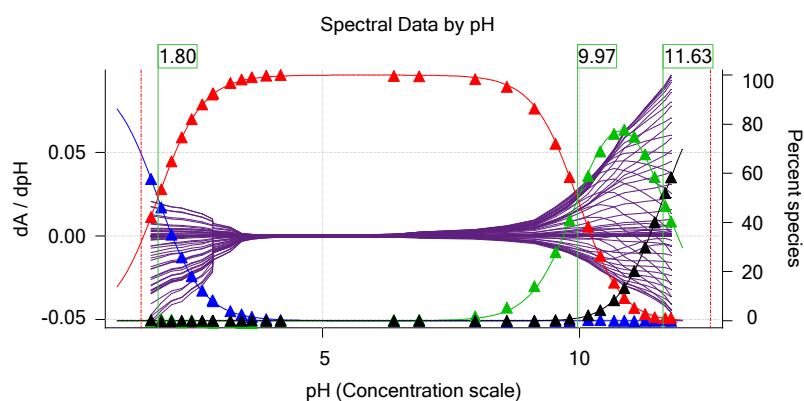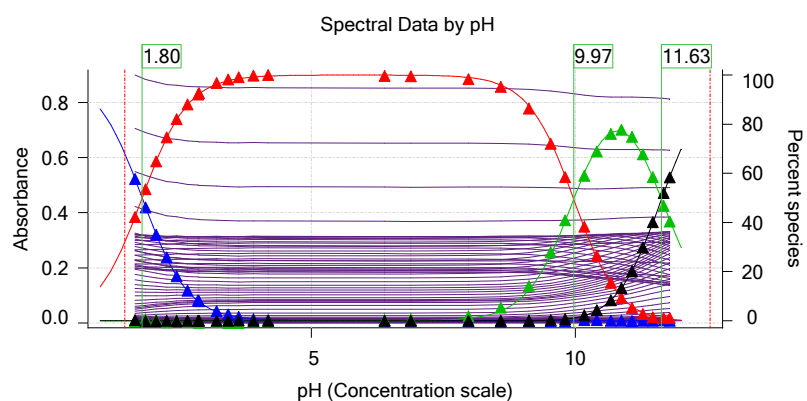

## UV-metric psKa Titration 2 of 3 171-28003 Points 34 to 72

### Results

pKa 1 **1.86**  
 pKa 2 **9.85**  
 pKa 3 **11.46**  
 RMSD **0.002 0.005 0.005 0.001**  
 Chi squared **0.0167**  
 PCA calculated number of pKas **3**  
 Average ionic strength **0.168 M**  
 Average temperature **25.0°C**  
 Analyte concentration range **34.0 µM to 32.0 µM**  
 Methanol weight % **49.3 %**  
 Dielectric constant **56.7**

Sample name: **D01**  
 Assay name: **UV-metric psKa**  
 Assay ID: **171-28003**  
 Filename: **C:\Sirius\_T3\Mehtap\20170928\_exp09\_uv\_pKa\171-28003\_D01\_UV-metric psKa.t3r**

Experiment start time: **9/28/2017 4:59:03 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Results (continued)

Water concentration **24.8 M**  
 Number of pKas source **Manual (3)**  
 Wavelength clipping **230.0 nm to 450.0 nm**  
 pH clipping **1.443 to 12.539**


## Warnings and errors

Errors **None**  
 Warnings **PCA calculation disagrees with predicted number of pKas**

## Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
<b>Assay Medium</b>				
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			


## Graphs

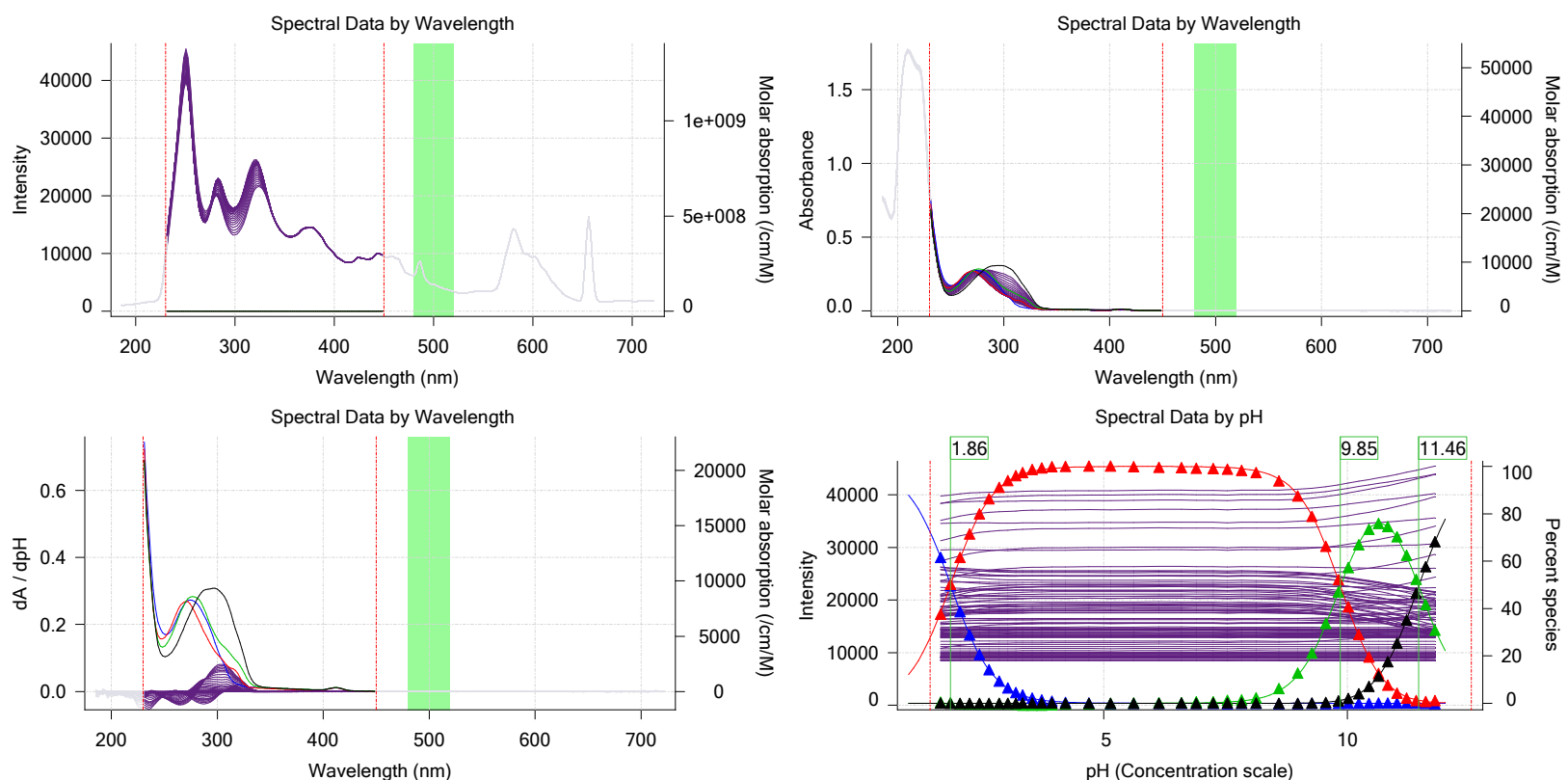

Sample name: **D01**  
 Assay name: **UV-metric psKa**  
 Assay ID: **171-28003**  
 Filename: **C:\Sirius\_T3\Mehtap\20170928\_exp09\_uv\_pKa\171-28003\_D01\_UV-metric psKa.t3r**

Experiment start time: **9/28/2017 4:59:03 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Graphs (continued)

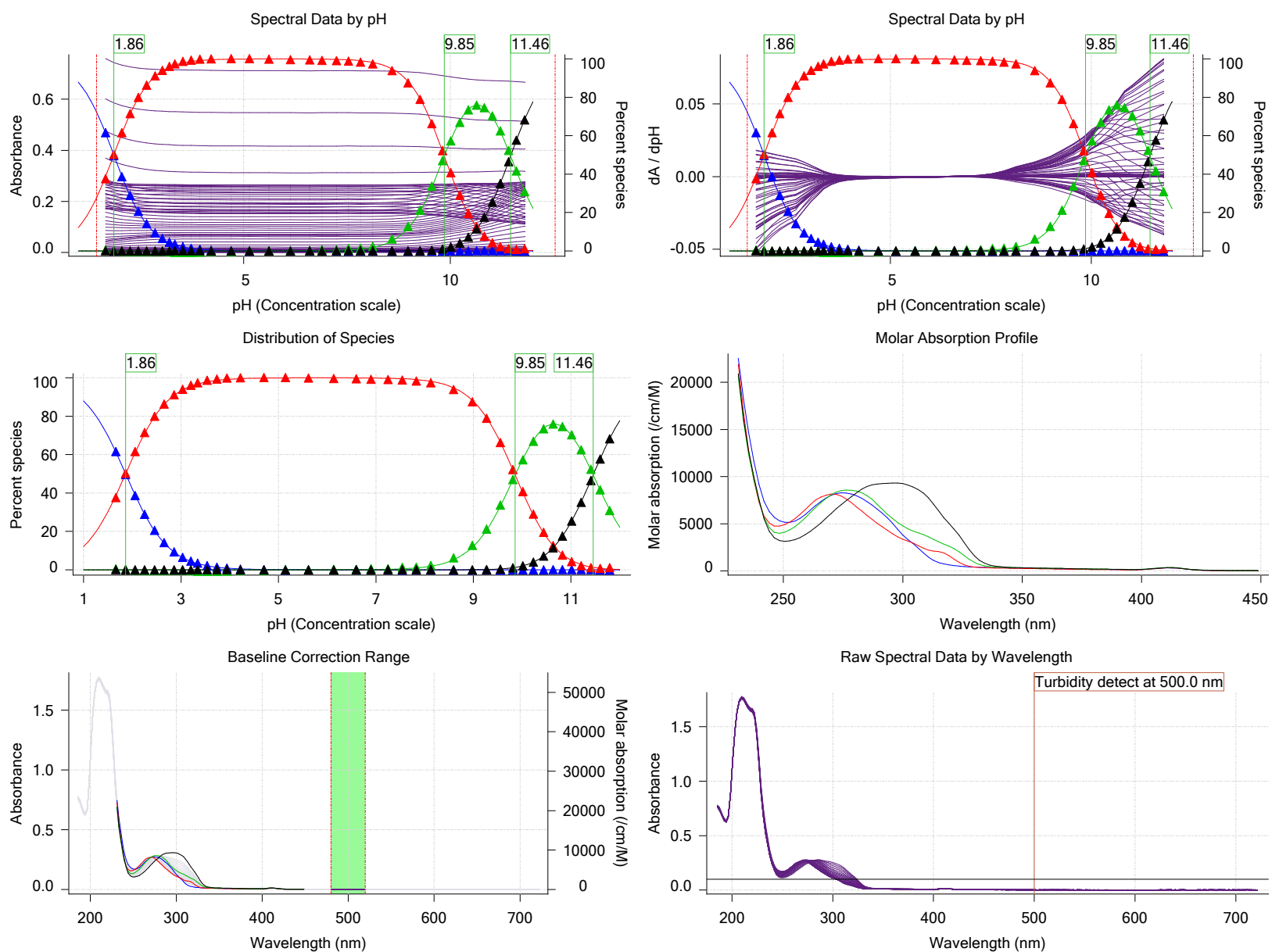

UV-metric psKa Titration 3 of 3 171-28003 Points 74 to 112

## Results

pKa 1 **1.99**  
 pKa 2 **9.84**  
 pKa 3 **11.42**  
 RMSD **0.002 0.007 0.009 0.003**  
 Chi squared **0.0640**  
 PCA calculated number of pKas **4**  
 Average ionic strength **0.176 M**  
 Average temperature **25.0°C**  
 Analyte concentration range **27.8 µM to 26.2 µM**  
 Methanol weight % **39.5 %**  
 Dielectric constant **61.2**

Sample name: **D01**  
 Assay name: **UV-metric psKa**  
 Assay ID: **171-28003**  
 Filename: **C:\Sirius\_T3\Mehtap\20170928\_exp09\_uv\_pKa\171-28003\_D01\_UV-metric psKa.t3r**

Experiment start time: **9/28/2017 4:59:03 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Results (continued)

Water concentration **30.3 M**

Number of pKas source **Manual (3)**  
 Wavelength clipping **230.0 nm to 450.0 nm**  
 pH clipping **1.459 to 12.548**


## Warnings and errors

Errors **None**  
 Warnings **None**

## Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
<b>Assay Medium</b>				
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			


## Graphs

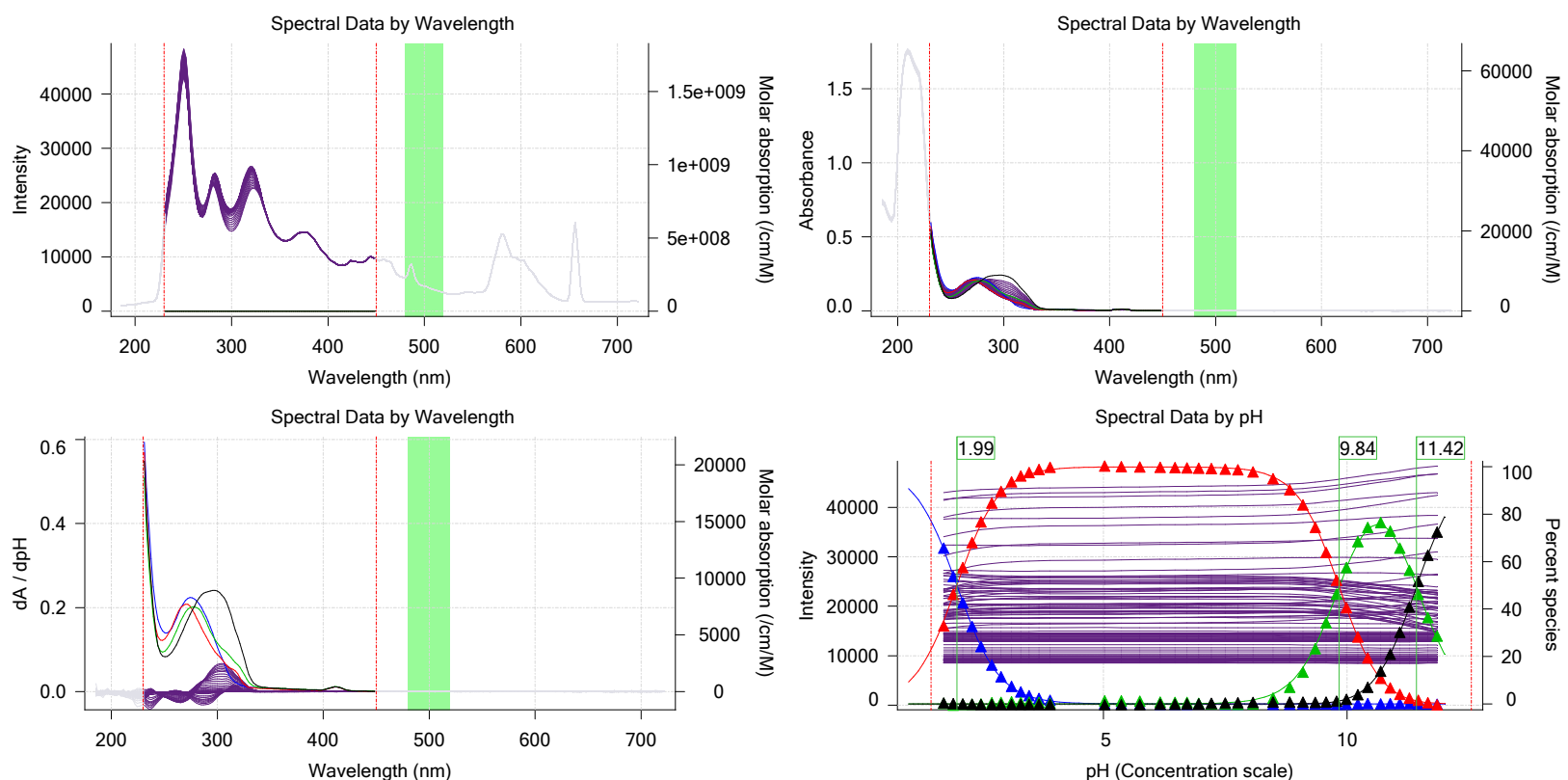

Sample name: **D01**  
 Assay name: **UV-metric psKa**  
 Assay ID: **171-28003**  
 Filename: **C:\Sirius\_T3\Mehtap\20170928\_exp09\_uv\_pKa\171-28003\_D01\_UV-metric psKa.t3r**

Experiment start time: **9/28/2017 4:59:03 PM**  
 Analyst: **Dorothy Leverse**  
 Instrument ID: **T311053**

## Graphs (continued)

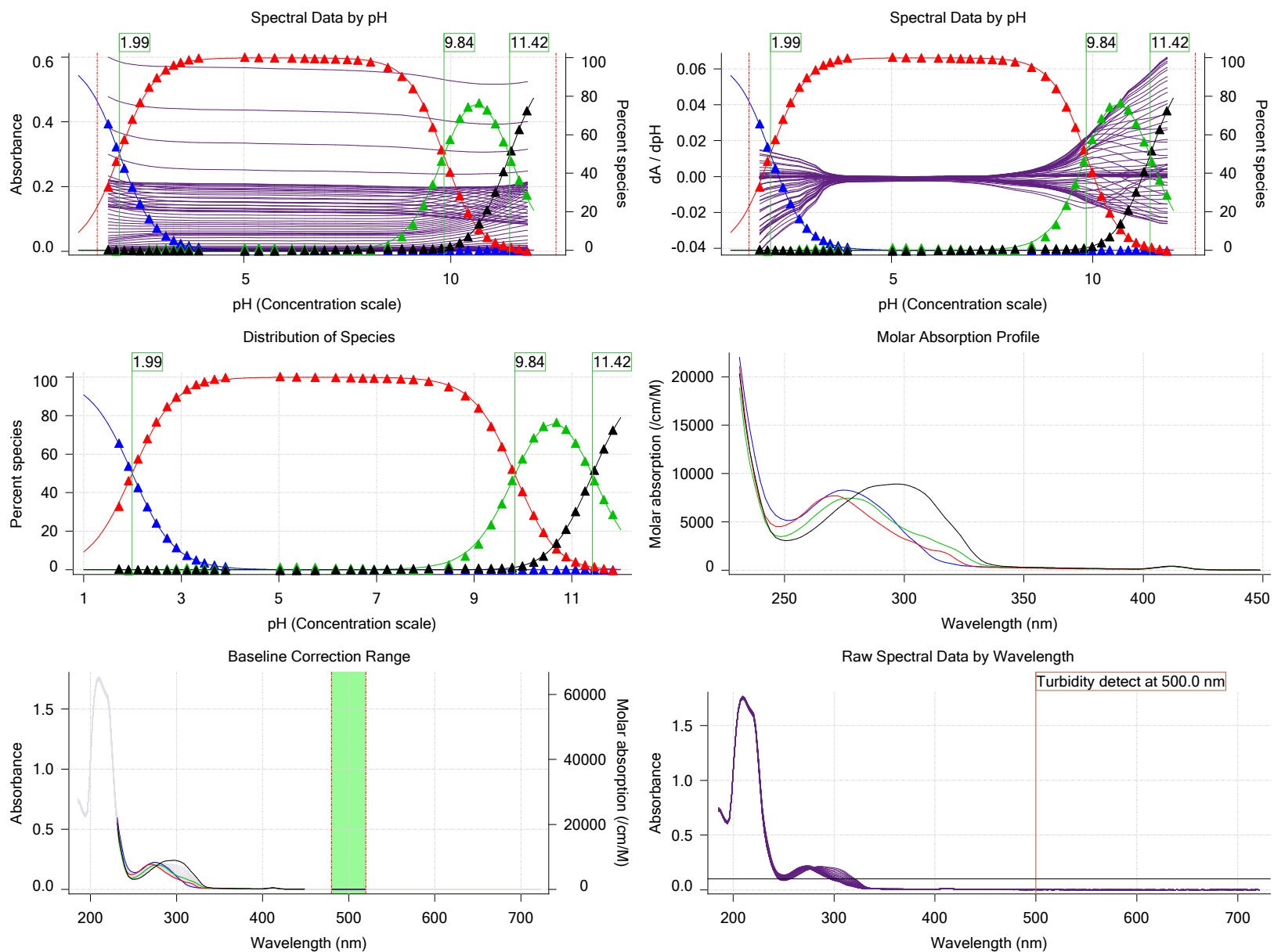

## Assay Model

Settings	Value	Date/Time changed	Imported from
Sample name	D01	9/22/2017 6:27:57 PM	User entered value
Sample by	Volume		Default value
Sample volume	0.0020 mL	9/22/2017 6:27:57 PM	User entered value
Solvent	DMSO		Default value
Sample concentration	0.032100 M	9/22/2017 6:27:57 PM	User entered value
Solubility	Unknown		Default value
Molecular weight	426.44	9/22/2017 6:28:08 PM	User entered value
Individual pKa ionic environments	No		Default value
Number of pKas	4	9/22/2017 6:27:57 PM	User entered value
Sample is a	Ampholyte	9/22/2017 6:27:57 PM	User entered value
pKa 1	2.10	9/22/2017 6:27:57 PM	User entered value
Type	Base	9/22/2017 6:27:57 PM	User entered value
pKa 2	2.99	9/22/2017 6:27:57 PM	User entered value
Type	Base	9/22/2017 6:27:57 PM	User entered value



Sample name: **D01**  
 Assay name: **UV-metric psKa**  
 Assay ID: **171-28003**  
 Filename: **C:\Sirius\_T3\Mehtap\20170928\_exp09\_uv\_pKa\171-28003\_D01\_UV-metric psKa.t3r**

Experiment start time: **9/28/2017 4:59:03 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Assay Model (continued)

Settings	Value	Date/Time changed	Imported from
pKa 3	8.77	9/22/2017 6:27:57 PM	User entered value
Type	Acid	9/22/2017 6:27:57 PM	User entered value
pKa 4	10.90	9/22/2017 6:27:57 PM	User entered value
Type	Acid	9/22/2017 6:27:57 PM	User entered value
logP (XH4 2+)	-10.00		Default value
logp (XH3 +)	-10.00		Default value
logP (neutral XH2)	-10.00	9/22/2017 6:27:57 PM	User entered value
logP (XH -)	-10.00		Default value
logP (X 2-)	-10.00		Default value

## Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
<b>General Settings</b>				
Analyst name	Dorothy Levorse			
Separate reference vial	Yes			
<b>Standard Experiment Settings</b>				
Number of titrations	3			
Minimum pH	2.000			
Maximum pH	12.000			
pH step between points of	0.200			
Minimum titrant addition	0.00002 mL			
Maximum titrant addition	0.10000 mL			
Argon flow rate	100%			
Start titration using	Cautious pH adjust			
<b>Advanced General Settings</b>				
Detect turbidity using	Spectrometer			
Monitor at a wavelength of	500.0 nm			
Absorbance threshold of	0.100			
Collect turbidity sensor data	No			
Stir after titrant addition for	5 seconds			
For titrant addition, stir at	15%			
<b>Titrant Pre-Dose</b>				
Titrant pre-dose	None			
<b>Assay Medium</b>				
Cosolvent in use	Yes			
Cosolvent type	Methanol			
Cosolvent volume	1.35 mL			
Cosolvent added	Automatic			
ISA water volume	0.15 mL			
Water added	Automatic			
After water addition, stir for	5 seconds			
At a speed of	15%			
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			
After medium addition, stir for	5 seconds			
<b>Sample Sonication</b>				
Sonicate	No			
<b>Sample Dissolution</b>				
Perform a dissolution stage	No			
<b>Carbonate purge</b>				
Perform a carbonate purge	No			
<b>Temperature Control</b>				
Wait for temperature	Yes			
Required start temperature	25.0°C			
Acceptable deviation	0.5°C			



Sample name: **D01**  
 Assay name: **UV-metric psKa**  
 Assay ID: **171-28003**  
 Filename: **C:\Sirius\_T3\Mehtap\20170928\_exp09\_uv\_pKa\171-28003\_D01\_UV-metric psKa.t3r**

Experiment start time: **9/28/2017 4:59:03 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
Time to wait	60 seconds			
Stir speed of	15%			
<b>Titration 1</b>				
Titrate from	Low to high pH			
Adjust to start pH	Yes			
After pH adjust stir for	10 seconds			
<b>Titration 2</b>				
Titrate from	Low to high pH			
Additional cosolvent volume	0.00 mL			
Add additional water	0.07 mL			
Additional water added	Automatic			
After pH adjust stir for	10 seconds			
<b>Titration 3</b>				
Titrate from	Low to high pH			
Additional cosolvent volume	0.00 mL			
Add additional water	0.17 mL			
Additional water added	Automatic			
After pH adjust stir for	10 seconds			
<b>Data Point Stability</b>				
Stir during data point collection	Yes			
For point collection, stir at	15%			
Delay before data point collection	0 seconds			
Number of points to average	20 points			
Time interval between points	0.50 seconds			
Required maximum standard deviation	0.00500 dpH/dt			
Stability timeout after	60 seconds			
<b>Experiment cleanup</b>				
Adjust pH to cleanup	To start pH			
And then stir for	60 seconds			
For cleaning, stir at	20%			
Then add water volume	0.25 mL			
And then stir for	30 seconds			

## Calibration Settings

Setting	Value	Date/Time changed	Imported from
Four-Plus alpha	0.105	9/28/2017 4:59:03 PM	C:\Sirius_T3\171-27006_Blank standardisation.t3r
Four-Plus S	1.0031	9/28/2017 4:59:03 PM	C:\Sirius_T3\171-27006_Blank standardisation.t3r
Four-Plus jH	0.7	9/28/2017 4:59:03 PM	C:\Sirius_T3\171-27006_Blank standardisation.t3r
Four-Plus jOH	-0.9	9/28/2017 4:59:03 PM	C:\Sirius_T3\171-27006_Blank standardisation.t3r
Base concentration factor	1.011	9/28/2017 4:59:03 PM	C:\Sirius_T3\KOH17122.t3r
Acid concentration factor	1.007	9/28/2017 4:59:03 PM	C:\Sirius_T3\171-27006_Blank standardisation.t3r

## Instrument Settings

Setting	Value	Batch Id	Install date
Instrument owner	Merck		
Instrument ID	T311053		
Instrument type	T3 Simulator		
Software version	1.1.3.0		
Dispenser module		T3DM1100253	3/31/2009 6:24:52 AM
Dispenser 0	Water		3/31/2009 6:25:05 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Water (0.15 M KCl)	8-18-17	9/26/2017 9:05:04 AM
Dispenser 2	Acid		3/31/2009 6:25:11 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		

Sample name: **D01**  
 Assay name: **UV-metric psKa**  
 Assay ID: **171-28003**  
 Filename: **C:\Sirius\_T3\Mehtap\20170928\_exp09\_uv\_pKa\171-28003\_D01\_UV-metric psKa.t3r**

Experiment start time: **9/28/2017 4:59:03 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Instrument Settings (continued)

Setting	Value	Batch Id	Install date
Titrant	Acid (0.5 M HCl)	166940	9/8/2017 9:21:27 AM
Dispenser 1	Base		3/31/2009 6:25:21 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Base (0.5 M KOH)	9-22-17	9/22/2017 4:02:42 PM
Dispenser 5	Cosolvent		3/31/2009 6:26:24 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Distribution valve 5	Distribution Valve		3/31/2009 6:28:19 AM
Firmware version	1.1.3		
Port A	Methanol (80%, 0.15 M KCl)	8-15-17	9/20/2017 4:38:16 PM
Port B	Cyclohexane		9/19/2017 2:15:02 PM
Dispenser 3	Buffer		8/3/2010 6:05:16 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Phosphate Buffer		9/12/2017 12:32:29 PM
Dispenser 6	Octanol		10/22/2010 11:52:43 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Octanol	9-14-17	9/14/2017 10:30:38 AM
Titration		T3TM1100153	3/31/2009 6:24:17 AM
Horizontal axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Probe I/O firmware version	1.1.1		
Electrode	T3 Electrode	T3E0769	8/15/2017 10:21:54 AM
E0 calibration	-2.72 mV		9/28/2017 4:59:27 PM
Filling solution	3M KCl	KCL095	9/28/2017 1:58:38 PM
Liquids			
Wash 1	50% IPA:50% Water		9/28/2017 1:57:12 PM
Wash 2	0.5% Triton X-100 in H2O		9/28/2017 1:57:15 PM
Buffer position 1	pH7 Wash		9/28/2017 1:57:18 PM
Buffer position 2	pH 7		9/28/2017 1:57:25 PM
Storage position			9/28/2017 1:57:49 PM
Wash water	9.6e+003 mL	9-27-17	9/27/2017 4:24:06 PM
Waste	4e+002 mL		9/27/2017 4:24:14 PM
Temperature controller			8/5/2010 7:35:13 AM
Turbidity detector			3/31/2009 6:24:45 AM
Spectrometer		072390	11/23/2010 12:22:28 PM
Dip probe		11086	
Wavelength coefficient A0	185.563		
Wavelength coefficient A1	2.17439		
Wavelength coefficient A2	-0.000285622		
Total lamp lit time	269:59:45		11/23/2010 12:22:28 PM
Calibrated on	9/26/2017 9:22:07 AM		
Integration time	11		
Scans averaged	10		
Autoloader		T3AL1100237	11/10/2015 10:34:13 AM
Left-right axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Front-back axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Configuration			
Alternate titration position	Titration position		
Alternate reference position	Reference position		
Maximum standard vial volume	3.50 mL		
Maximum alternate vial volume	25.00 mL		
Automatic action idle period	5 minute(s)		



## Assay Settings

Sample name: **D01**  
Assay name: **UV-metric psKa**  
Assay ID: **17I-28003**  
Filename: **C:\Sirius\_T3\Mehtap\20170928\_exp09\_uv\_pKa\17I-28003\_D01\_UV-metric psKa.t3r**

Experiment start time: **9/28/2017 4:59:03 PM**  
Analyst: **Dorothy Levorse**  
Instrument ID: **T311053**

### Instrument Settings (continued)

Setting	Value	Batch Id	Install date
Titrant tube volume	1.3 mL		
Syringe flush count	3.50		
Flowing wash pump volume	20.0 mL		
Flowing wash stir duration	5 s		
Flowing wash stir speed	30%		
Solvent wash stir duration	5 s		
Solvent wash stir speed	30%		
Surfactant wash stir duration	5 s		
Surfactant wash stir speed	30%		
E0 calibration minimum number of points	10		
E0 calibration maximum standard deviation	0.01500		
E0 calibration timeout period	60 s		
E0 calibration stir duration	5 s		
E0 calibration preparation stir speed	30%		
E0 calibration buffer wash stir duration	5 s		
E0 calibration buffer wash stir speed	30%		
E0 calibration reading stir speed	0%		
Spectrometer calibration stir duration	5 s		
Spectrometer calibration stir speed	30%		
Spectrometer calibration wash pump volume	20.0 mL		
Spectrometer calibration wash stir duration	5 s		
Spectrometer calibration wash stir speed	30%		
Overhead dispense height	10000		

### Refinement Settings

Setting	Value	Default value
Turbidity detection method	Spectrometer	Spectrometer
Turbidity wavelength to assess	500.0 nm	500.0 nm
Turbidity maximum absorbance	0.100	0.100
Turbidity probe threshold	50.00	50.00
Exclude turbid points	Yes	Yes
Low intensity warning threshold	100	100
Minimum absorbance change threshold	0.100	0.100
Eigenvector autocorrelation threshold	0.80	0.80
Maximum RMSD severe warning	0.250	0.250
Maximum RMSD warning	0.050	0.050

### Tray Information

Title  
Location A5