

Sample name: **M02**
 Assay name: **UV-metric psKa**
 Assay ID: **17I-22023**
 Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22023_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 7:46:07 PM**
 Analyst: **Dorothy Levorse**
 Instrument ID: **T311053**




Yasuda-Shedlovsky result

Extrapolation type	pKa 0%	SD	Intercept	Slope	R ²	Ionic strength	Temperature
Yasuda-Shedlovsky	5.04	±0.08	9.32	-199.4852	0.9883	0.165 M	25.0°C

Component assay results

Titration	Methanol weight%	Direction	Result type	Dielectric constant	[H ₂ O]	Ionic strength	Temperature	psKa 1
17I-22023 Points 4 to 42	49.51 %	Up	UV-metric pKa	56.6	24.7 M	0.157 M	25.0°C	✓ 4.42
17I-22023 Points 44 to 81	40.01 %	Up	UV-metric pKa	61.0	30.0 M	0.165 M	25.0°C	✓ 4.54
17I-22023 Points 83 to 127	30.14 %	Up	UV-metric pKa	65.4	35.8 M	0.172 M	25.0°C	✓ 4.74

Graphs

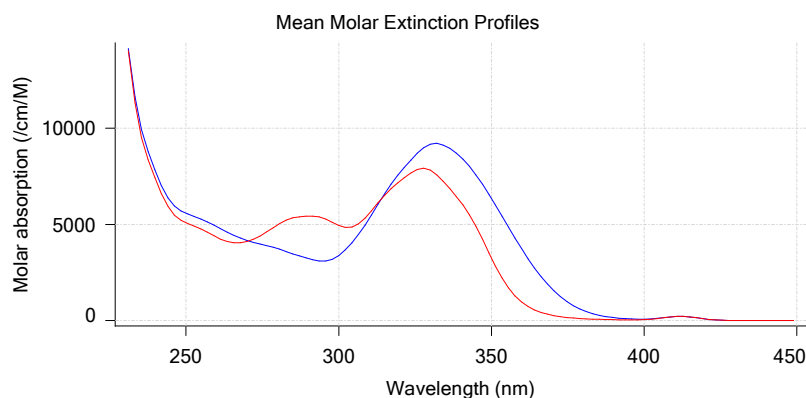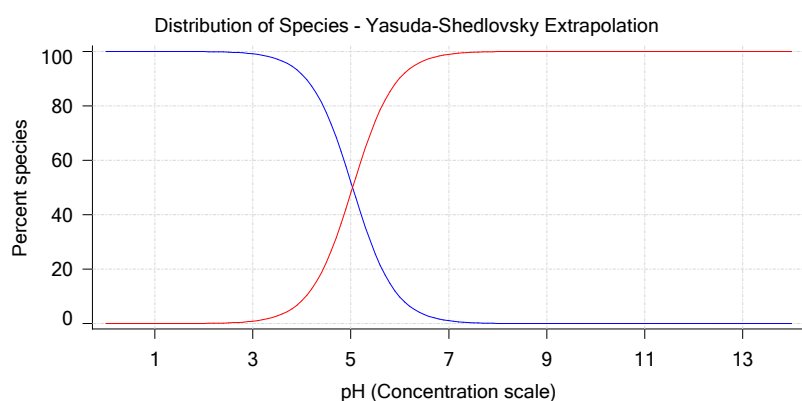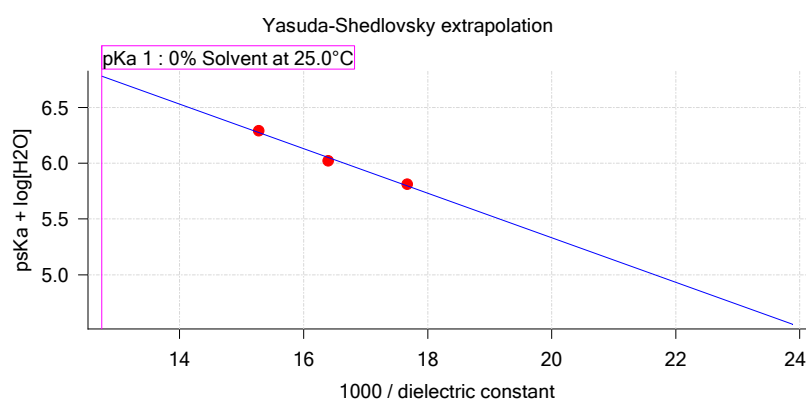

UV-metric psKa Titration 1 of 3 17I-22023 Points 4 to 42

Results

pKa 1	4.42
RMSD	0.022 0.022
Chi squared	0.0893
PCA calculated number of pKas	4
Average ionic strength	0.157 M
Average temperature	25.0°C
Analyte concentration range	58.2 µM to 54.8 µM
Methanol weight %	49.5 %
Dielectric constant	56.6
Water concentration	24.7 M
Number of pKas source	Predicted
Wavelength clipping	230.0 nm to 450.0 nm

Sample name: **M02**
 Assay name: **UV-metric psKa**
 Assay ID: **17I-22023**
 Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22023_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 7:46:07 PM**
 Analyst: **Dorothy Levorse**
 Instrument ID: **T311053**

Results (continued)

pH clipping 1.477 to 12.527


Warnings and errors

Errors None
 Warnings PCA calculation disagrees with predicted number of pKas

Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
Assay Medium				
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			


Graphs

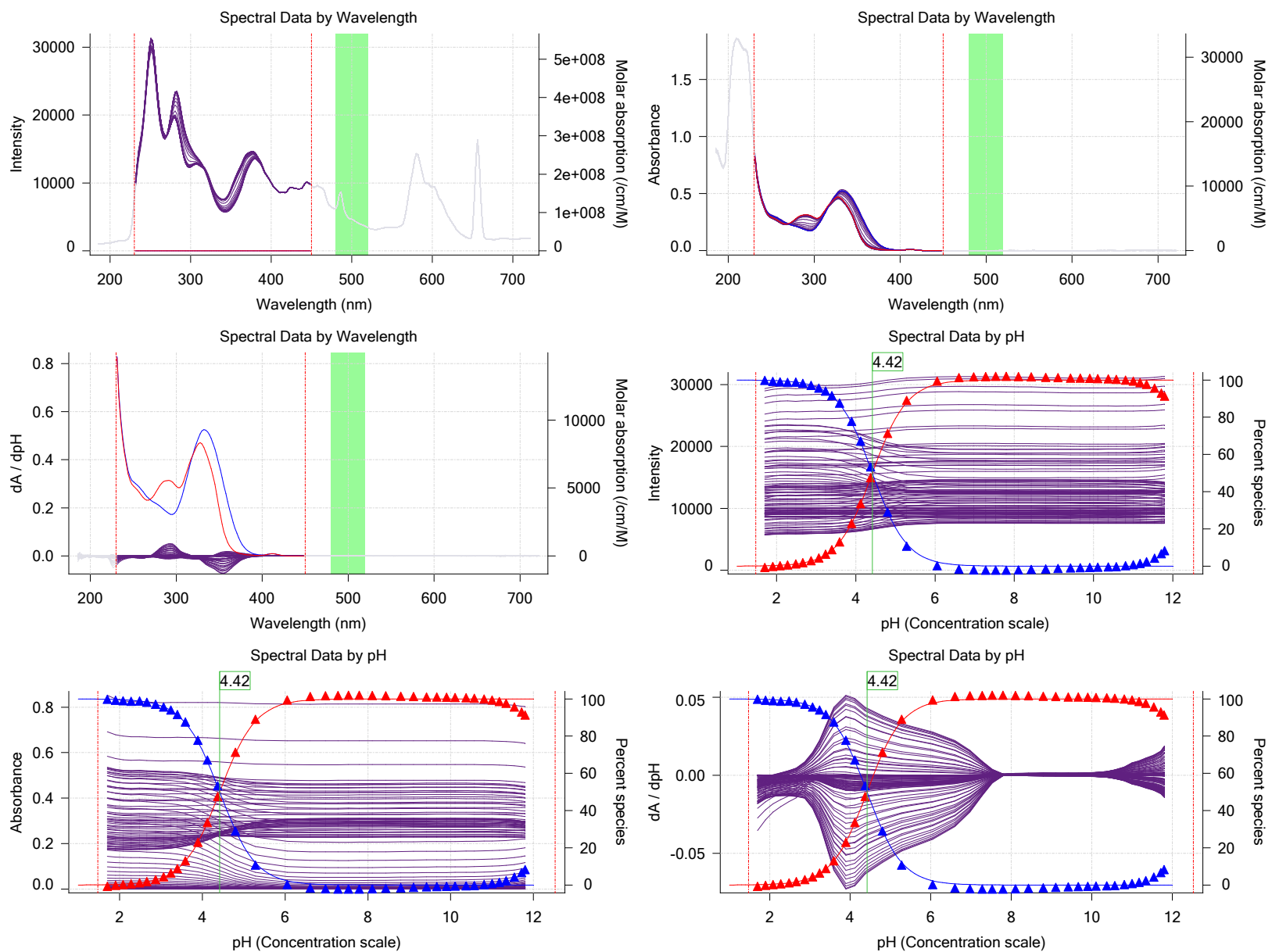

Sample name: **M02**
Assay name: **UV-metric psKa**
Assay ID: **171-22023**
Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\171-22023_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 7:46:07 PM**
Analyst: **Dorothy Levorse**
Instrument ID: **T311053**

Graphs (continued)

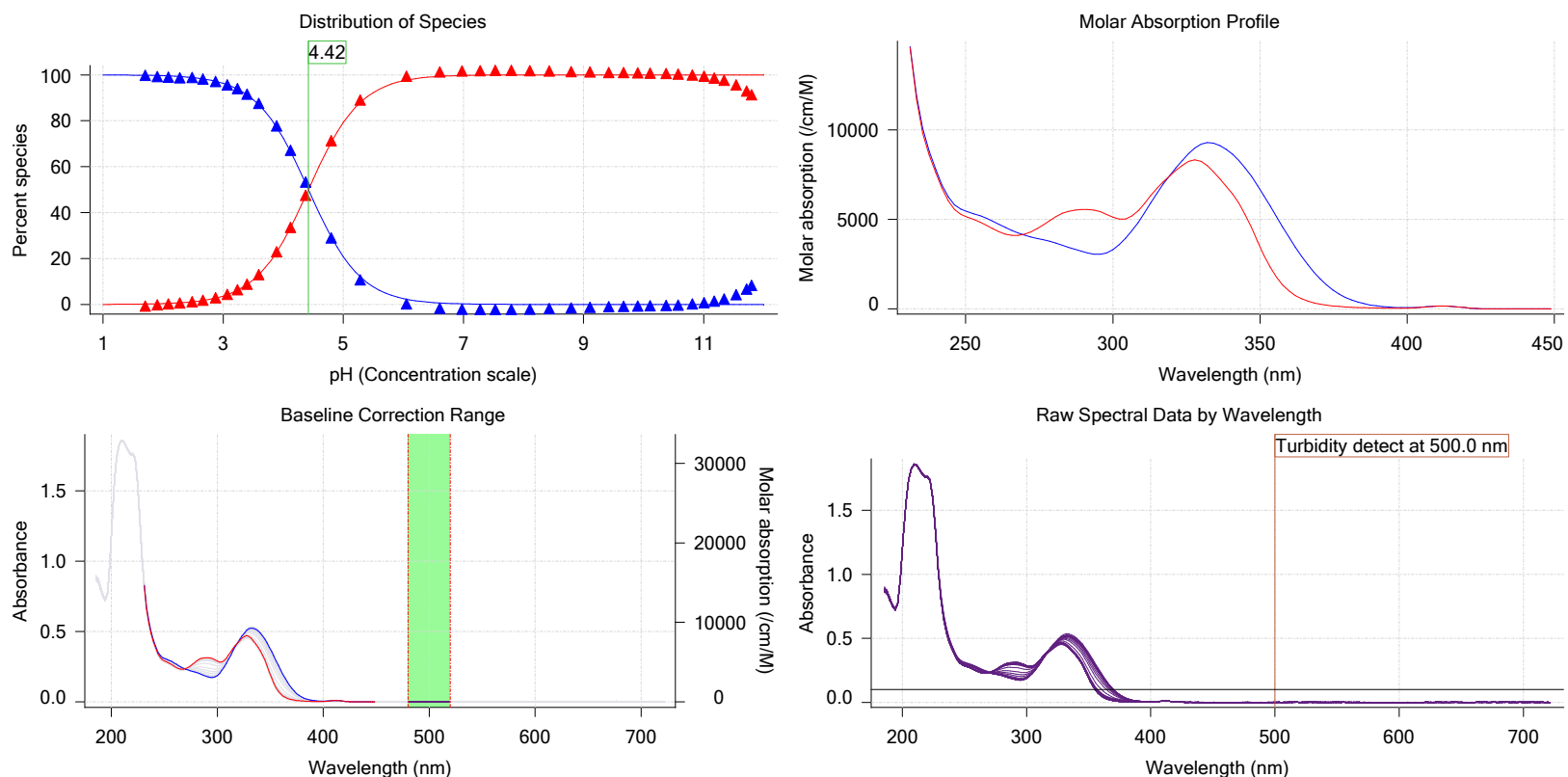

UV-metric psKa Titration 2 of 3 171-22023 Points 44 to 81

Results

pKa 1	4.54
RMSD	0.025 0.026
Chi squared	0.0964
PCA calculated number of pKas	3
Average ionic strength	0.165 M
Average temperature	25.0°C
Analyte concentration range	47.8 µM to 45.2 µM
Methanol weight %	40.0 %
Dielectric constant	61.0
Water concentration	30.0 M
Number of pKas source	Predicted
Wavelength clipping	230.0 nm to 450.0 nm
pH clipping	1.486 to 12.532

Warnings and errors

Errors: None
Warnings: PCA calculation disagrees with predicted number of pKas

Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			

Assay Medium


Sample name: **M02**
 Assay name: **UV-metric psKa**
 Assay ID: **17I-22023**
 Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22023_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 7:46:07 PM**
 Analyst: **Dorothy Levorse**
 Instrument ID: **T311053**

Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			


Graphs

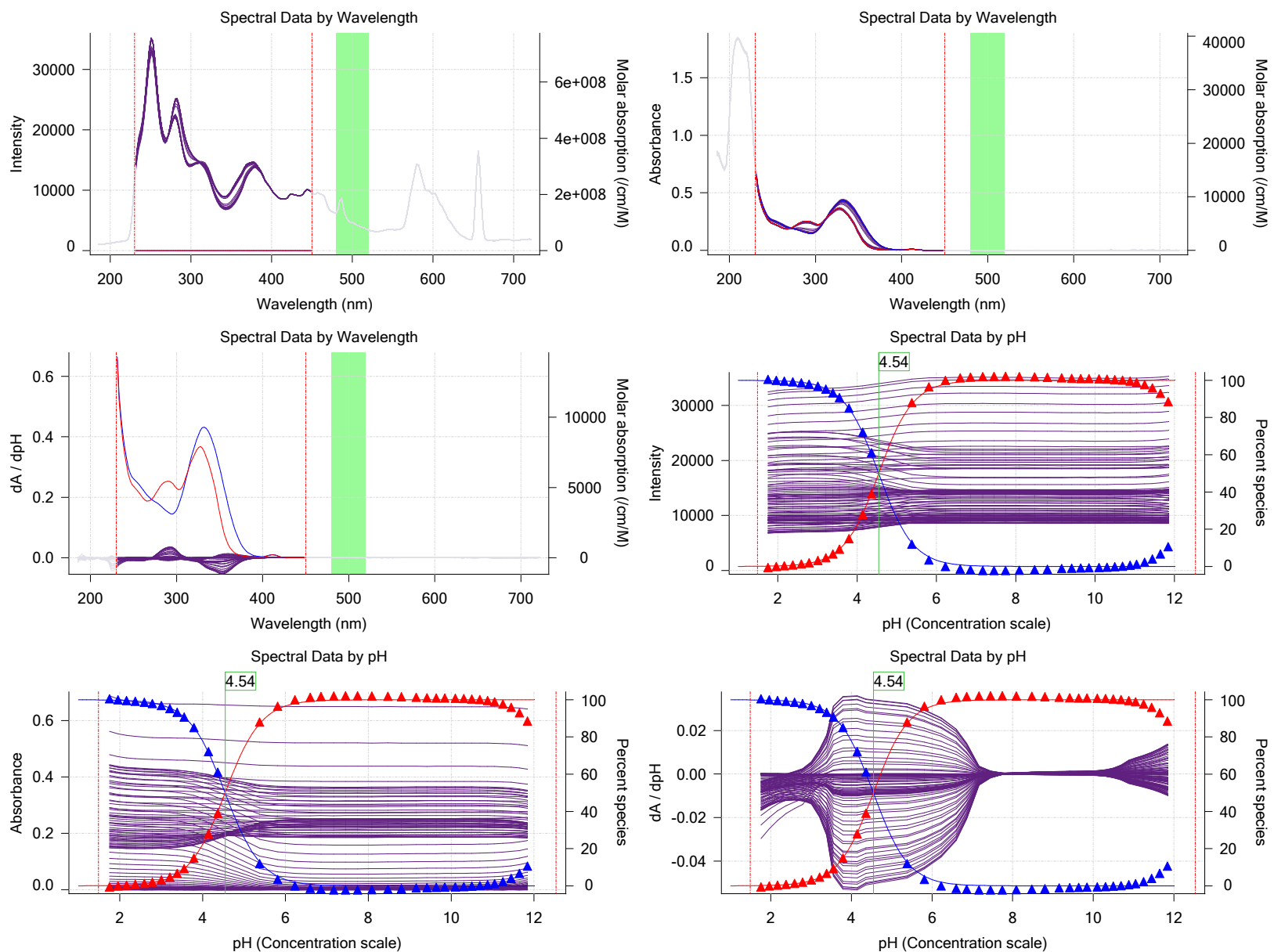

Sample name: **M02**
Assay name: **UV-metric psKa**
Assay ID: **171-22023**
Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\171-22023_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 7:46:07 PM**
Analyst: **Dorothy Leverse**
Instrument ID: **T311053**

Graphs (continued)

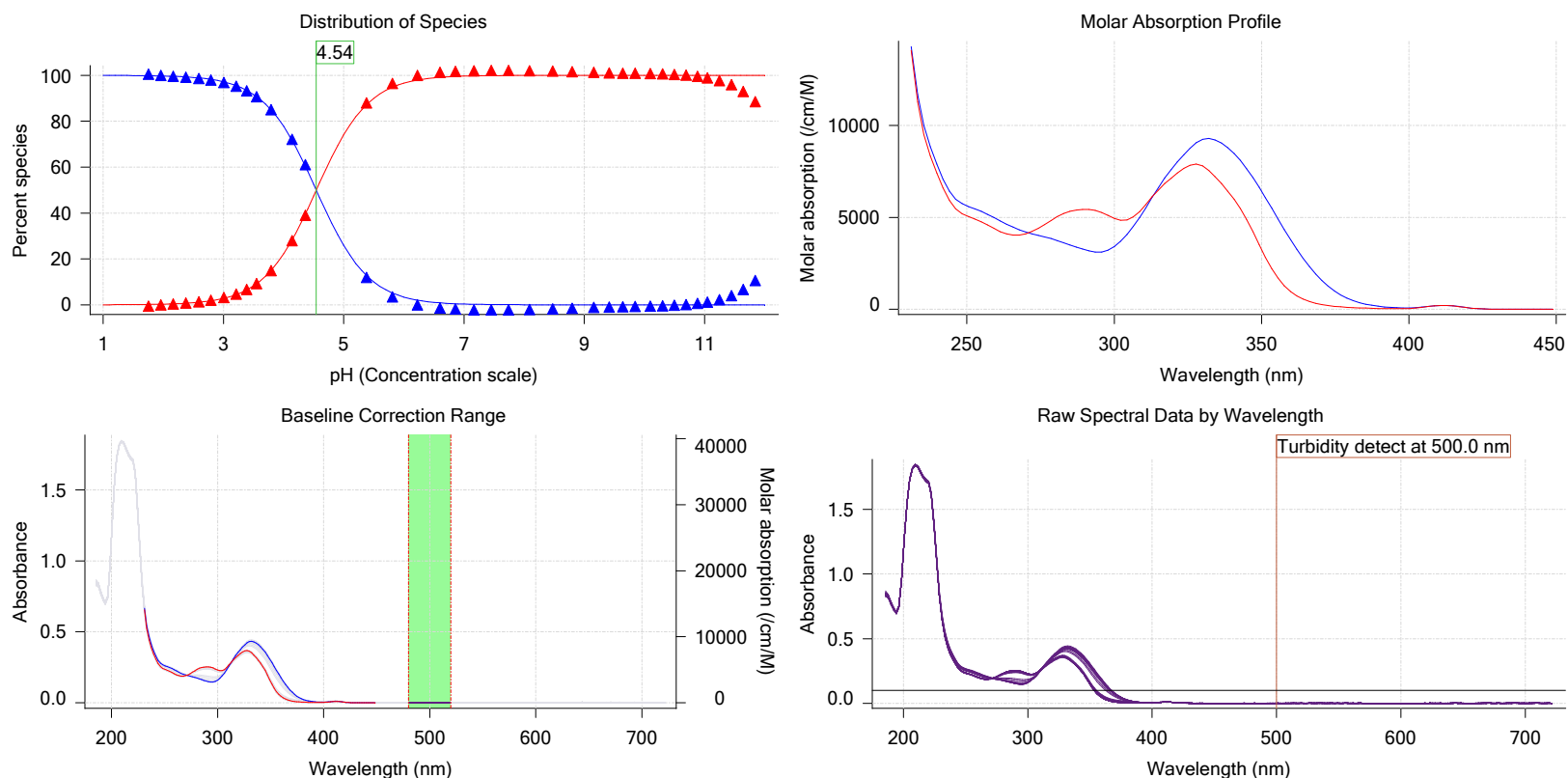

UV-metric psKa Titration 3 of 3 171-22023 Points 83 to 127

Results

pKa 1	4.74
RMSD	0.034 0.036
Chi squared	0.1293
PCA calculated number of pKas	4
Average ionic strength	0.172 M
Average temperature	25.0°C
Analyte concentration range	36.8 µM to 34.8 µM
Methanol weight %	30.1 %
Dielectric constant	65.4
Water concentration	35.8 M
Number of pKas source	Predicted
Wavelength clipping	230.0 nm to 450.0 nm
pH clipping	1.495 to 12.537

Warnings and errors

Errors: None
Warnings: PCA calculation disagrees with predicted number of pKas

Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			

Assay Medium


Sample name: **M02**
 Assay name: **UV-metric psKa**
 Assay ID: **17I-22023**
 Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22023_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 7:46:07 PM**
 Analyst: **Dorothy Levorse**
 Instrument ID: **T311053**

Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			


Graphs

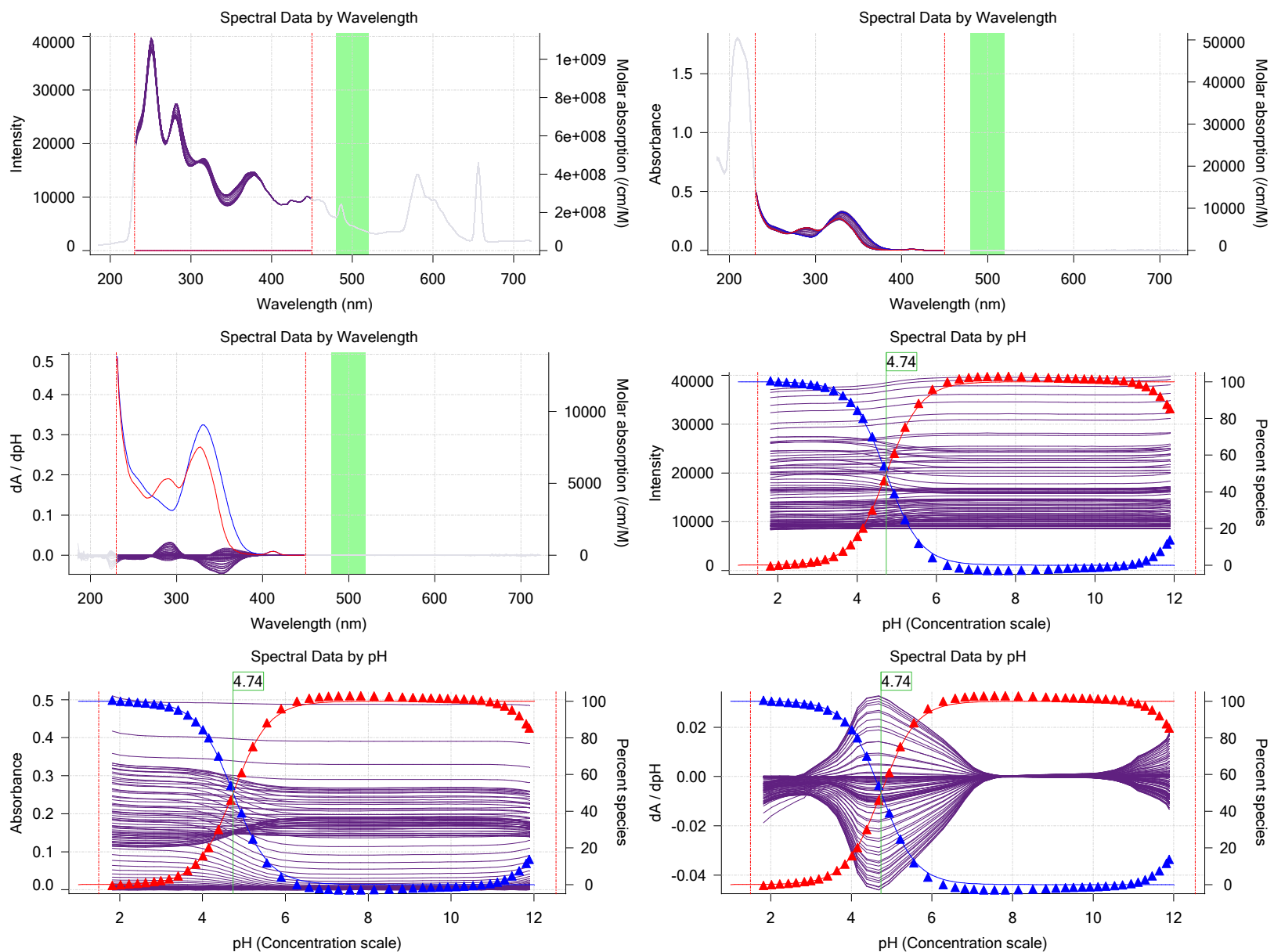

Sample name: **M02**
Assay name: **UV-metric psKa**
Assay ID: **171-22023**
Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\171-22023_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 7:46:07 PM**
Analyst: **Dorothy Leverse**
Instrument ID: **T311053**

Graphs (continued)

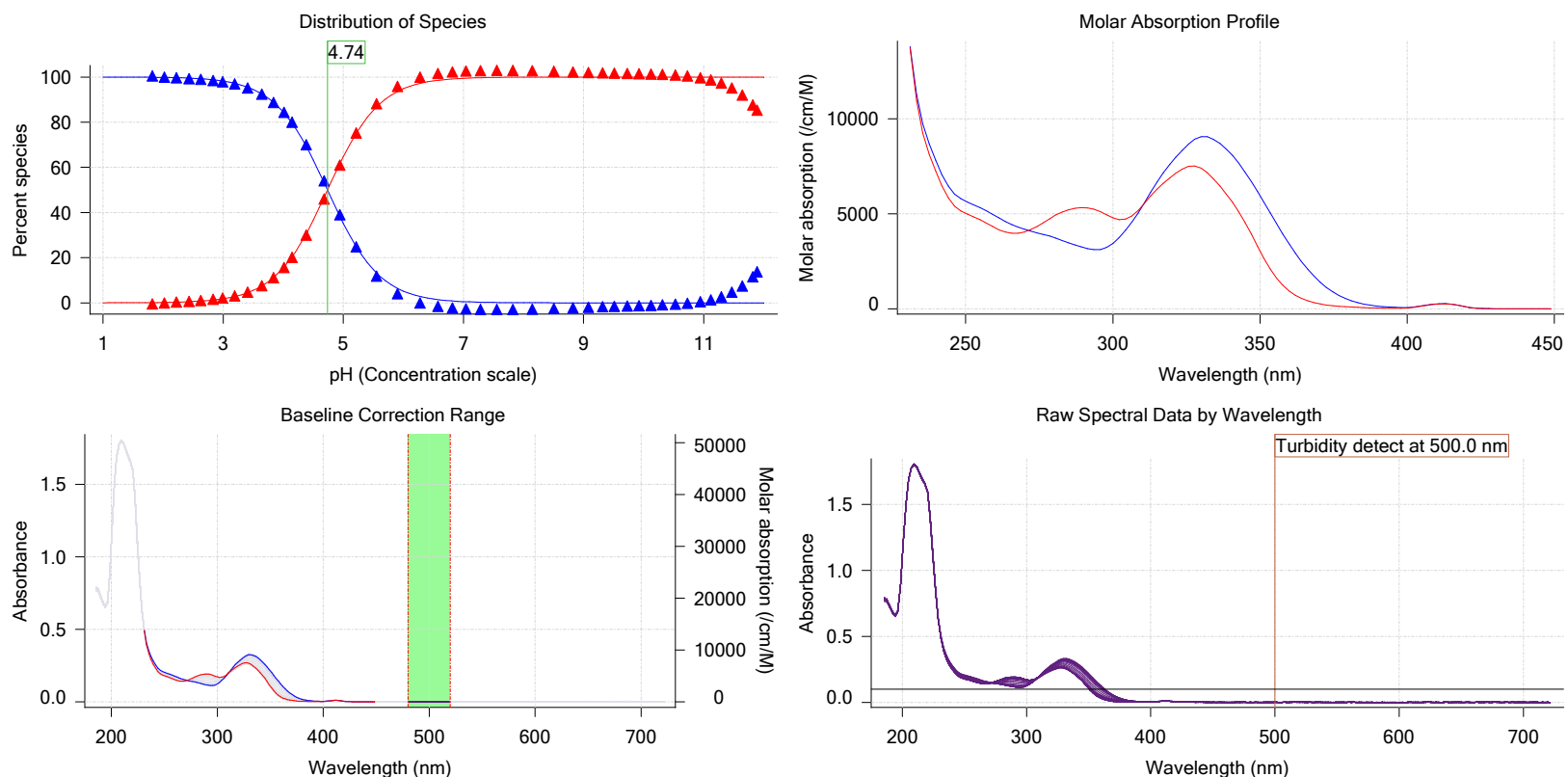

Assay Model

Settings	Value	Date/Time changed	Imported from
Sample name	M02	9/22/2017 10:27:56 AM	User entered value
Sample by	Volume		Default value
Sample volume	0.0020 mL	9/22/2017 10:27:56 AM	User entered value
Solvent	DMSO		Default value
Sample concentration	0.046400 M	9/22/2017 10:27:56 AM	User entered value
Solubility	Unknown		Default value
Molecular weight	289.26	9/22/2017 10:28:05 AM	User entered value
Individual pKa ionic environments	No		Default value
Number of pKas	1	9/22/2017 10:27:56 AM	User entered value
Sample is a	Base	9/22/2017 10:27:56 AM	User entered value
pKa 1	5.60	9/22/2017 10:27:56 AM	User entered value
logp (XH +)	-10.00		Default value
logP (neutral X)	-10.00	9/22/2017 10:27:56 AM	User entered value

Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
General Settings				
Analyst name	Dorothy Leverse			
Separate reference vial	Yes			
Standard Experiment Settings				
Number of titrations	3			
Minimum pH	2.000			
Maximum pH	12.000			
pH step between points of	0.200			
Minimum titrant addition	0.00002 mL			
Maximum titrant addition	0.10000 mL			
Argon flow rate	100%			

Sample name: **M02**
 Assay name: **UV-metric psKa**
 Assay ID: **17I-22023**
 Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22023_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 7:46:07 PM**
 Analyst: **Dorothy Levorse**
 Instrument ID: **T311053**

Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
Start titration using	Cautious pH adjust			
Advanced General Settings				
Detect turbidity using	Spectrometer			
Monitor at a wavelength of	500.0 nm			
Absorbance threshold of	0.100			
Collect turbidity sensor data	No			
Stir after titrant addition for	5 seconds			
For titrant addition, stir at	15%			
Titration Pre-Dose				
Titration pre-dose	None			
Assay Medium				
Cosolvent in use	Yes			
Cosolvent type	Methanol			
Cosolvent volume	1.15 mL			
Cosolvent added	Automatic			
ISA water volume	0.35 mL			
Water added	Automatic			
After water addition, stir for	5 seconds			
At a speed of	15%			
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			
After medium addition, stir for	5 seconds			
Sample Sonication				
Sonicate	No			
Sample Dissolution				
Perform a dissolution stage	No			
Carbonate purge				
Perform a carbonate purge	No			
Temperature Control				
Wait for temperature	Yes			
Required start temperature	25.0°C			
Acceptable deviation	0.5°C			
Time to wait	60 seconds			
Stir speed of	15%			
Titration 1				
Titrate from	Low to high pH			
Adjust to start pH	Yes			
After pH adjust stir for	10 seconds			
Titration 2				
Titrate from	Low to high pH			
Additional cosolvent volume	0.00 mL			
Add additional water	0.15 mL			
Additional water added	Automatic			
After pH adjust stir for	10 seconds			
Titration 3				
Titrate from	Low to high pH			
Additional cosolvent volume	0.00 mL			
Add additional water	0.34 mL			
Additional water added	Automatic			
After pH adjust stir for	10 seconds			
Data Point Stability				
Stir during data point collection	Yes			
For point collection, stir at	15%			
Delay before data point collection	0 seconds			
Number of points to average	20 points			
Time interval between points	0.50 seconds			

Sample name: **M02**
 Assay name: **UV-metric psKa**
 Assay ID: **17I-22023**
 Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22023_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 7:46:07 PM**
 Analyst: **Dorothy Levorse**
 Instrument ID: **T311053**

Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
Required maximum standard deviation	0.00500 dpH/dt			
Stability timeout after	60 seconds			
Experiment cleanup				
Adjust pH to cleanup	To start pH			
And then stir for	60 seconds			
For cleaning, stir at	20%			
Then add water volume	0.25 mL			
And then stir for	30 seconds			

Calibration Settings

Setting	Value	Date/Time changed	Imported from
Four-Plus alpha	0.105	9/22/2017 7:46:07 PM	C:\Sirius_T3\17I-22019_Blank standardisation.t3r
Four-Plus S	1.0019	9/22/2017 7:46:07 PM	C:\Sirius_T3\17I-22019_Blank standardisation.t3r
Four-Plus jH	0.3	9/22/2017 7:46:07 PM	C:\Sirius_T3\17I-22019_Blank standardisation.t3r
Four-Plus jOH	-1.0	9/22/2017 7:46:07 PM	C:\Sirius_T3\17I-22019_Blank standardisation.t3r
Base concentration factor	1.000	9/22/2017 7:46:07 PM	C:\Sirius_T3\17I-22016_KHP_Base standardisation using KHP.t3r
Acid concentration factor	0.995	9/22/2017 7:46:07 PM	C:\Sirius_T3\17I-22019_Blank standardisation.t3r

Instrument Settings

Setting	Value	Batch Id	Install date
Instrument owner	Merck		
Instrument ID	T311053		
Instrument type	T3 Simulator		
Software version	1.1.3.0		
Dispenser module		T3DM1100253	3/31/2009 6:24:52 AM
Dispenser 0	Water		3/31/2009 6:25:05 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Water (0.15 M KCl)	8-18-17	9/18/2017 9:13:04 AM
Dispenser 2	Acid		3/31/2009 6:25:11 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Acid (0.5 M HCl)	166940	9/8/2017 9:21:27 AM
Dispenser 1	Base		3/31/2009 6:25:21 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Base (0.5 M KOH)	9-22-17	9/22/2017 4:02:42 PM
Dispenser 5	Cosolvent		3/31/2009 6:26:24 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Distribution valve 5	Distribution Valve		3/31/2009 6:28:19 AM
Firmware version	1.1.3		
Port A	Methanol (80%, 0.15 M KCl)	8-15-17	9/20/2017 4:38:16 PM
Port B	Cyclohexane		9/19/2017 2:15:02 PM
Dispenser 3	Buffer		8/3/2010 6:05:16 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Phosphate Buffer		9/12/2017 12:32:29 PM
Dispenser 6	Octanol		10/22/2010 11:52:43 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Octanol	9-14-17	9/14/2017 10:30:38 AM
Titrator		T3TM1100153	3/31/2009 6:24:17 AM
Horizontal axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		



Assay Settings

Sample name: **M02** Experiment start time: **9/22/2017 7:46:07 PM**
Assay name: **UV-metric psKa** Analyst: **Dorothy Levorse**
Assay ID: **17I-22023** Instrument ID: **T311053**
Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22023_M02_UV-metric psKa.t3r**

Instrument Settings (continued)

Setting	Value	Batch Id	Install date
Probe I/O firmware version	1.1.1		
Electrode	T3 Electrode	T3E0769	8/15/2017 10:21:54 AM
E0 calibration	-8.13 mV		9/22/2017 7:46:31 PM
Filling solution	3M KCl	KCL095	9/18/2017 9:17:15 AM
Liquids			
Wash 1	50% IPA:50% Water		9/22/2017 10:05:34 AM
Wash 2	0.5% Triton X-100 in H2O		9/22/2017 10:05:36 AM
Buffer position 1	pH7 Wash		9/22/2017 10:05:39 AM
Buffer position 2	pH 7		9/22/2017 10:05:42 AM
Storage position			9/22/2017 10:06:32 AM
Wash water	6.1e+002 mL	9-18-17	9/18/2017 8:54:32 AM
Waste	9.5e+003 mL		9/18/2017 8:54:39 AM
Temperature controller			8/5/2010 7:35:13 AM
Turbidity detector			3/31/2009 6:24:45 AM
Spectrometer		072390	11/23/2010 12:22:28 PM
Dip probe		11086	
Wavelength coefficient A0	185.563		
Wavelength coefficient A1	2.17439		
Wavelength coefficient A2	-0.000285622		
Total lamp lit time	207:39:18		11/23/2010 12:22:28 PM
Calibrated on	9/18/2017 9:35:14 AM		
Integration time	11		
Scans averaged	10		
Autoloader		T3AL1100237	11/10/2015 10:34:13 AM
Left-right axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Front-back axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Configuration			
Alternate titration position	Titration position		
Alternate reference position	Reference position		
Maximum standard vial volume	3.50 mL		
Maximum alternate vial volume	25.00 mL		
Automatic action idle period	5 minute(s)		
Titrant tube volume	1.3 mL		
Syringe flush count	3.50		
Flowing wash pump volume	20.0 mL		
Flowing wash stir duration	5 s		
Flowing wash stir speed	30%		
Solvent wash stir duration	5 s		
Solvent wash stir speed	30%		
Surfactant wash stir duration	5 s		
Surfactant wash stir speed	30%		
E0 calibration minimum number of points	10		
E0 calibration maximum standard deviation	0.01500		
E0 calibration timeout period	60 s		
E0 calibration stir duration	5 s		
E0 calibration preparation stir speed	30%		
E0 calibration buffer wash stir duration	5 s		
E0 calibration buffer wash stir speed	30%		
E0 calibration reading stir speed	0%		
Spectrometer calibration stir duration	5 s		
Spectrometer calibration stir speed	30%		
Spectrometer calibration wash pump volume	20.0 mL		
Spectrometer calibration wash stir duration	5 s		
Spectrometer calibration wash stir speed	30%		
Overhead dispense height	10000		



Assay Settings

Sample name: **M02**
Assay name: **UV-metric psKa**
Assay ID: **17I-22023**
Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22023_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 7:46:07 PM**
Analyst: **Dorothy Levorse**
Instrument ID: **T311053**

Refinement Settings

Setting	Value	Default value
Turbidity detection method	Spectrometer	Spectrometer
Turbidity wavelength to assess	500.0 nm	500.0 nm
Turbidity maximum absorbance	0.100	0.100
Turbidity probe threshold	50.00	50.00
Exclude turbid points	Yes	Yes
Low intensity warning threshold	100	100
Minimum absorbance change threshold	0.100	0.100
Eigenvector autocorrelation threshold	0.80	0.80
Maximum RMSD severe warning	0.250	0.250
Maximum RMSD warning	0.050	0.050

Tray Information

Title
Location B3