



Sample name: **M11**
Assay name: **UV-metric pKa**
Assay ID: **17I-16007**
Filename: **C:\Sirius_T3\17I-16007_M11_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 2:03:46 AM**
Analyst: **Dorothy Leverse**
Instrument ID: **T311053**

Results

pKa 1 **3.89**
RMSD **0.001 0.002**
Chi squared **0.0091**
PCA calculated number of pKas **2**
Average ionic strength **0.158 M**
Average temperature **24.9°C**
Analyte concentration range **80.0 µM to 72.5 µM**

Number of pKas source **Predicted**
Wavelength clipping **230.0 nm to 450.0 nm**
pH clipping **1.271 to 12.704**


Warnings and errors

Errors None
Warnings PCA calculation disagrees with predicted number of pKas

Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
Assay Medium				
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			









Graphs

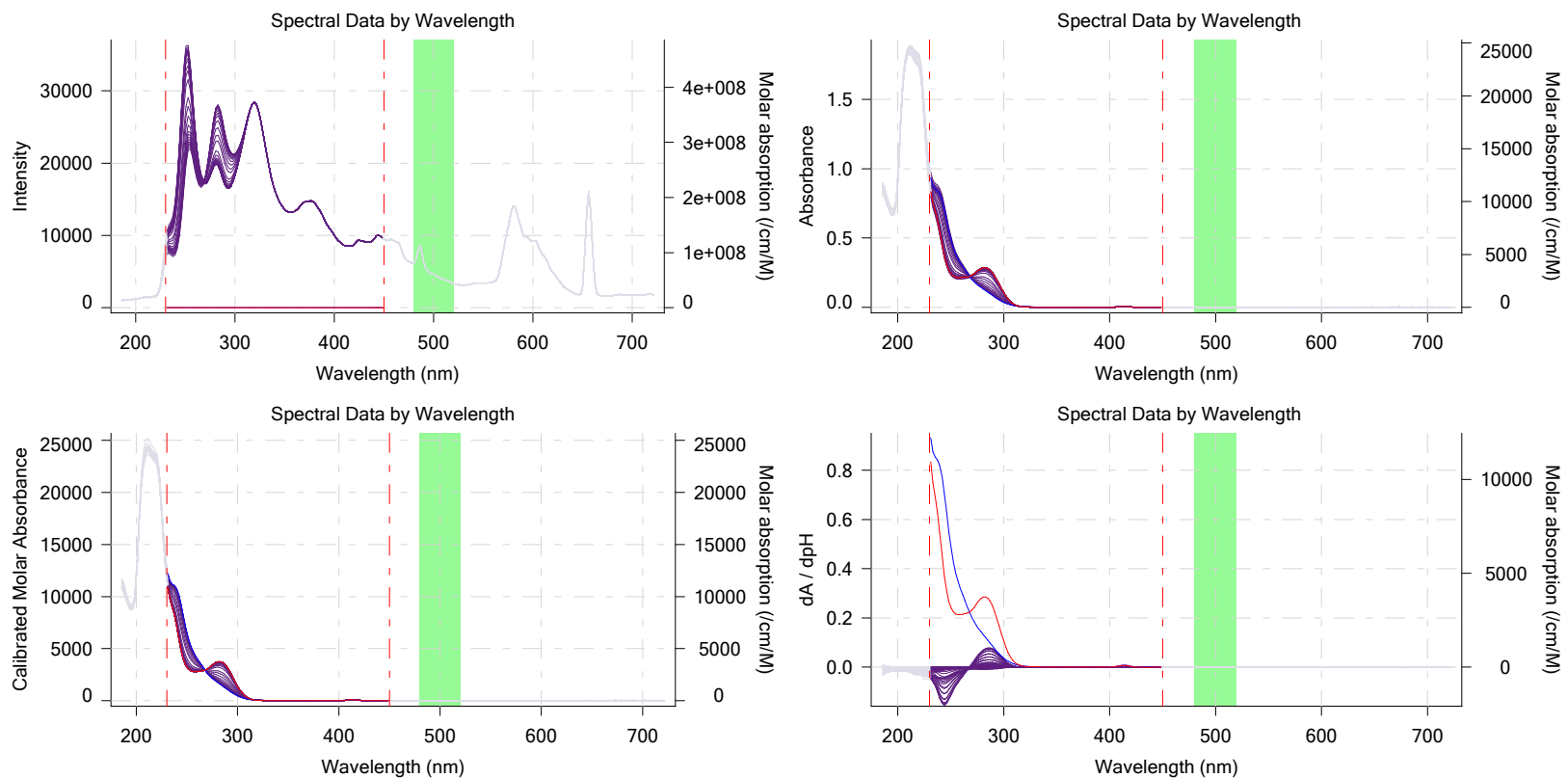

Sample name: **M11**
 Assay name: **UV-metric pKa**
 Assay ID: **17I-16007**
 Filename: **C:\Sirius_T3\17I-16007_M11_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 2:03:46 AM**
 Analyst: **Dorothy Leverse**
 Instrument ID: **T311053**


Graphs (continued)

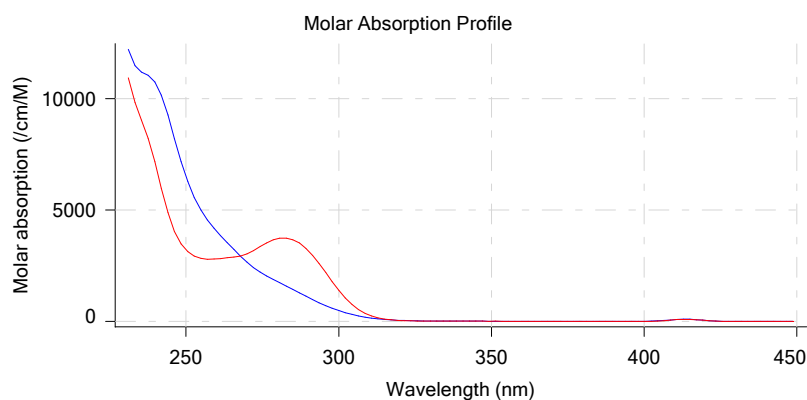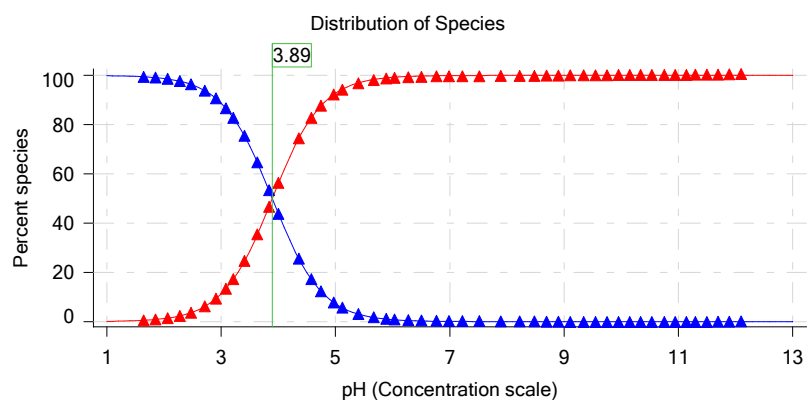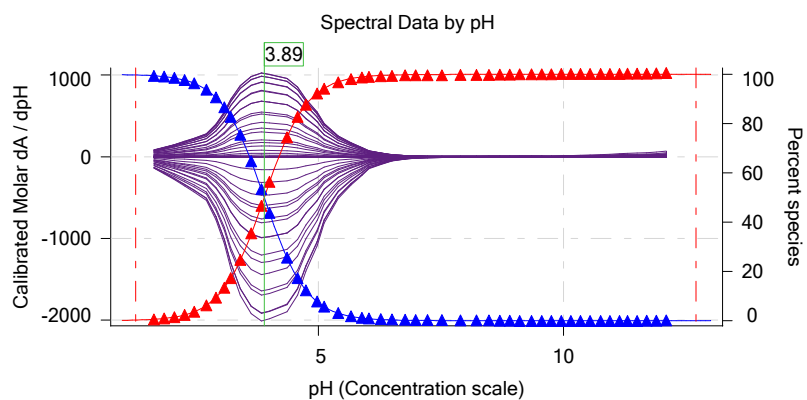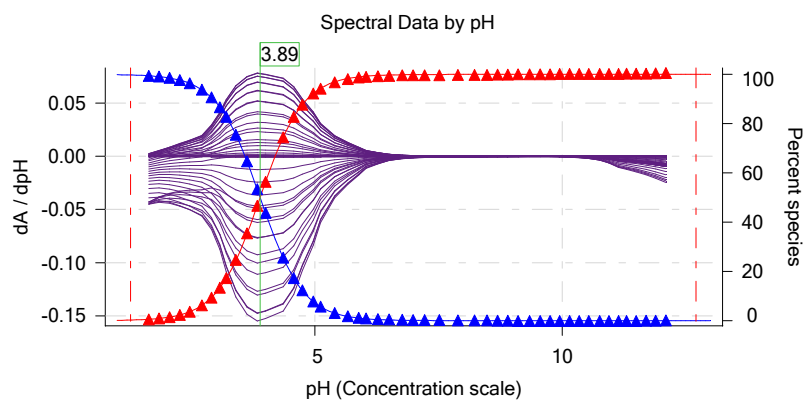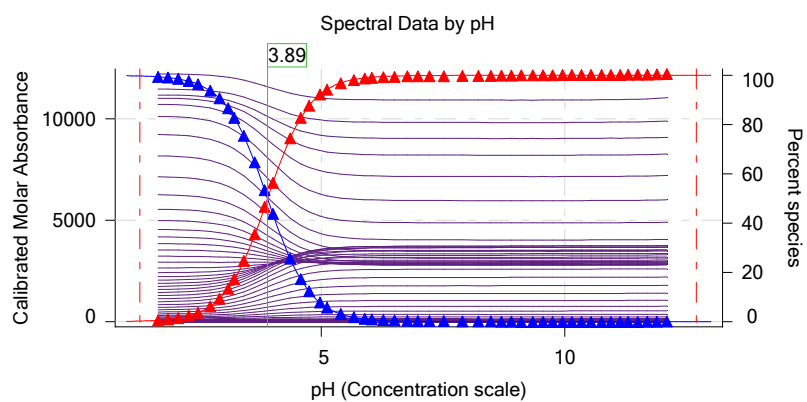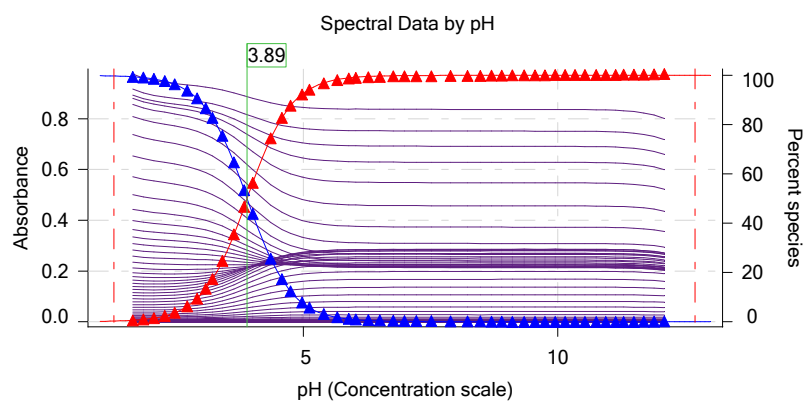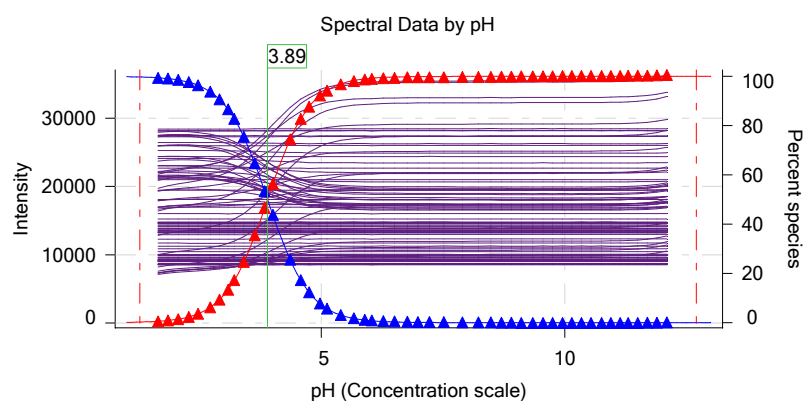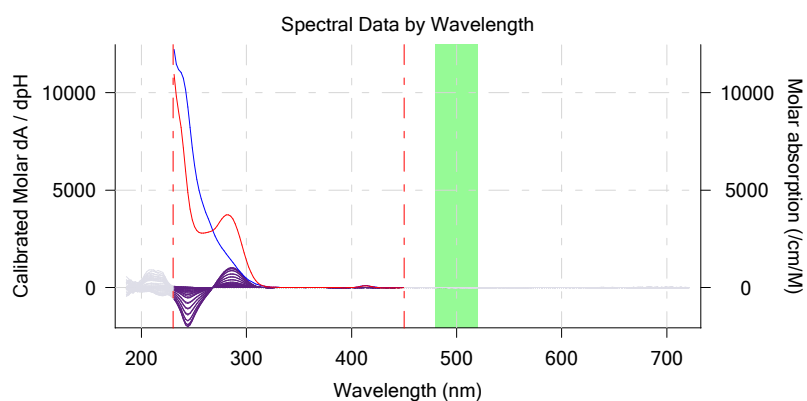

Sample name: **M11**
 Assay name: **UV-metric pKa**
 Assay ID: **171-16007**
 Filename: **C:\Sirius_T3\171-16007_M11_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 2:03:46 AM**
 Analyst: **Dorothy Levorse**
 Instrument ID: **T311053**

Graphs (continued)

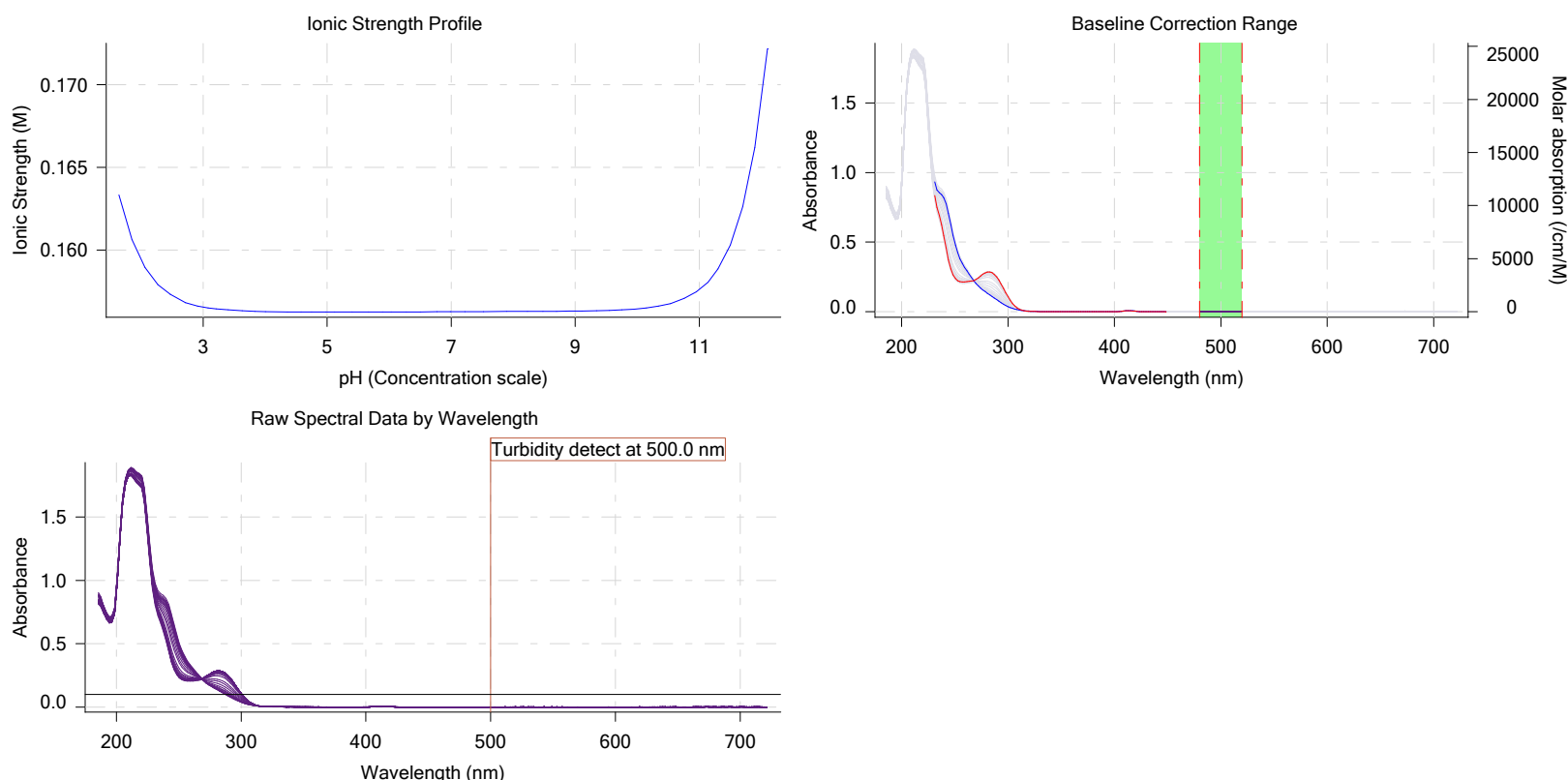

Assay Model

Settings	Value	Date/Time changed	Imported from
Sample name	M11	9/15/2017 3:37:20 PM	User entered value
Sample by	Volume		Default value
Sample volume	0.0020 mL	9/15/2017 3:37:20 PM	User entered value
Solvent	DMSO		Default value
Sample concentration	0.063900 M	9/15/2017 3:37:20 PM	User entered value
Solubility	Unknown		Default value
Molecular weight	211.22	9/15/2017 3:37:26 PM	User entered value
Individual pKa ionic environments	No		Default value
Number of pKas	1	9/15/2017 3:37:20 PM	User entered value
Sample is a	Base	9/15/2017 3:37:20 PM	User entered value
pKa 1	3.90	9/15/2017 3:37:20 PM	User entered value
logp (XH +)	-10.00		Default value
logP (neutral X)	-10.00	9/15/2017 3:37:20 PM	User entered value

Events

Time	Event	Water	Acid	Base	Buffer	pH	dpH/dt	pH R-squared	pH SD
3:13.7	Dark spectrum								
3:15.1	Reference spectrum								
3:42.8	Volume reset due to vial change								
5:13.2	Initial pH = 7.60								
6:26.2	Data point 4	1.50000 mL	0.07107 mL	0.00000 mL	0.02500 mL	1.771	-0.00634	0.42213	0.00049
6:55.1	Data point 5	1.50000 mL	0.07107 mL	0.02554 mL	0.02500 mL	1.972	0.01699	0.67702	0.00102
7:12.1	Data point 6	1.50000 mL	0.07107 mL	0.04238 mL	0.02500 mL	2.180	0.00317	0.05047	0.00070
7:29.0	Data point 7	1.50000 mL	0.07107 mL	0.05282 mL	0.02500 mL	2.390	0.01919	0.90293	0.00100
7:45.8	Data point 8	1.50000 mL	0.07107 mL	0.05917 mL	0.02500 mL	2.584	0.00343	0.41839	0.00026
8:02.5	Data point 9	1.50000 mL	0.07107 mL	0.06329 mL	0.02500 mL	2.829	0.00801	0.70833	0.00047

Sample name: **M11**
 Assay name: **UV-metric pKa**
 Assay ID: **17I-16007**
 Filename: **C:\Sirius_T3\17I-16007_M11_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 2:03:46 AM**
 Analyst: **Dorothy Leverse**
 Instrument ID: **T311053**

Events (continued)

Time	Event	Water	Acid	Base	Buffer	pH	dpH/dt	pH R-squared	pH SD	dpH/dt time
8:34.6	Data point 10	1.50000 mL	0.07107 mL	0.06569 mL	0.02500 mL	3.023	0.00374	0.38424	0.00030	10.0 s
8:51.2	Data point 11	1.50000 mL	0.07107 mL	0.06719 mL	0.02500 mL	3.196	0.00233	0.18069	0.00027	10.0 s
9:07.7	Data point 12	1.50000 mL	0.07107 mL	0.06820 mL	0.02500 mL	3.325	0.00157	0.09557	0.00025	10.0 s
9:29.5	Data point 13	1.50000 mL	0.07107 mL	0.06898 mL	0.02500 mL	3.520	0.00386	0.36858	0.00031	10.0 s
9:46.1	Data point 14	1.50000 mL	0.07107 mL	0.06947 mL	0.02500 mL	3.745	0.01139	0.83977	0.00061	10.0 s
10:02.8	Data point 15	1.50000 mL	0.07107 mL	0.06978 mL	0.02500 mL	3.950	0.01971	0.90933	0.00102	10.0 s
10:19.4	Data point 16	1.50000 mL	0.07107 mL	0.06999 mL	0.02500 mL	4.117	0.04048	0.93619	0.00206	10.0 s
10:41.1	Data point 17	1.50000 mL	0.07107 mL	0.07020 mL	0.02500 mL	4.474	0.06358	0.95154	0.00322	10.0 s
11:02.8	Data point 18	1.50000 mL	0.07107 mL	0.07030 mL	0.02500 mL	4.691	0.08976	0.92854	0.00460	10.5 s
11:19.9	Data point 19	1.50000 mL	0.07107 mL	0.07034 mL	0.02500 mL	4.861	0.09432	0.96465	0.00475	12.0 s
11:38.5	Data point 20	1.50000 mL	0.07107 mL	0.07039 mL	0.02500 mL	5.086	0.09966	0.99236	0.00494	15.0 s
12:00.1	Data point 21	1.50000 mL	0.07107 mL	0.07041 mL	0.02500 mL	5.234	0.09524	0.97979	0.00475	16.5 s
12:28.3	Data point 22	1.50000 mL	0.07107 mL	0.07046 mL	0.02500 mL	5.516	0.07663	0.88801	0.00404	10.0 s
12:50.0	Data point 23	1.50000 mL	0.07107 mL	0.07051 mL	0.02500 mL	5.782	0.04123	0.62718	0.00257	10.0 s
13:11.5	Data point 24	1.50000 mL	0.07107 mL	0.07056 mL	0.02500 mL	6.005	0.07077	0.69366	0.00419	10.0 s
13:28.0	Data point 25	1.50000 mL	0.07107 mL	0.07060 mL	0.02500 mL	6.149	0.08018	0.68398	0.00479	10.0 s
13:54.7	Data point 26	1.50000 mL	0.07107 mL	0.07067 mL	0.02500 mL	6.391	0.03834	0.63083	0.00238	10.0 s
14:21.6	Data point 27	1.50000 mL	0.07107 mL	0.07077 mL	0.02500 mL	6.617	0.03854	0.62565	0.00240	10.0 s
14:53.0	Data point 28	1.50000 mL	0.07107 mL	0.07086 mL	0.02500 mL	6.885	0.04543	0.76022	0.00257	10.0 s
15:19.8	Data point 29	1.50000 mL	0.07107 mL	0.07093 mL	0.02500 mL	7.101	0.07978	0.82147	0.00435	10.0 s
15:46.5	Data point 30	1.50000 mL	0.07107 mL	0.07100 mL	0.02500 mL	7.338	0.08191	0.85827	0.00436	11.0 s
16:14.2	Data point 31	1.50000 mL	0.07107 mL	0.07107 mL	0.02500 mL	7.634	0.09095	0.82314	0.00495	12.5 s
16:43.5	Data point 32	1.50000 mL	0.07107 mL	0.07114 mL	0.02500 mL	8.011	0.08253	0.83263	0.00446	14.0 s
17:14.3	Data point 33	1.50000 mL	0.07107 mL	0.07121 mL	0.02500 mL	8.343	0.07113	0.60187	0.00453	12.0 s
17:43.1	Data point 34	1.50000 mL	0.07107 mL	0.07128 mL	0.02500 mL	8.598	0.08000	0.78085	0.00447	11.5 s
18:11.1	Data point 35	1.50000 mL	0.07107 mL	0.07135 mL	0.02500 mL	8.814	0.08312	0.82481	0.00452	10.5 s
18:38.2	Data point 36	1.50000 mL	0.07107 mL	0.07143 mL	0.02500 mL	9.018	0.06528	0.87242	0.00345	10.0 s
19:04.8	Data point 37	1.50000 mL	0.07107 mL	0.07152 mL	0.02500 mL	9.225	0.04080	0.78610	0.00229	10.0 s
19:31.7	Data point 38	1.50000 mL	0.07107 mL	0.07166 mL	0.02500 mL	9.461	0.00943	0.35538	0.00078	10.0 s
19:58.5	Data point 39	1.50000 mL	0.07107 mL	0.07185 mL	0.02500 mL	9.670	0.00384	0.12824	0.00053	10.0 s
20:25.2	Data point 40	1.50000 mL	0.07107 mL	0.07206 mL	0.02500 mL	9.870	-0.00009	0.00012	0.00042	10.0 s
20:57.0	Data point 41	1.50000 mL	0.07107 mL	0.07237 mL	0.02500 mL	10.071	-0.00999	0.73722	0.00057	10.0 s
21:28.9	Data point 42	1.50000 mL	0.07107 mL	0.07279 mL	0.02500 mL	10.266	-0.01491	0.87516	0.00079	10.0 s
22:01.0	Data point 43	1.50000 mL	0.07107 mL	0.07342 mL	0.02500 mL	10.460	-0.01523	0.82876	0.00083	10.0 s
22:32.9	Data point 44	1.50000 mL	0.07107 mL	0.07441 mL	0.02500 mL	10.656	-0.01123	0.84422	0.00060	10.0 s
23:15.4	Data point 45	1.50000 mL	0.07107 mL	0.07646 mL	0.02500 mL	10.874	-0.01429	0.94540	0.00073	10.0 s
23:47.5	Data point 46	1.50000 mL	0.07107 mL	0.07874 mL	0.02500 mL	11.066	-0.01387	0.90225	0.00072	10.0 s
24:14.6	Data point 47	1.50000 mL	0.07107 mL	0.08154 mL	0.02500 mL	11.258	-0.01233	0.79498	0.00068	10.0 s
24:31.3	Data point 48	1.50000 mL	0.07107 mL	0.08615 mL	0.02500 mL	11.418	-0.01128	0.81642	0.00062	10.0 s
25:03.6	Data point 49	1.50000 mL	0.07107 mL	0.09426 mL	0.02500 mL	11.614	-0.01041	0.82694	0.00057	10.0 s
25:36.0	Data point 50	1.50000 mL	0.07107 mL	0.10767 mL	0.02500 mL	11.812	-0.00887	0.79144	0.00049	10.0 s
26:03.4	Data point 51	1.50000 mL	0.07107 mL	0.12843 mL	0.02500 mL	12.004	-0.00564	0.47854	0.00040	10.0 s
26:36.3	Data point 52	1.50000 mL	0.07107 mL	0.16470 mL	0.02500 mL	12.204	-0.00466	0.59069	0.00030	10.0 s
28:36.5	Assay volumes	1.75000 mL	0.24516 mL	0.16470 mL	0.02500 mL					

Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
General Settings				
Analyst name	Dorothy Leverse			
Separate reference vial	Yes			
Standard Experiment Settings				
Number of titrations	1			
Minimum pH	1.800			
Maximum pH	12.200			
pH step between points of	0.200			
Minimum titrant addition	0.00002 mL			
Maximum titrant addition	0.10000 mL			

Sample name: **M11**
 Assay name: **UV-metric pKa**
 Assay ID: **171-16007**
 Filename: **C:\Sirius_T3\171-16007_M11_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 2:03:46 AM**
 Analyst: **Dorothy Levorse**
 Instrument ID: **T311053**

Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
Argon flow rate	100%			
Start titration using	Cautious pH adjust			
Advanced General Settings				
Detect turbidity using	Spectrometer			
Monitor at a wavelength of	500.0 nm			
Absorbance threshold of	0.100			
Collect turbidity sensor data	No			
Stir after titrant addition for	5 seconds			
For titrant addition, stir at	15%			
Titrant Pre-Dose				
Titrant pre-dose	None			
Assay Medium				
Cosolvent in use	No			
ISA water volume	1.50 mL			
Water added	Automatic			
After water addition, stir for	5 seconds			
At a speed of	15%			
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			
After medium addition, stir for	5 seconds			
Sample Sonication				
Sonicate	No			
Sample Dissolution				
Perform a dissolution stage	No			
Carbonate purge				
Perform a carbonate purge	No			
Temperature Control				
Wait for temperature	Yes			
Required start temperature	25.0°C			
Acceptable deviation	0.5°C			
Time to wait	60 seconds			
Stir speed of	15%			
Titration 1				
Titrate from	Low to high pH			
Adjust to start pH	Yes			
After pH adjust stir for	10 seconds			
Data Point Stability				
Stir during data point collection	Yes			
For point collection, stir at	15%			
Delay before data point collection	0 seconds			
Number of points to average	20 points			
Time interval between points	0.50 seconds			
Required maximum standard deviation	0.00500 dpH/dt			
Stability timeout after	60 seconds			
Experiment cleanup				
Adjust pH to cleanup	To start pH			
And then stir for	60 seconds			
For cleaning, stir at	20%			
Then add water volume	0.25 mL			
And then stir for	30 seconds			

Calibration Settings

Setting	Value	Date/Time changed	Imported from
Four-Plus alpha	0.112	9/16/2017 2:03:46 AM	C:\Sirius_T3\HC117115.t3r
Four-Plus S	1.0006	9/16/2017 2:03:46 AM	C:\Sirius_T3\HC117115.t3r



Assay Settings

Sample name: **M11**
Assay name: **UV-metric pKa**
Assay ID: **17I-16007**
Filename: **C:\Sirius_T3\17I-16007_M11_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 2:03:46 AM**
Analyst: **Dorothy Levorse**
Instrument ID: **T311053**

Calibration Settings (continued)

Setting	Value	Date/Time changed	Imported from
Four-Plus jH	0.7	9/16/2017 2:03:46 AM	C:\Sirius_T3\HCl17I15.t3r
Four-Plus jOH	-0.6	9/16/2017 2:03:46 AM	C:\Sirius_T3\HCl17I15.t3r
Base concentration factor	1.015	9/16/2017 2:03:46 AM	C:\Sirius_T3\KOH17I11.t3r
Acid concentration factor	1.003	9/16/2017 2:03:46 AM	C:\Sirius_T3\HCl17I15.t3r

Instrument Settings

Setting	Value	Batch Id	Install date
Instrument owner	Merck		
Instrument ID	T311053		
Instrument type	T3 Simulator		
Software version	1.1.3.0		
Dispenser module		T3DM1100253	3/31/2009 5:24:52 AM
Dispenser 0	Water		3/31/2009 5:25:05 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Water (0.15 M KCl)	8-18-17	9/8/2017 8:22:43 AM
Dispenser 2	Acid		3/31/2009 5:25:11 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Acid (0.5 M HCl)	166940	9/8/2017 8:21:27 AM
Dispenser 1	Base		3/31/2009 5:25:21 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Base (0.5 M KOH)	01/06/17	9/8/2017 8:20:03 AM
Dispenser 5	Cosolvent		3/31/2009 5:26:24 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Distribution valve 5	Distribution Valve		3/31/2009 5:28:19 AM
Firmware version	1.1.3		
Port A	Methanol (80%, 0.15 M KCl)	8-15-17	9/13/2017 11:23:11 AM
Dispenser 3	Buffer		8/3/2010 5:05:16 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Phosphate Buffer		9/12/2017 11:32:29 AM
Dispenser 6	Octanol		10/22/2010 10:52:43 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Octanol	9-14-17	9/14/2017 9:30:38 AM
Titration		T3TM1100153	3/31/2009 5:24:17 AM
Horizontal axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Probe I/O firmware version	1.1.1		
Electrode	T3 Electrode	T3E0769	8/15/2017 9:21:54 AM
E0 calibration	-8.30 mV		9/16/2017 2:04:10 AM
Filling solution	3M KCl	KCL095	9/13/2017 8:16:19 AM
Liquids			
Wash 1	50% IPA:50% Water		9/15/2017 8:38:18 AM
Wash 2	0.5% Triton X-100 in H2O		9/15/2017 8:38:22 AM
Buffer position 1	pH7 Wash		9/15/2017 8:38:24 AM
Buffer position 2	pH 7		9/15/2017 8:38:27 AM
Storage position			9/15/2017 8:38:55 AM
Wash water	3.3e+003 mL	9-11-17	9/11/2017 3:28:43 PM
Waste	6.8e+003 mL		9/11/2017 3:28:49 PM
Temperature controller			8/5/2010 6:35:13 AM
Turbidity detector			3/31/2009 5:24:45 AM
Spectrometer		072390	11/23/2010 11:22:28 AM



Assay Settings

Sample name: **M11**
Assay name: **UV-metric pKa**
Assay ID: **171-16007**
Filename: **C:\Sirius_T3\171-16007_M11_UV-metric pKa.t3r**

Experiment start time: **9/16/2017 2:03:46 AM**
Analyst: **Dorothy Levorse**
Instrument ID: **T311053**

Instrument Settings (continued)

Setting	Value	Batch Id	Install date
Dip probe		11086	
Wavelength coefficient A0	185.563		
Wavelength coefficient A1	2.17439		
Wavelength coefficient A2	-0.000285622		
Total lamp lit time	114:03:31		11/23/2010 11:22:28 AM
Calibrated on	9/6/2017 8:33:02 AM		
Integration time	11		
Scans averaged	10		
Autoloader		T3AL1100237	11/10/2015 9:34:13 AM
Left-right axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Front-back axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Configuration			
Alternate titration position	Titration position		
Alternate reference position	Reference position		
Maximum standard vial volume	3.50 mL		
Maximum alternate vial volume	25.00 mL		
Automatic action idle period	5 minute(s)		
Titrant tube volume	1.3 mL		
Syringe flush count	3.50		
Flowing wash pump volume	20.0 mL		
Flowing wash stir duration	5 s		
Flowing wash stir speed	30%		
Solvent wash stir duration	5 s		
Solvent wash stir speed	30%		
Surfactant wash stir duration	5 s		
Surfactant wash stir speed	30%		
E0 calibration minimum number of points	10		
E0 calibration maximum standard deviation	0.01500		
E0 calibration timeout period	60 s		
E0 calibration stir duration	5 s		
E0 calibration preparation stir speed	30%		
E0 calibration buffer wash stir duration	5 s		
E0 calibration buffer wash stir speed	30%		
E0 calibration reading stir speed	0%		
Spectrometer calibration stir duration	5 s		
Spectrometer calibration stir speed	30%		
Spectrometer calibration wash pump volume	20.0 mL		
Spectrometer calibration wash stir duration	5 s		
Spectrometer calibration wash stir speed	30%		
Overhead dispense height	10000		

Refinement Settings

Setting	Value	Default value
Turbidity detection method	Spectrometer	Spectrometer
Turbidity wavelength to assess	500.0 nm	500.0 nm
Turbidity maximum absorbance	0.100	0.100
Turbidity probe threshold	50.00	50.00
Exclude turbid points	Yes	Yes
Low intensity warning threshold	100	100
Minimum absorbance change threshold	0.100	0.100
Eigenvector autocorrelation threshold	0.80	0.80
Maximum RMSD severe warning	0.250	0.250
Maximum RMSD warning	0.050	0.050