

Sample name: **M16**  
 Assay name: **UV-metric pKa**  
 Assay ID: **17K-10009**  
 Filename: **C:\Sirius\_T3\Mehtap\20171110\_exp17\_pKa\17K-10009\_M16\_UV-metric pKa.t3r**

Experiment start time: **11/10/2017 4:58:04 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Results

pKa 1 **4.71**  
 pKa 2 **8.96**  
 RMSD **0.010 0.003 0.005**  
 Chi squared **0.1967**  
 PCA calculated number of pKas **3**  
 Average ionic strength **0.158 M**  
 Average temperature **24.9°C**  
 Analyte concentration range **62.6 µM to 57.2 µM**

Number of pKas source **Predicted**  
 Wavelength clipping **230.0 nm to 450.0 nm**  
 pH clipping **1.279 to 12.748**


## Warnings and errors

Errors **None**  
 Warnings **PCA calculation disagrees with predicted number of pKas**

## Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
<b>Assay Medium</b>				
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			









## Graphs

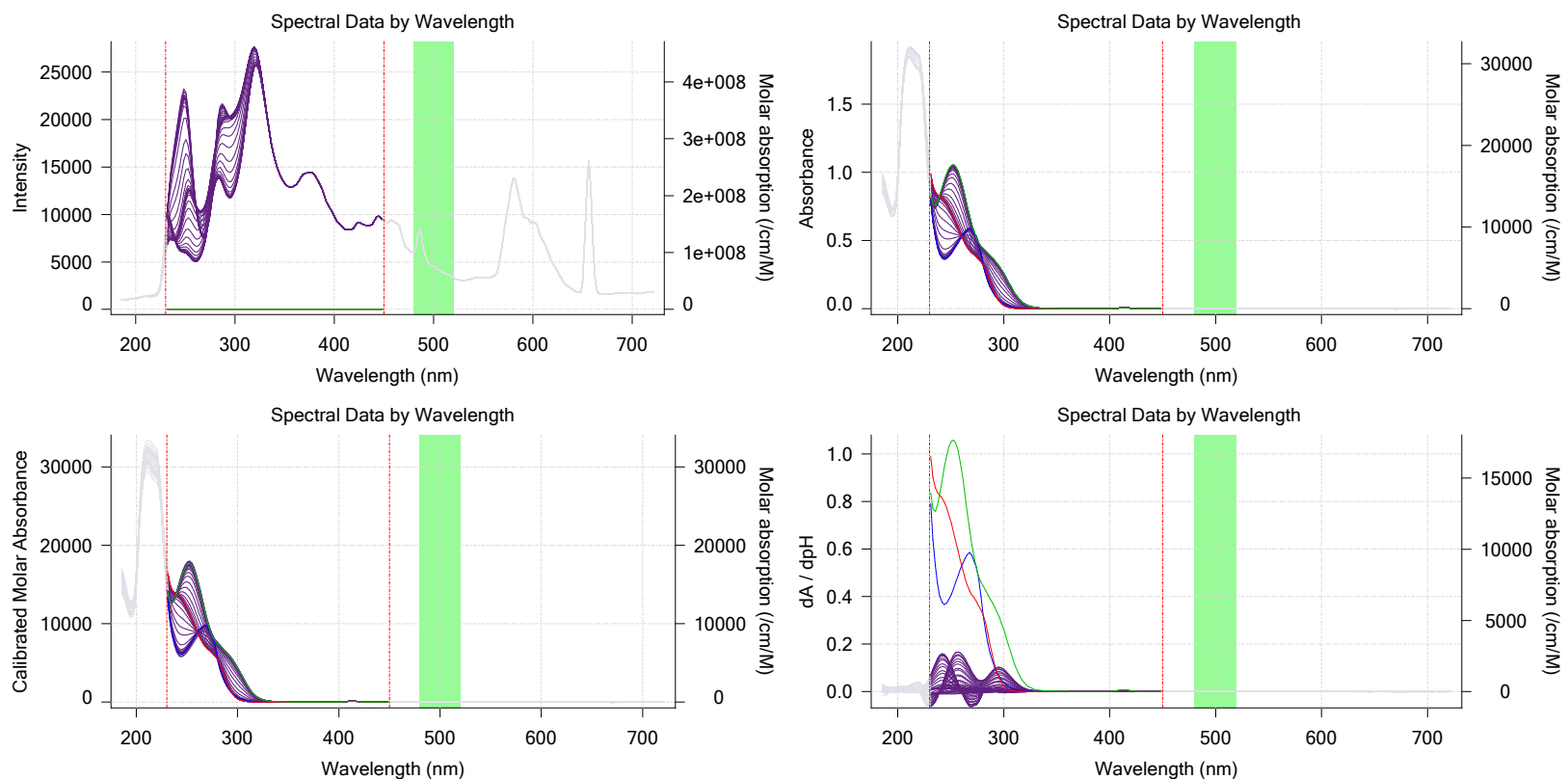

Sample name: **M16**  
 Assay name: **UV-metric pKa**  
 Assay ID: **17K-10009**  
 Filename: **C:\Sirius\_T3\Mehtap\20171110\_exp17\_pKa\17K-10009\_M16\_UV-metric pKa.t3r**

Experiment start time: **11/10/2017 4:58:04 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**


## Graphs (continued)

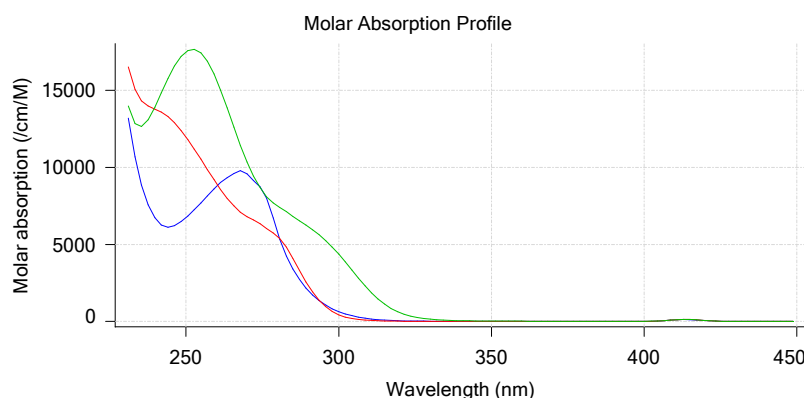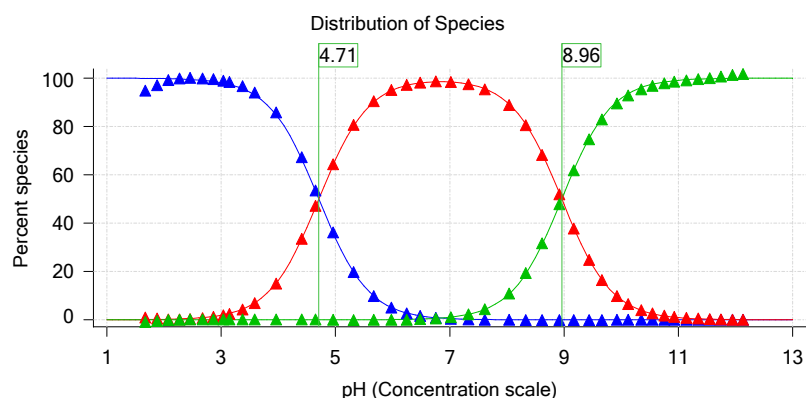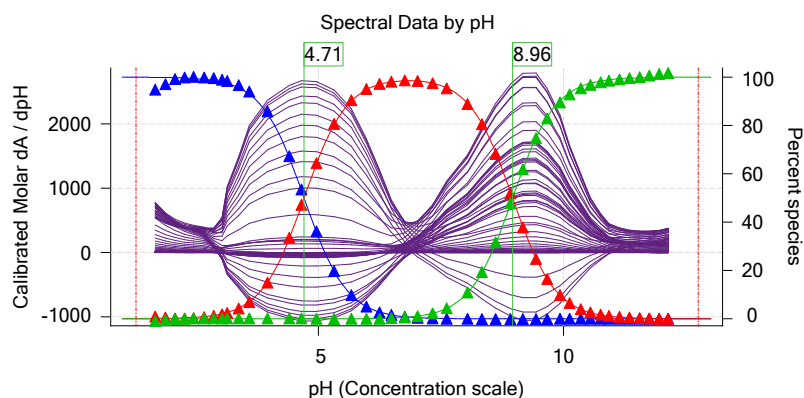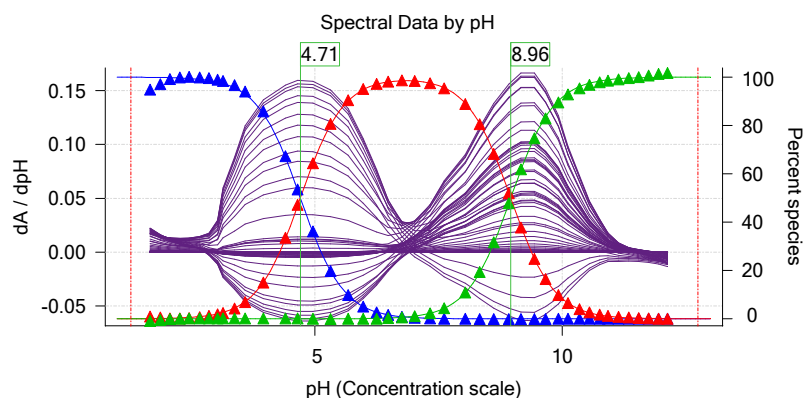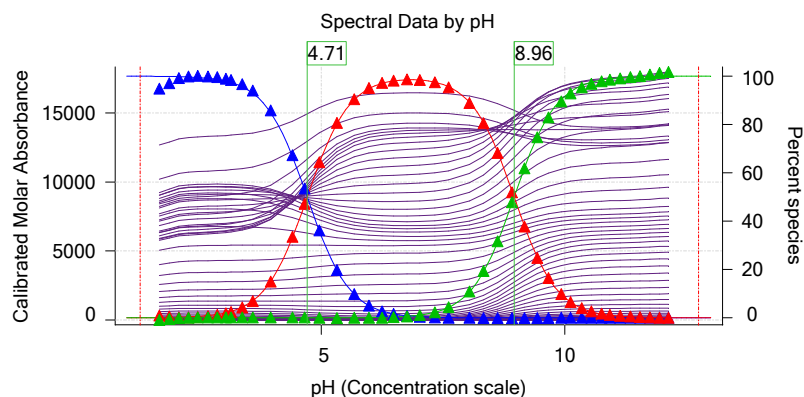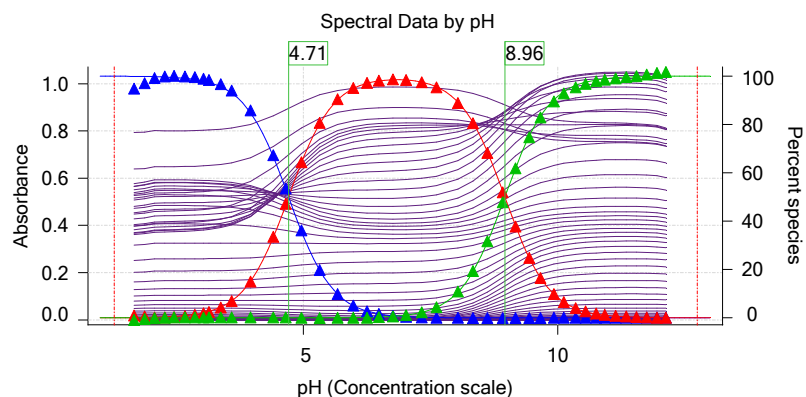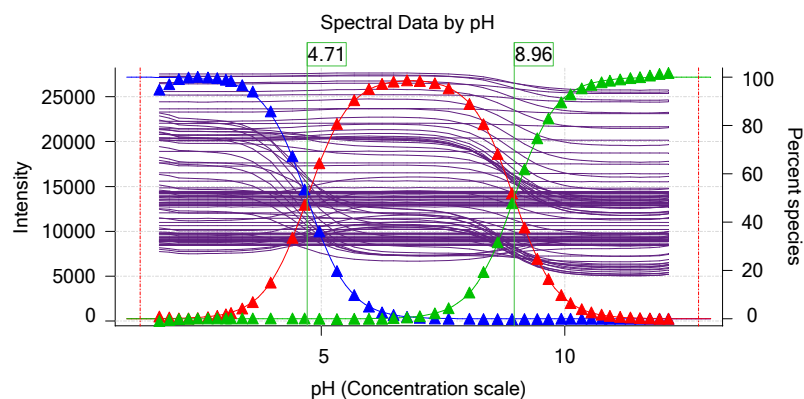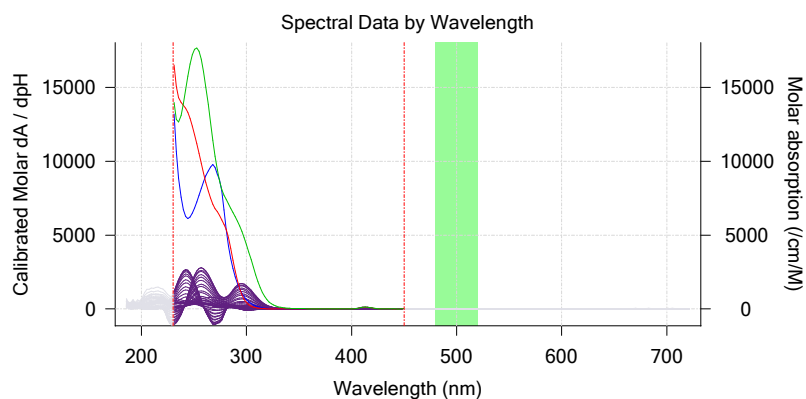

Sample name: **M16**  
 Assay name: **UV-metric pKa**  
 Assay ID: **17K-10009**  
 Filename: **C:\Sirius\_T3\Mehtap\20171110\_exp17\_pKa\17K-10009\_M16\_UV-metric pKa.t3r**

Experiment start time: **11/10/2017 4:58:04 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Graphs (continued)

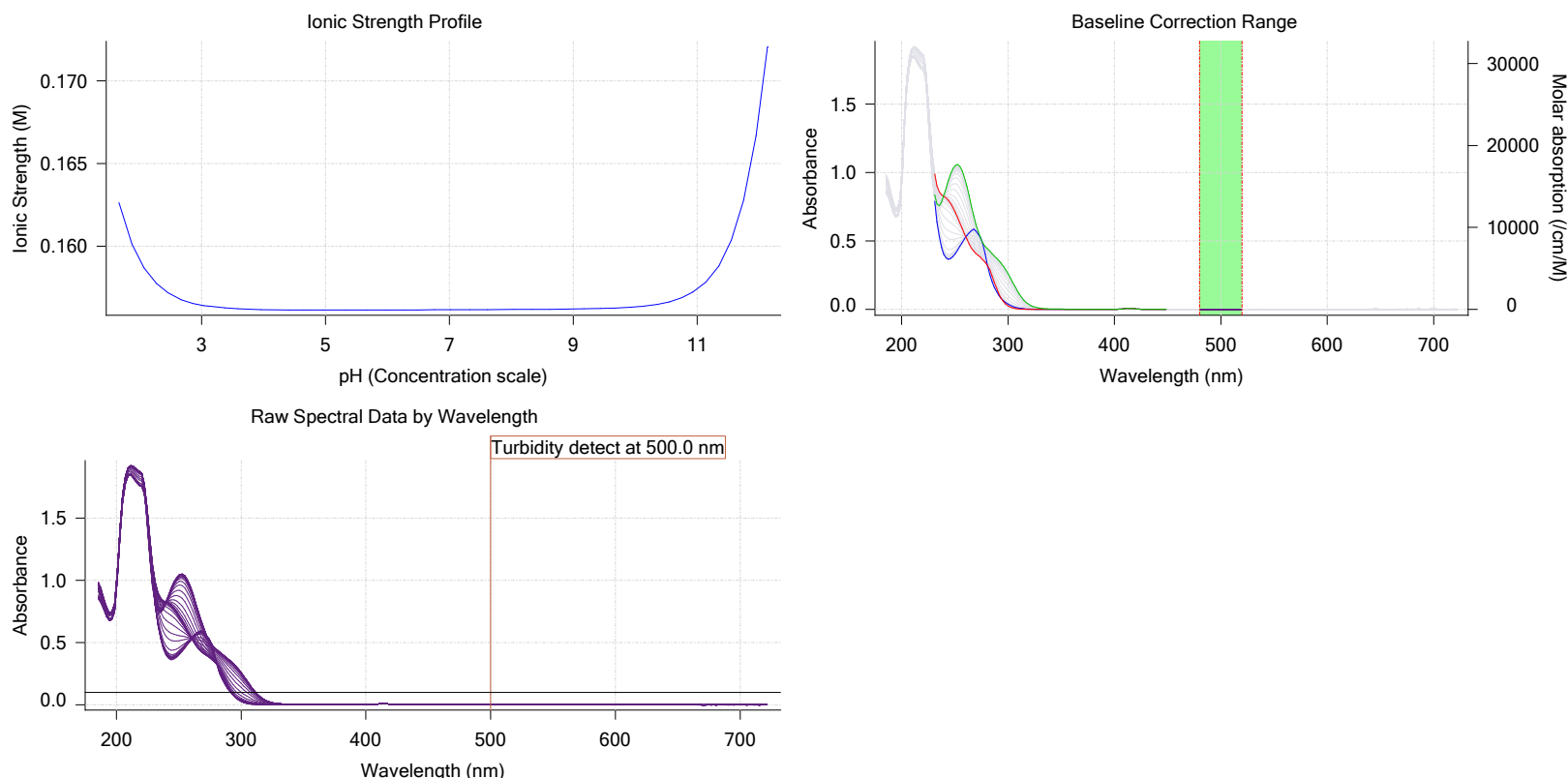

## Assay Model

Settings	Value	Date/Time changed	Imported from
Sample name	M16	11/9/2017 5:28:23 PM	User entered value
Sample by	Volume		Default value
Sample volume	0.0020 mL	11/9/2017 5:28:23 PM	User entered value
Solvent	DMSO		Default value
Sample concentration	0.050000 M	11/9/2017 5:28:23 PM	User entered value
Solubility	Unknown		Default value
Molecular weight	210.23	11/9/2017 5:28:29 PM	User entered value
Individual pKa ionic environments	No		Default value
Number of pKas	2	11/9/2017 5:28:23 PM	User entered value
Sample is a	Ampholyte	11/9/2017 5:28:23 PM	User entered value
pKa 1	3.86	11/9/2017 5:28:23 PM	User entered value
Type	Base	11/9/2017 5:28:23 PM	User entered value
pKa 2	9.07	11/9/2017 5:28:23 PM	User entered value
Type	Acid	11/9/2017 5:28:23 PM	User entered value
logp (XH <sub>2</sub> +)	-10.00		Default value
logP (neutral XH)	-10.00	11/9/2017 5:28:23 PM	User entered value
logP (X -)	-10.00		Default value

## Events

Time	Event	Water	Acid	Base	Buffer	pH	dpH/dt	pH R-squared	pH SD
2:58.5	Dark spectrum								
2:59.8	Reference spectrum								
3:27.5	Volume reset due to vial change								
4:57.9	Initial pH = 7.05								
6:10.8	Data point 4	1.50000 mL	0.07098 mL	0.00000 mL	0.02500 mL	1.779	-0.00596	0.52455	0.00041
6:39.4	Data point 5	1.50000 mL	0.07098 mL	0.02561 mL	0.02500 mL	1.980	0.00947	0.75147	0.00054

Sample name: **M16**  
 Assay name: **UV-metric pKa**  
 Assay ID: **17K-10009**  
 Filename: **C:\Sirius\_T3\Mehtap\20171110\_exp17\_pKa\17K-10009\_M16\_UV-metric pKa.t3r**

Experiment start time: **11/10/2017 4:58:04 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Events (continued)

Time	Event	Water	Acid	Base	Buffer	pH	dpH/dt	pH R-squared	pH SD	dpH/dt time
6:56.3	Data point 6	1.50000 mL	0.07098 mL	0.04160 mL	0.02500 mL	2.170	0.00067	0.00223	0.00070	10.0 s
7:13.0	Data point 7	1.50000 mL	0.07098 mL	0.05191 mL	0.02500 mL	2.368	0.01704	0.87061	0.00090	10.0 s
7:29.7	Data point 8	1.50000 mL	0.07098 mL	0.05835 mL	0.02500 mL	2.551	0.00897	0.81075	0.00049	10.0 s
7:46.3	Data point 9	1.50000 mL	0.07098 mL	0.06258 mL	0.02500 mL	2.766	0.01047	0.85459	0.00056	10.0 s
8:02.9	Data point 10	1.50000 mL	0.07098 mL	0.06517 mL	0.02500 mL	2.956	0.01383	0.81033	0.00076	10.0 s
8:19.6	Data point 11	1.50000 mL	0.07098 mL	0.06684 mL	0.02500 mL	3.129	0.01221	0.90712	0.00063	10.0 s
8:36.3	Data point 12	1.50000 mL	0.07098 mL	0.06797 mL	0.02500 mL	3.235	0.00390	0.31117	0.00035	10.0 s
9:08.2	Data point 13	1.50000 mL	0.07098 mL	0.06912 mL	0.02500 mL	3.468	0.01408	0.86584	0.00075	10.0 s
9:34.9	Data point 14	1.50000 mL	0.07098 mL	0.06968 mL	0.02500 mL	3.683	0.02123	0.92825	0.00109	10.0 s
9:51.5	Data point 15	1.50000 mL	0.07098 mL	0.07001 mL	0.02500 mL	4.056	0.03786	0.94143	0.00193	10.0 s
10:13.2	Data point 16	1.50000 mL	0.07098 mL	0.07030 mL	0.02500 mL	4.506	0.07198	0.93124	0.00368	10.0 s
10:39.9	Data point 17	1.50000 mL	0.07098 mL	0.07044 mL	0.02500 mL	4.752	0.09876	0.98446	0.00491	13.0 s
11:09.6	Data point 18	1.50000 mL	0.07098 mL	0.07053 mL	0.02500 mL	5.059	0.09986	0.97933	0.00498	17.5 s
11:43.7	Data point 19	1.50000 mL	0.07098 mL	0.07060 mL	0.02500 mL	5.421	0.09859	0.96157	0.00496	22.0 s
12:22.3	Data point 20	1.50000 mL	0.07098 mL	0.07067 mL	0.02500 mL	5.775	0.10059	0.99523	0.00497	16.0 s
12:55.1	Data point 21	1.50000 mL	0.07098 mL	0.07074 mL	0.02500 mL	6.083	0.09577	0.92782	0.00491	12.5 s
13:24.2	Data point 22	1.50000 mL	0.07098 mL	0.07081 mL	0.02500 mL	6.354	0.09132	0.95413	0.00461	11.5 s
13:52.5	Data point 23	1.50000 mL	0.07098 mL	0.07088 mL	0.02500 mL	6.590	0.09421	0.92472	0.00483	11.0 s
14:25.3	Data point 24	1.50000 mL	0.07098 mL	0.07098 mL	0.02500 mL	6.866	0.09713	0.96447	0.00488	11.5 s
14:58.7	Data point 25	1.50000 mL	0.07098 mL	0.07107 mL	0.02500 mL	7.130	0.09005	0.93588	0.00459	12.0 s
15:32.5	Data point 26	1.50000 mL	0.07098 mL	0.07117 mL	0.02500 mL	7.438	0.09571	0.97126	0.00479	14.5 s
16:03.7	Data point 27	1.50000 mL	0.07098 mL	0.07124 mL	0.02500 mL	7.726	0.09962	0.99139	0.00494	19.0 s
16:39.4	Data point 28	1.50000 mL	0.07098 mL	0.07131 mL	0.02500 mL	8.152	0.09514	0.97607	0.00475	18.5 s
17:09.3	Data point 29	1.50000 mL	0.07098 mL	0.07135 mL	0.02500 mL	8.446	0.09076	0.93138	0.00464	17.0 s
17:37.9	Data point 30	1.50000 mL	0.07098 mL	0.07140 mL	0.02500 mL	8.733	0.09473	0.95105	0.00479	17.5 s
18:07.0	Data point 31	1.50000 mL	0.07098 mL	0.07147 mL	0.02500 mL	9.034	0.09309	0.85884	0.00495	12.5 s
18:31.1	Data point 32	1.50000 mL	0.07098 mL	0.07154 mL	0.02500 mL	9.287	0.09508	0.94816	0.00482	11.0 s
18:48.5	Data point 33	1.50000 mL	0.07098 mL	0.07164 mL	0.02500 mL	9.556	0.06851	0.93049	0.00350	10.0 s
19:15.4	Data point 34	1.50000 mL	0.07098 mL	0.07178 mL	0.02500 mL	9.783	0.02516	0.80757	0.00138	10.0 s
19:31.9	Data point 35	1.50000 mL	0.07098 mL	0.07197 mL	0.02500 mL	10.046	0.00319	0.17404	0.00038	10.0 s
20:03.7	Data point 36	1.50000 mL	0.07098 mL	0.07230 mL	0.02500 mL	10.241	0.00193	0.08188	0.00033	10.0 s
20:20.2	Data point 37	1.50000 mL	0.07098 mL	0.07277 mL	0.02500 mL	10.472	-0.00473	0.58499	0.00031	10.0 s
20:52.1	Data point 38	1.50000 mL	0.07098 mL	0.07361 mL	0.02500 mL	10.670	-0.00845	0.80992	0.00046	10.0 s
21:08.6	Data point 39	1.50000 mL	0.07098 mL	0.07479 mL	0.02500 mL	10.878	-0.00469	0.49129	0.00033	10.0 s
21:25.2	Data point 40	1.50000 mL	0.07098 mL	0.07667 mL	0.02500 mL	11.050	-0.00470	0.56451	0.00031	10.0 s
21:41.8	Data point 41	1.50000 mL	0.07098 mL	0.07947 mL	0.02500 mL	11.263	-0.00565	0.58631	0.00036	10.0 s
21:58.5	Data point 42	1.50000 mL	0.07098 mL	0.08403 mL	0.02500 mL	11.470	-0.00439	0.45909	0.00032	10.0 s
22:15.2	Data point 43	1.50000 mL	0.07098 mL	0.09144 mL	0.02500 mL	11.673	-0.00494	0.46894	0.00036	10.0 s
22:32.1	Data point 44	1.50000 mL	0.07098 mL	0.10339 mL	0.02500 mL	11.871	-0.00319	0.22828	0.00033	10.0 s
22:49.0	Data point 45	1.50000 mL	0.07098 mL	0.12260 mL	0.02500 mL	12.070	0.00174	0.14798	0.00022	10.0 s
23:06.2	Data point 46	1.50000 mL	0.07098 mL	0.14995 mL	0.02500 mL	12.248	-0.00044	0.00685	0.00026	10.0 s
25:11.4	Assay volumes	1.75000 mL	0.23676 mL	0.14995 mL	0.02500 mL					

## Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
<b>General Settings</b>				
Analyst name	Dorothy Levorse			
Separate reference vial	Yes			
<b>Standard Experiment Settings</b>				
Number of titrations	1			
Minimum pH	1.800			
Maximum pH	12.200			
pH step between points of	0.200			
Minimum titrant addition	0.00002 mL			
Maximum titrant addition	0.10000 mL			
Argon flow rate	100%			
Start titration using	Cautious pH adjust			

Sample name: **M16**  
 Assay name: **UV-metric pKa**  
 Assay ID: **17K-10009**  
 Filename: **C:\Sirius\_T3\Mehtap\20171110\_exp17\_pKa\17K-10009\_M16\_UV-metric pKa.t3r**

Experiment start time: **11/10/2017 4:58:04 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
<b>Advanced General Settings</b>				
Detect turbidity using	Spectrometer			
Monitor at a wavelength of	500.0 nm			
Absorbance threshold of	0.100			
Collect turbidity sensor data	No			
Stir after titrant addition for	5 seconds			
For titrant addition, stir at	15%			
<b>Titrant Pre-Dose</b>				
Titrant pre-dose	None			
<b>Assay Medium</b>				
Cosolvent in use	No			
ISA water volume	1.50 mL			
Water added	Automatic			
After water addition, stir for	5 seconds			
At a speed of	15%			
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			
After medium addition, stir for	5 seconds			
<b>Sample Sonication</b>				
Sonicate	No			
<b>Sample Dissolution</b>				
Perform a dissolution stage	No			
<b>Carbonate purge</b>				
Perform a carbonate purge	No			
<b>Temperature Control</b>				
Wait for temperature	Yes			
Required start temperature	25.0°C			
Acceptable deviation	0.5°C			
Time to wait	60 seconds			
Stir speed of	15%			
<b>Titration 1</b>				
Titrate from	Low to high pH			
Adjust to start pH	Yes			
After pH adjust stir for	10 seconds			
<b>Data Point Stability</b>				
Stir during data point collection	Yes			
For point collection, stir at	15%			
Delay before data point collection	0 seconds			
Number of points to average	20 points			
Time interval between points	0.50 seconds			
Required maximum standard deviation	0.00500 dpH/dt			
Stability timeout after	60 seconds			
<b>Experiment cleanup</b>				
Adjust pH to cleanup	To start pH			
And then stir for	60 seconds			
For cleaning, stir at	20%			
Then add water volume	0.25 mL			
And then stir for	30 seconds			

## Calibration Settings

Setting	Value	Date/Time changed	Imported from
Four-Plus alpha	0.081	11/10/2017 4:58:04 PM	C:\Sirius_T3\HCl17K10.t3r
Four-Plus S	1.0039	11/10/2017 4:58:04 PM	C:\Sirius_T3\HCl17K10.t3r
Four-Plus jH	1.1	11/10/2017 4:58:04 PM	C:\Sirius_T3\HCl17K10.t3r
Four-Plus jOH	-0.5	11/10/2017 4:58:04 PM	C:\Sirius_T3\HCl17K10.t3r



Sample name: **M16**  
 Assay name: **UV-metric pKa**  
 Assay ID: **17K-10009**  
 Filename: **C:\Sirius\_T3\Mehtap\20171110\_exp17\_pKa\17K-10009\_M16\_UV-metric pKa.t3r**

Experiment start time: **11/10/2017 4:58:04 PM**  
 Analyst: **Dorothy Levorse**  
 Instrument ID: **T311053**

## Calibration Settings (continued)

Setting	Value	Date/Time changed	Imported from
Base concentration factor	1.008	11/10/2017 4:58:04 PM	C:\Sirius_T3\KOH17K09.t3r
Acid concentration factor	0.999	11/10/2017 4:58:04 PM	C:\Sirius_T3\HCl17K10.t3r

## Instrument Settings

Setting	Value	Batch Id	Install date
Instrument owner	Merck		
Instrument ID	T311053		
Instrument type	T3 Simulator		
Software version	1.1.3.0		
Dispenser module		T3DM1100253	3/31/2009 5:24:52 AM
Dispenser 0	Water		3/31/2009 5:25:05 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Water (0.15 M KCl)	10-10-2017	11/8/2017 11:33:30 AM
Dispenser 2	Acid		3/31/2009 5:25:11 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Acid (0.5 M HCl)	11-8-17	11/8/2017 11:32:21 AM
Dispenser 1	Base		3/31/2009 5:25:21 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Base (0.5 M KOH)	10-30-17	10/30/2017 8:01:46 AM
Dispenser 5	Cosolvent		3/31/2009 5:26:24 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Distribution valve 5	Distribution Valve		3/31/2009 5:28:19 AM
Firmware version	1.1.3		
Port A	Methanol (80%, 0.15 M KCl)	9-26-17	11/1/2017 10:56:16 AM
Port B	Cyclohexane		10/19/2017 2:11:05 PM
Port C	MeCN (50%, 0.15 M KCl)	10-30-17	10/30/2017 8:02:00 AM
Dispenser 3	Buffer		8/3/2010 5:05:16 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Phosphate Buffer		11/8/2017 11:32:27 AM
Dispenser 6	Octanol		10/22/2010 10:52:43 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Octanol	9-14-17	10/13/2017 7:46:59 AM
Titration		T3TM1100153	3/31/2009 5:24:17 AM
Horizontal axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Probe I/O firmware version	1.1.1		
Electrode	T3 Electrode	T3E0769	8/15/2017 9:21:54 AM
E0 calibration	-4.93 mV		11/10/2017 4:58:28 PM
Filling solution	3M KCl	KCL095	11/10/2017 10:16:10 AM
Liquids			
Wash 1	50% IPA:50% Water		11/10/2017 10:14:45 AM
Wash 2	0.5% Triton X-100 in H2O		11/10/2017 10:14:49 AM
Buffer position 1	pH7 Wash		11/10/2017 10:14:51 AM
Buffer position 2	pH 7		11/10/2017 10:14:54 AM
Storage position			11/10/2017 10:15:25 AM
Wash water	9.5e+003 mL	11-10-17	11/10/2017 10:14:37 AM
Waste	7.9e+003 mL		10/13/2017 8:58:05 AM
Temperature controller			8/5/2010 6:35:13 AM
Turbidity detector			3/31/2009 5:24:45 AM
Spectrometer		072390	11/23/2010 11:22:28 AM



## Assay Settings

Sample name: **M16** Experiment start time: **11/10/2017 4:58:04 PM**  
Assay name: **UV-metric pKa** Analyst: **Dorothy Levorse**  
Assay ID: **17K-10009** Instrument ID: **T311053**  
Filename: **C:\Sirius\_T3\Mehtap\20171110\_exp17\_pKa\17K-10009\_M16\_UV-metric pKa.t3r**

## Instrument Settings (continued)

Setting	Value	Batch Id	Install date
Dip probe		11086	
Wavelength coefficient A0	185.563		
Wavelength coefficient A1	2.17439		
Wavelength coefficient A2	-0.000285622		
Total lamp lit time	518:34:28		11/23/2010 11:22:28 AM
Calibrated on	11/8/2017 1:14:37 PM		
Integration time	10		
Scans averaged	10		
Autoloader		T3AL1100237	11/10/2015 9:34:13 AM
Left-right axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Front-back axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Configuration			
Alternate titration position	Titration position		
Alternate reference position	Reference position		
Maximum standard vial volume	3.50 mL		
Maximum alternate vial volume	25.00 mL		
Automatic action idle period	5 minute(s)		
Titrant tube volume	1.3 mL		
Syringe flush count	3.50		
Flowing wash pump volume	20.0 mL		
Flowing wash stir duration	5 s		
Flowing wash stir speed	30%		
Solvent wash stir duration	5 s		
Solvent wash stir speed	30%		
Surfactant wash stir duration	5 s		
Surfactant wash stir speed	30%		
E0 calibration minimum number of points	10		
E0 calibration maximum standard deviation	0.01500		
E0 calibration timeout period	60 s		
E0 calibration stir duration	5 s		
E0 calibration preparation stir speed	30%		
E0 calibration buffer wash stir duration	5 s		
E0 calibration buffer wash stir speed	30%		
E0 calibration reading stir speed	0%		
Spectrometer calibration stir duration	5 s		
Spectrometer calibration stir speed	30%		
Spectrometer calibration wash pump volume	20.0 mL		
Spectrometer calibration wash stir duration	5 s		
Spectrometer calibration wash stir speed	30%		
Overhead dispense height	10000		

## Refinement Settings

Setting	Value	Default value
Turbidity detection method	Spectrometer	Spectrometer
Turbidity wavelength to assess	500.0 nm	500.0 nm
Turbidity maximum absorbance	0.100	0.100
Turbidity probe threshold	50.00	50.00
Exclude turbid points	Yes	Yes
Low intensity warning threshold	100	100
Minimum absorbance change threshold	0.100	0.100
Eigenvector autocorrelation threshold	0.80	0.80
Maximum RMSD severe warning	0.250	0.250
Maximum RMSD warning	0.050	0.050