

Sample name: **M02**
Assay name: **UV-metric psKa**
Assay ID: **17I-22024**
Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22024_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 8:59:52 PM**
Analyst: **Dorothy Levorse**
Instrument ID: **T311053**




Yasuda-Shedlovsky result

| Extrapolation type | pKa 0% | SD | Intercept | Slope | R ² | Ionic strength | Temperature |
|--------------------|--------|-------|-----------|-----------|----------------|----------------|-------------|
| Yasuda-Shedlovsky | 5.02 | ±0.09 | 9.21 | -192.3384 | 0.9849 | 0.165 M | 25.0°C |

Component assay results

| Titration | Methanol weight% | Direction | Result type | Dielectric constant | [H2O] | Ionic strength | Temperature | psKa 1 |
|----------------------------|------------------|-----------|---------------|---------------------|--------|----------------|-------------|--------|
| 17I-22024 Points 4 to 36 | 49.54 % | Up | UV-metric pKa | 56.6 | 24.7 M | 0.157 M | 25.0°C | ✓ 4.43 |
| 17I-22024 Points 38 to 78 | 40.02 % | Up | UV-metric pKa | 61.0 | 30.0 M | 0.165 M | 25.0°C | ✓ 4.55 |
| 17I-22024 Points 80 to 123 | 30.16 % | Up | UV-metric pKa | 65.4 | 35.8 M | 0.172 M | 25.0°C | ✓ 4.73 |

Graphs

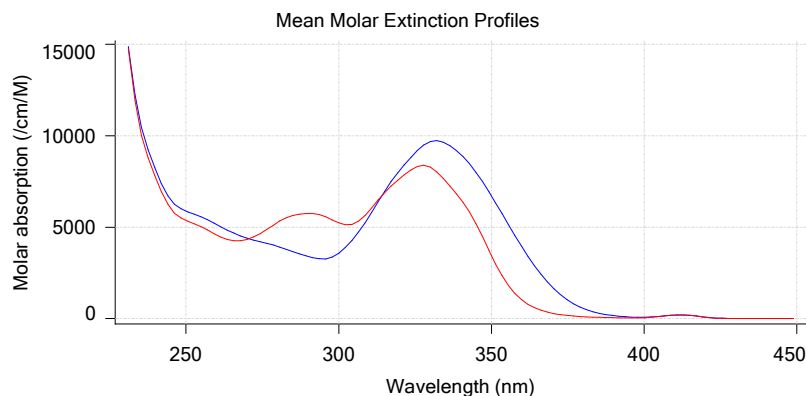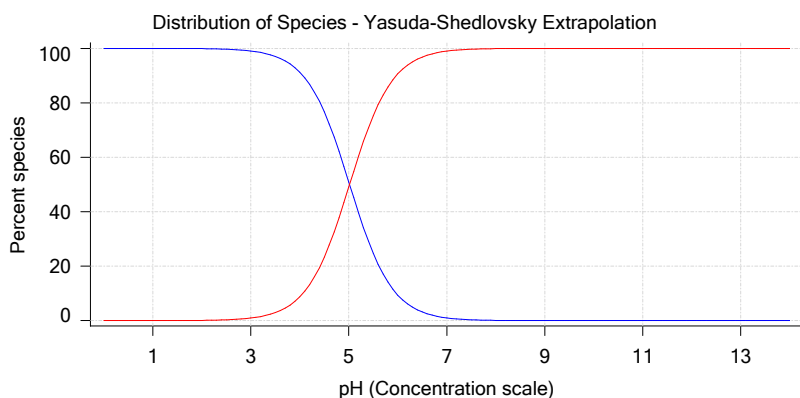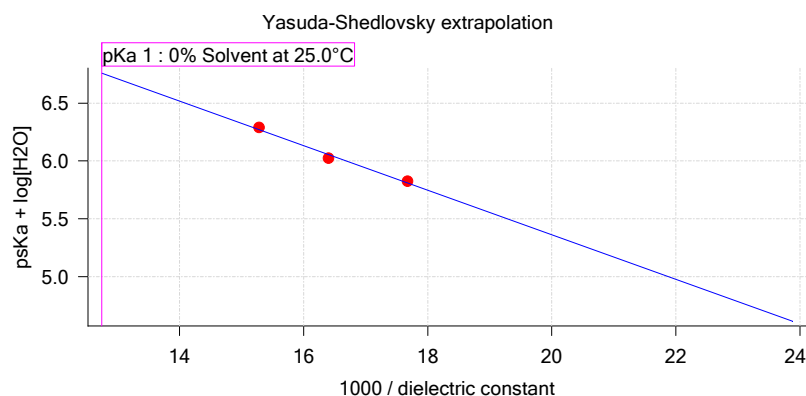

UV-metric psKa Titration 1 of 3 17I-22024 Points 4 to 36

Results

| | |
|-------------------------------|-----------------------------|
| pKa 1 | 4.43 |
| RMSD | 0.023 0.024 |
| Chi squared | 0.1028 |
| PCA calculated number of pKas | 2 |
| Average ionic strength | 0.157 M |
| Average temperature | 25.0°C |
| Analyte concentration range | 58.2 µM to 54.8 µM |
| Methanol weight % | 49.5 % |
| Dielectric constant | 56.6 |
| Water concentration | 24.7 M |
| Number of pKas source | Predicted |
| Wavelength clipping | 230.0 nm to 450.0 nm |

Sample name: **M02**
 Assay name: **UV-metric psKa**
 Assay ID: **17I-22024**
 Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22024_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 8:59:52 PM**
 Analyst: **Dorothy Leverse**
 Instrument ID: **T311053**

Results (continued)

pH clipping 1.475 to 12.517


Warnings and errors

Errors: None
 Warnings: PCA calculation disagrees with predicted number of pKas

Assay Settings

| Setting | Value | Original Value | Date/Time changed | Imported from |
|-----------------------------|------------------|----------------|-------------------|---------------|
| Buffer in use | Yes | | | |
| Buffer type | Phosphate Buffer | | | |
| Assay Medium | | | | |
| Volume of buffer introduced | 0.025000 mL | | | |
| Add buffer manually | Manual | | | |


Graphs

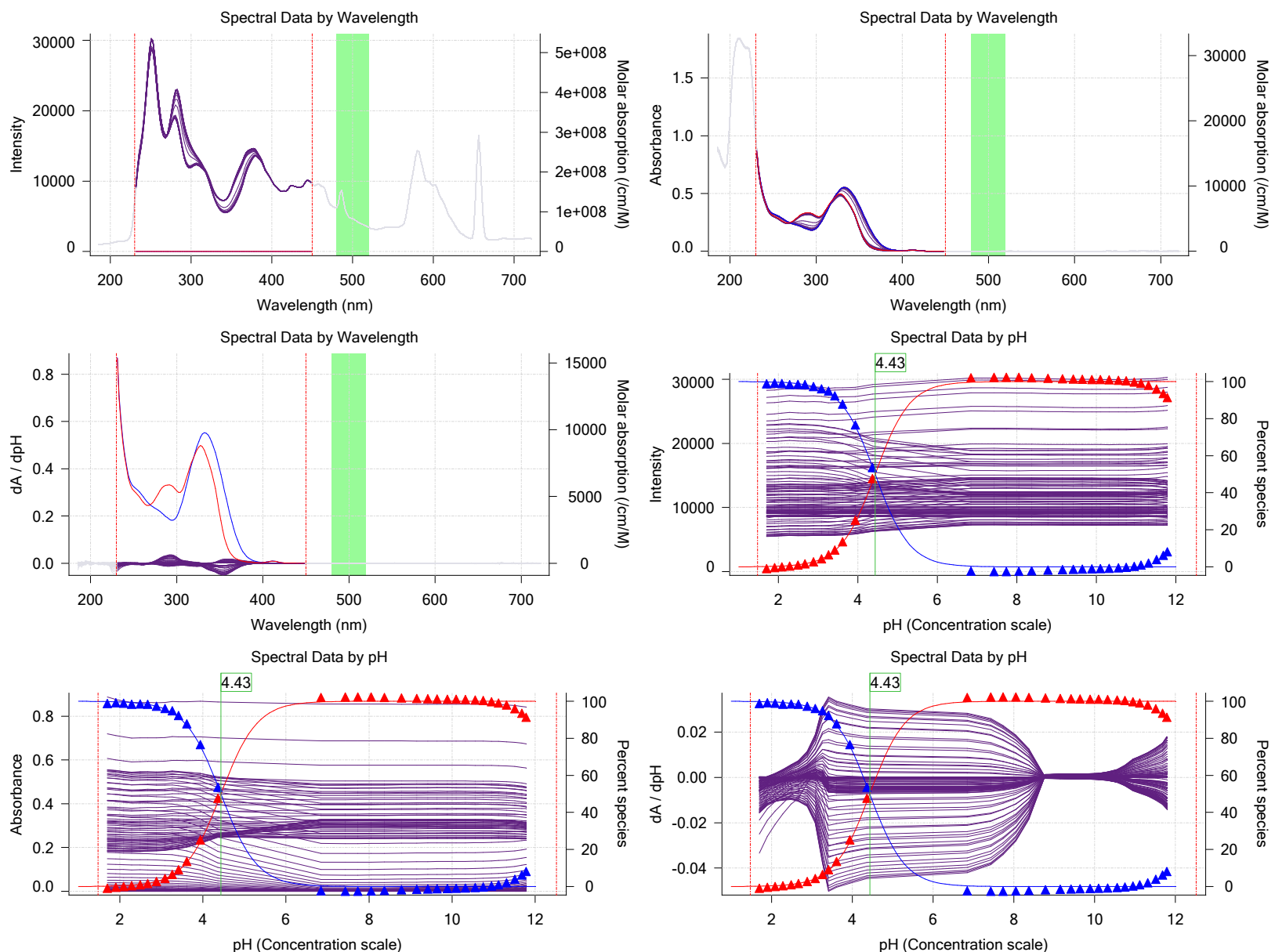

Sample name: **M02**
Assay name: **UV-metric psKa**
Assay ID: **171-22024**
Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\171-22024_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 8:59:52 PM**
Analyst: **Dorothy Leverse**
Instrument ID: **T311053**

Graphs (continued)

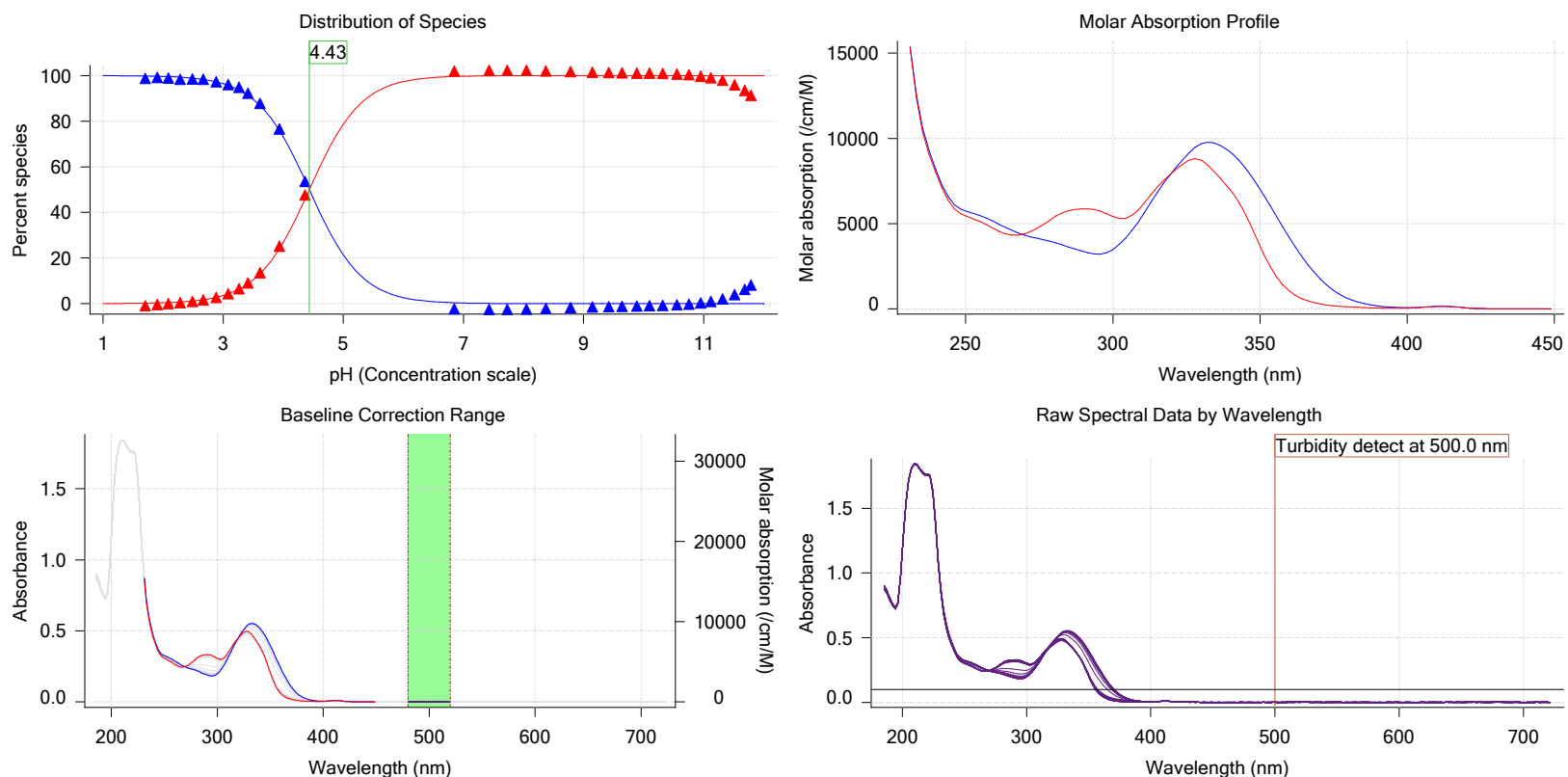

UV-metric psKa Titration 2 of 3 171-22024 Points 38 to 78

Results

| | |
|-------------------------------|-----------------------------|
| pKa 1 | 4.55 |
| RMSD | 0.027 0.029 |
| Chi squared | 0.1091 |
| PCA calculated number of pKas | 4 |
| Average ionic strength | 0.165 M |
| Average temperature | 25.0°C |
| Analyte concentration range | 47.8 µM to 45.2 µM |
| Methanol weight % | 40.0 % |
| Dielectric constant | 61.0 |
| Water concentration | 30.0 M |
| Number of pKas source | Predicted |
| Wavelength clipping | 230.0 nm to 450.0 nm |
| pH clipping | 1.495 to 12.529 |

Warnings and errors

Errors: None
Warnings: PCA calculation disagrees with predicted number of pKas

Assay Settings

| Setting | Value | Original Value | Date/Time changed | Imported from |
|---------------|------------------|----------------|-------------------|---------------|
| Buffer in use | Yes | | | |
| Buffer type | Phosphate Buffer | | | |

Assay Medium


Sample name: **M02**
 Assay name: **UV-metric psKa**
 Assay ID: **17I-22024**
 Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22024_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 8:59:52 PM**
 Analyst: **Dorothy Leverse**
 Instrument ID: **T311053**

Assay Settings (continued)

| Setting | Value | Original Value | Date/Time changed | Imported from |
|-----------------------------|-------------|----------------|-------------------|---------------|
| Volume of buffer introduced | 0.025000 mL | | | |
| Add buffer manually | Manual | | | |


Graphs

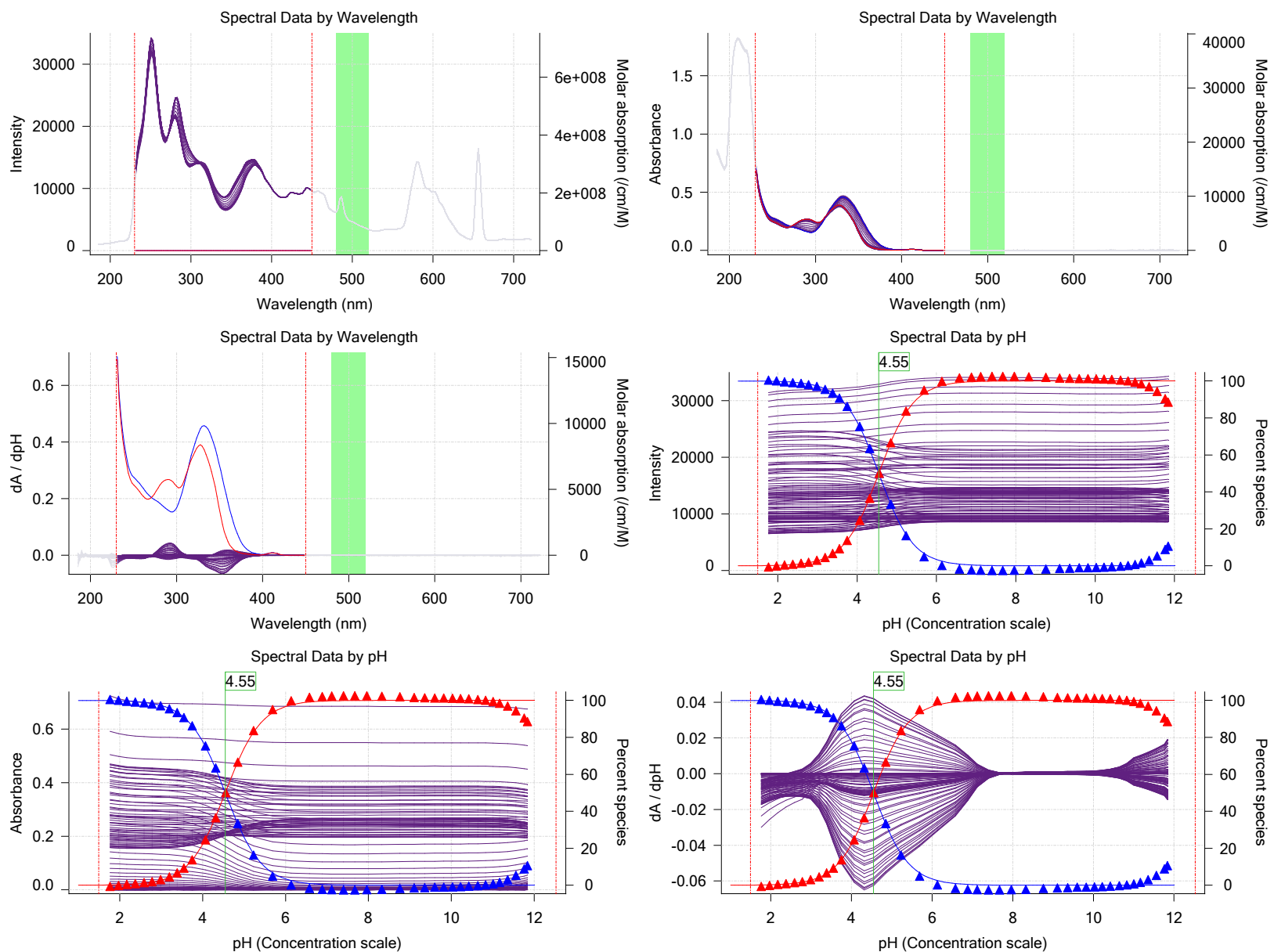

Sample name: **M02**
Assay name: **UV-metric psKa**
Assay ID: **171-22024**
Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\171-22024_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 8:59:52 PM**
Analyst: **Dorothy Leverse**
Instrument ID: **T311053**

Graphs (continued)

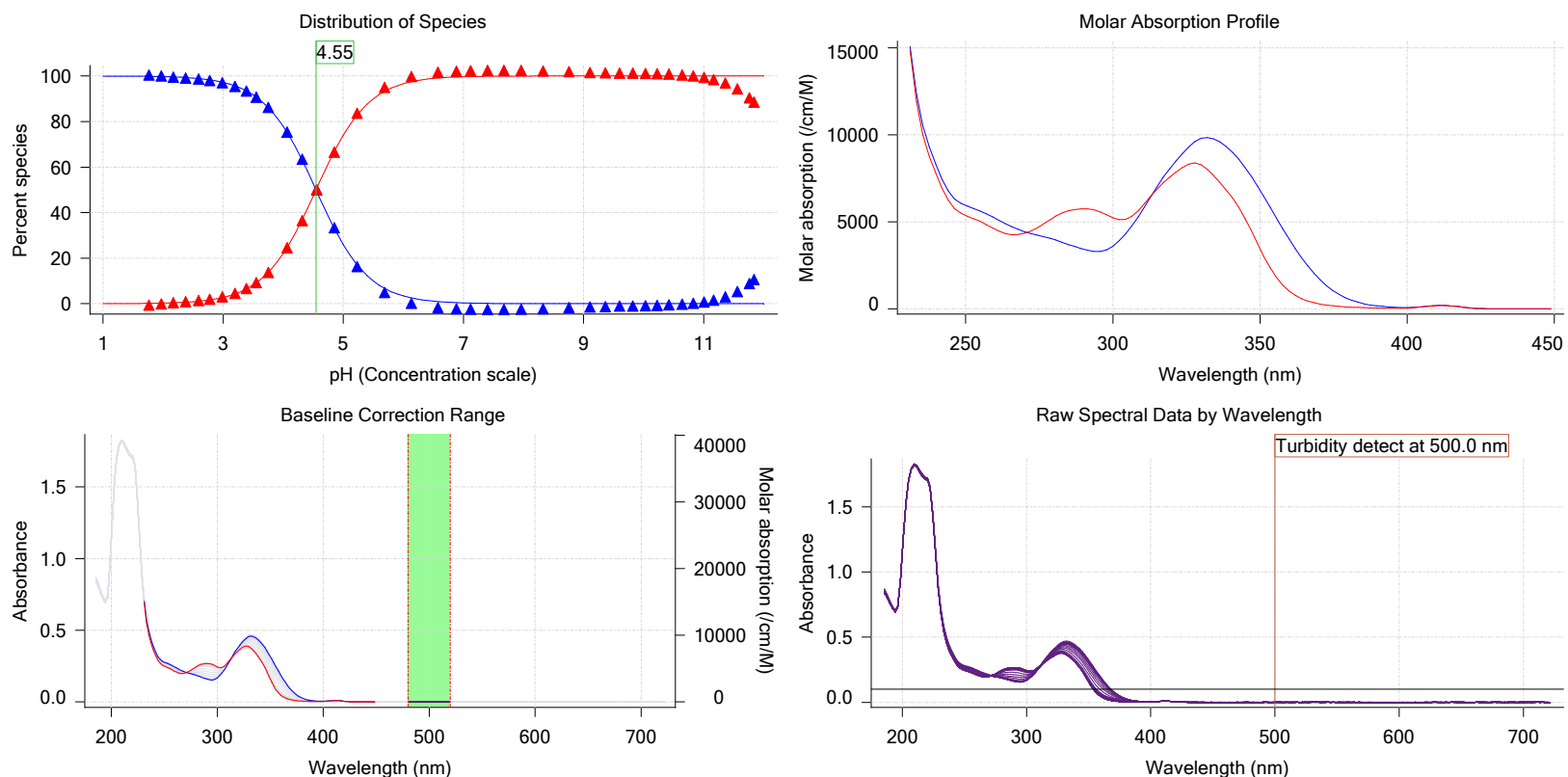

UV-metric psKa Titration 3 of 3 171-22024 Points 80 to 123

Results

| | |
|-------------------------------|-----------------------------|
| pKa 1 | 4.73 |
| RMSD | 0.032 0.033 |
| Chi squared | 0.1231 |
| PCA calculated number of pKas | 3 |
| Average ionic strength | 0.172 M |
| Average temperature | 25.0°C |
| Analyte concentration range | 36.8 µM to 34.8 µM |
| Methanol weight % | 30.2 % |
| Dielectric constant | 65.4 |
| Water concentration | 35.8 M |
| Number of pKas source | Predicted |
| Wavelength clipping | 230.0 nm to 450.0 nm |
| pH clipping | 1.500 to 12.512 |

Warnings and errors

Errors: None
Warnings: PCA calculation disagrees with predicted number of pKas

Assay Settings

| Setting | Value | Original Value | Date/Time changed | Imported from |
|---------------|------------------|----------------|-------------------|---------------|
| Buffer in use | Yes | | | |
| Buffer type | Phosphate Buffer | | | |

Assay Medium


Sample name: **M02**
 Assay name: **UV-metric psKa**
 Assay ID: **171-22024**
 Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\171-22024_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 8:59:52 PM**
 Analyst: **Dorothy Levorse**
 Instrument ID: **T311053**

Assay Settings (continued)

| Setting | Value | Original Value | Date/Time changed | Imported from |
|-----------------------------|-------------|----------------|-------------------|---------------|
| Volume of buffer introduced | 0.025000 mL | | | |
| Add buffer manually | Manual | | | |


Graphs

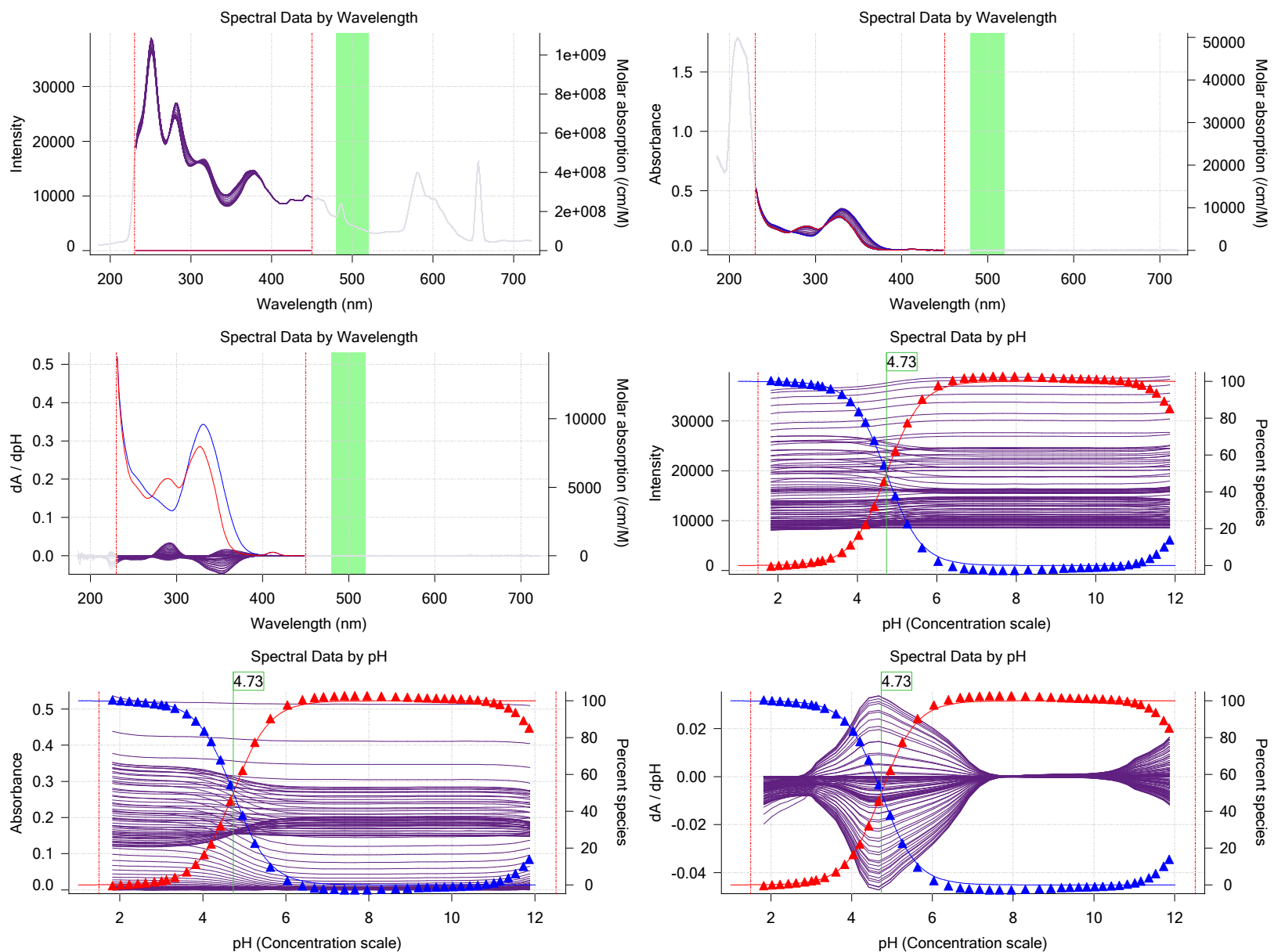

Sample name: **M02**
Assay name: **UV-metric psKa**
Assay ID: **171-22024**
Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\171-22024_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 8:59:52 PM**
Analyst: **Dorothy Leverse**
Instrument ID: **T311053**

Graphs (continued)

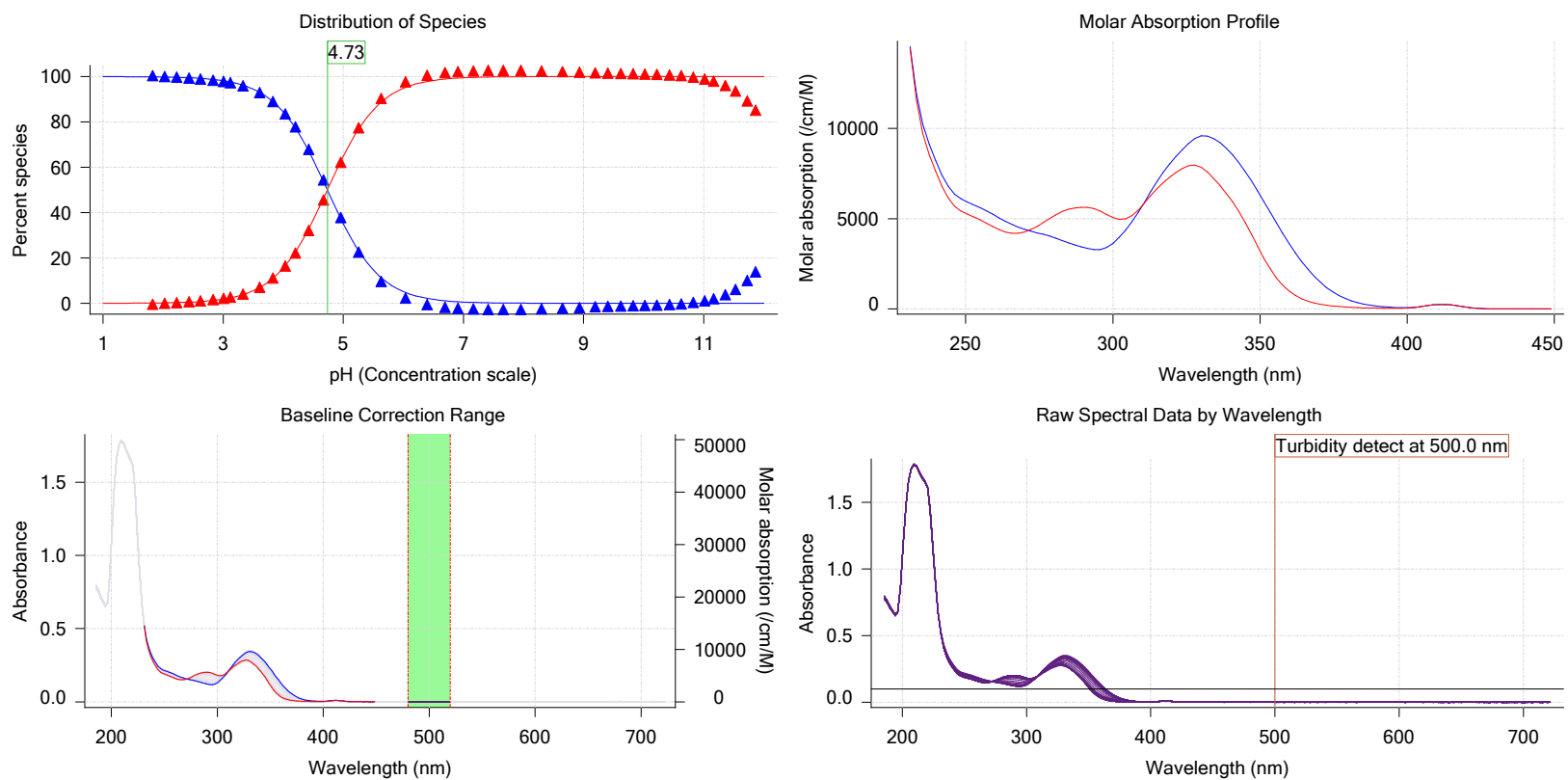

Assay Model

| Settings | Value | Date/Time changed | Imported from |
|-----------------------------------|------------|-----------------------|--------------------|
| Sample name | M02 | 9/22/2017 10:27:56 AM | User entered value |
| Sample by | Volume | | Default value |
| Sample volume | 0.0020 mL | 9/22/2017 10:27:56 AM | User entered value |
| Solvent | DMSO | | Default value |
| Sample concentration | 0.046400 M | 9/22/2017 10:27:56 AM | User entered value |
| Solubility | Unknown | | Default value |
| Molecular weight | 289.26 | 9/22/2017 10:28:05 AM | User entered value |
| Individual pKa ionic environments | No | | Default value |
| Number of pKas | 1 | 9/22/2017 10:27:56 AM | User entered value |
| Sample is a | Base | 9/22/2017 10:27:56 AM | User entered value |
| pKa 1 | 5.60 | 9/22/2017 10:27:56 AM | User entered value |
| logp (XH +) | -10.00 | | Default value |
| logP (neutral X) | -10.00 | 9/22/2017 10:27:56 AM | User entered value |

Assay Settings

| Setting | Value | Original Value | Date/Time changed | Imported from |
|-------------------------------------|-----------------|----------------|-------------------|---------------|
| General Settings | | | | |
| Analyst name | Dorothy Leverse | | | |
| Separate reference vial | Yes | | | |
| Standard Experiment Settings | | | | |
| Number of titrations | 3 | | | |
| Minimum pH | 2.000 | | | |
| Maximum pH | 12.000 | | | |
| pH step between points of | 0.200 | | | |
| Minimum titrant addition | 0.00002 mL | | | |
| Maximum titrant addition | 0.10000 mL | | | |
| Argon flow rate | 100% | | | |

Sample name: **M02**
 Assay name: **UV-metric psKa**
 Assay ID: **17I-22024**
 Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22024_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 8:59:52 PM**
 Analyst: **Dorothy Levorse**
 Instrument ID: **T311053**

Assay Settings (continued)

| Setting | Value | Original Value | Date/Time changed | Imported from |
|------------------------------------|--------------------|----------------|-------------------|---------------|
| Start titration using | Cautious pH adjust | | | |
| Advanced General Settings | | | | |
| Detect turbidity using | Spectrometer | | | |
| Monitor at a wavelength of | 500.0 nm | | | |
| Absorbance threshold of | 0.100 | | | |
| Collect turbidity sensor data | No | | | |
| Stir after titrant addition for | 5 seconds | | | |
| For titrant addition, stir at | 15% | | | |
| Titration Pre-Dose | | | | |
| Titration pre-dose | None | | | |
| Assay Medium | | | | |
| Cosolvent in use | Yes | | | |
| Cosolvent type | Methanol | | | |
| Cosolvent volume | 1.15 mL | | | |
| Cosolvent added | Automatic | | | |
| ISA water volume | 0.35 mL | | | |
| Water added | Automatic | | | |
| After water addition, stir for | 5 seconds | | | |
| At a speed of | 15% | | | |
| Buffer in use | Yes | | | |
| Buffer type | Phosphate Buffer | | | |
| Volume of buffer introduced | 0.025000 mL | | | |
| Add buffer manually | Manual | | | |
| After medium addition, stir for | 5 seconds | | | |
| Sample Sonication | | | | |
| Sonicate | No | | | |
| Sample Dissolution | | | | |
| Perform a dissolution stage | No | | | |
| Carbonate purge | | | | |
| Perform a carbonate purge | No | | | |
| Temperature Control | | | | |
| Wait for temperature | Yes | | | |
| Required start temperature | 25.0°C | | | |
| Acceptable deviation | 0.5°C | | | |
| Time to wait | 60 seconds | | | |
| Stir speed of | 15% | | | |
| Titration 1 | | | | |
| Titrate from | Low to high pH | | | |
| Adjust to start pH | Yes | | | |
| After pH adjust stir for | 10 seconds | | | |
| Titration 2 | | | | |
| Titrate from | Low to high pH | | | |
| Additional cosolvent volume | 0.00 mL | | | |
| Add additional water | 0.15 mL | | | |
| Additional water added | Automatic | | | |
| After pH adjust stir for | 10 seconds | | | |
| Titration 3 | | | | |
| Titrate from | Low to high pH | | | |
| Additional cosolvent volume | 0.00 mL | | | |
| Add additional water | 0.34 mL | | | |
| Additional water added | Automatic | | | |
| After pH adjust stir for | 10 seconds | | | |
| Data Point Stability | | | | |
| Stir during data point collection | Yes | | | |
| For point collection, stir at | 15% | | | |
| Delay before data point collection | 0 seconds | | | |
| Number of points to average | 20 points | | | |
| Time interval between points | 0.50 seconds | | | |

Sample name: **M02** Experiment start time: **9/22/2017 8:59:52 PM**
 Assay name: **UV-metric psKa** Analyst: **Dorothy Levorse**
 Assay ID: **17I-22024** Instrument ID: **T311053**
 Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22024_M02_UV-metric psKa.t3r**

Assay Settings (continued)

| Setting | Value | Original Value | Date/Time changed | Imported from |
|-------------------------------------|----------------|----------------|-------------------|---------------|
| Required maximum standard deviation | 0.00500 dpH/dt | | | |
| Stability timeout after | 60 seconds | | | |
| Experiment cleanup | | | | |
| Adjust pH to cleanup | To start pH | | | |
| And then stir for | 60 seconds | | | |
| For cleaning, stir at | 20% | | | |
| Then add water volume | 0.25 mL | | | |
| And then stir for | 30 seconds | | | |

Calibration Settings

| Setting | Value | Date/Time changed | Imported from |
|---------------------------|--------|----------------------|---|
| Four-Plus alpha | 0.105 | 9/22/2017 8:59:52 PM | C:\Sirius_T3\17I-22019_Blank standardisation.t3r |
| Four-Plus S | 1.0019 | 9/22/2017 8:59:52 PM | C:\Sirius_T3\17I-22019_Blank standardisation.t3r |
| Four-Plus jH | 0.3 | 9/22/2017 8:59:52 PM | C:\Sirius_T3\17I-22019_Blank standardisation.t3r |
| Four-Plus jOH | -1.0 | 9/22/2017 8:59:52 PM | C:\Sirius_T3\17I-22019_Blank standardisation.t3r |
| Base concentration factor | 1.000 | 9/22/2017 8:59:52 PM | C:\Sirius_T3\17I-22016_KHP_Base standardisation using KHP.t3r |
| Acid concentration factor | 0.995 | 9/22/2017 8:59:52 PM | C:\Sirius_T3\17I-22019_Blank standardisation.t3r |

Instrument Settings

| Setting | Value | Batch Id | Install date |
|----------------------------------|----------------------------|-------------|------------------------|
| Instrument owner | Merck | | |
| Instrument ID | T311053 | | |
| Instrument type | T3 Simulator | | |
| Software version | 1.1.3.0 | | |
| Dispenser module | | T3DM1100253 | 3/31/2009 6:24:52 AM |
| Dispenser 0 | Water | | 3/31/2009 6:25:05 AM |
| Syringe volume | 2.5 mL | | |
| Firmware version | 1.2.1(r2) | | |
| Titrant | Water (0.15 M KCl) | 8-18-17 | 9/18/2017 9:13:04 AM |
| Dispenser 2 | Acid | | 3/31/2009 6:25:11 AM |
| Syringe volume | 0.5 mL | | |
| Firmware version | 1.2.1(r2) | | |
| Titrant | Acid (0.5 M HCl) | 166940 | 9/8/2017 9:21:27 AM |
| Dispenser 1 | Base | | 3/31/2009 6:25:21 AM |
| Syringe volume | 0.5 mL | | |
| Firmware version | 1.2.1(r2) | | |
| Titrant | Base (0.5 M KOH) | 9-22-17 | 9/22/2017 4:02:42 PM |
| Dispenser 5 | Cosolvent | | 3/31/2009 6:26:24 AM |
| Syringe volume | 2.5 mL | | |
| Firmware version | 1.2.1(r2) | | |
| Distribution valve 5 | Distribution Valve | | 3/31/2009 6:28:19 AM |
| Firmware version | 1.1.3 | | |
| Port A | Methanol (80%, 0.15 M KCl) | 8-15-17 | 9/20/2017 4:38:16 PM |
| Port B | Cyclohexane | | 9/19/2017 2:15:02 PM |
| Dispenser 3 | Buffer | | 8/3/2010 6:05:16 AM |
| Syringe volume | 0.5 mL | | |
| Firmware version | 1.2.1(r2) | | |
| Titrant | Phosphate Buffer | | 9/12/2017 12:32:29 PM |
| Dispenser 6 | Octanol | | 10/22/2010 11:52:43 AM |
| Syringe volume | 0.5 mL | | |
| Firmware version | 1.2.1(r2) | | |
| Titrant | Octanol | 9-14-17 | 9/14/2017 10:30:38 AM |
| Titrator | | T3TM1100153 | 3/31/2009 6:24:17 AM |
| Horizontal axis firmware version | 1.17 AI1DI2DO2 Stepper 2 | | |
| Vertical axis firmware version | 1.17 AI1DI2DO2 Stepper 2 | | |
| Chassis I/O firmware version | 1.11 AI1DI0DO4 Norgren I/O | | |



Assay Settings

Sample name: **M02** Experiment start time: **9/22/2017 8:59:52 PM**
Assay name: **UV-metric psKa** Analyst: **Dorothy Levorse**
Assay ID: **171-22024** Instrument ID: **T311053**
Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\171-22024_M02_UV-metric psKa.t3r**

Instrument Settings (continued)

| Setting | Value | Batch Id | Install date |
|---|----------------------------|-------------|------------------------|
| Probe I/O firmware version | 1.1.1 | | |
| Electrode | T3 Electrode | T3E0769 | 8/15/2017 10:21:54 AM |
| E0 calibration | -8.79 mV | | 9/22/2017 9:00:16 PM |
| Filling solution | 3M KCl | KCL095 | 9/18/2017 9:17:15 AM |
| Liquids | | | |
| Wash 1 | 50% IPA:50% Water | | 9/22/2017 10:05:34 AM |
| Wash 2 | 0.5% Triton X-100 in H2O | | 9/22/2017 10:05:36 AM |
| Buffer position 1 | pH7 Wash | | 9/22/2017 10:05:39 AM |
| Buffer position 2 | pH 7 | | 9/22/2017 10:05:42 AM |
| Storage position | | | 9/22/2017 10:06:32 AM |
| Wash water | 4.6e+002 mL | 9-18-17 | 9/18/2017 8:54:32 AM |
| Waste | 9.6e+003 mL | | 9/18/2017 8:54:39 AM |
| Temperature controller | | | 8/5/2010 7:35:13 AM |
| Turbidity detector | | | 3/31/2009 6:24:45 AM |
| Spectrometer | | 072390 | 11/23/2010 12:22:28 PM |
| Dip probe | | 11086 | |
| Wavelength coefficient A0 | 185.563 | | |
| Wavelength coefficient A1 | 2.17439 | | |
| Wavelength coefficient A2 | -0.000285622 | | |
| Total lamp lit time | 207:39:18 | | 11/23/2010 12:22:28 PM |
| Calibrated on | 9/18/2017 9:35:14 AM | | |
| Integration time | 11 | | |
| Scans averaged | 10 | | |
| Autoloader | | T3AL1100237 | 11/10/2015 10:34:13 AM |
| Left-right axis firmware version | 1.17 AI1DI2DO2 Stepper 2 | | |
| Front-back axis firmware version | 1.17 AI1DI2DO2 Stepper 2 | | |
| Vertical axis firmware version | 1.17 AI1DI2DO2 Stepper 2 | | |
| Chassis I/O firmware version | 1.11 AI1DI0DO4 Norgren I/O | | |
| Configuration | | | |
| Alternate titration position | Titration position | | |
| Alternate reference position | Reference position | | |
| Maximum standard vial volume | 3.50 mL | | |
| Maximum alternate vial volume | 25.00 mL | | |
| Automatic action idle period | 5 minute(s) | | |
| Titrant tube volume | 1.3 mL | | |
| Syringe flush count | 3.50 | | |
| Flowing wash pump volume | 20.0 mL | | |
| Flowing wash stir duration | 5 s | | |
| Flowing wash stir speed | 30% | | |
| Solvent wash stir duration | 5 s | | |
| Solvent wash stir speed | 30% | | |
| Surfactant wash stir duration | 5 s | | |
| Surfactant wash stir speed | 30% | | |
| E0 calibration minimum number of points | 10 | | |
| E0 calibration maximum standard deviation | 0.01500 | | |
| E0 calibration timeout period | 60 s | | |
| E0 calibration stir duration | 5 s | | |
| E0 calibration preparation stir speed | 30% | | |
| E0 calibration buffer wash stir duration | 5 s | | |
| E0 calibration buffer wash stir speed | 30% | | |
| E0 calibration reading stir speed | 0% | | |
| Spectrometer calibration stir duration | 5 s | | |
| Spectrometer calibration stir speed | 30% | | |
| Spectrometer calibration wash pump volume | 20.0 mL | | |
| Spectrometer calibration wash stir duration | 5 s | | |
| Spectrometer calibration wash stir speed | 30% | | |
| Overhead dispense height | 10000 | | |



Assay Settings

Sample name: **M02**
Assay name: **UV-metric psKa**
Assay ID: **17I-22024**
Filename: **C:\Sirius_T3\Mehtap\20170922_exp06_M02_M14-M16_D01-D03\17I-22024_M02_UV-metric psKa.t3r**

Experiment start time: **9/22/2017 8:59:52 PM**
Analyst: **Dorothy Levorse**
Instrument ID: **T311053**

Refinement Settings

| Setting | Value | Default value |
|---------------------------------------|--------------|---------------|
| Turbidity detection method | Spectrometer | Spectrometer |
| Turbidity wavelength to assess | 500.0 nm | 500.0 nm |
| Turbidity maximum absorbance | 0.100 | 0.100 |
| Turbidity probe threshold | 50.00 | 50.00 |
| Exclude turbid points | Yes | Yes |
| Low intensity warning threshold | 100 | 100 |
| Minimum absorbance change threshold | 0.100 | 0.100 |
| Eigenvector autocorrelation threshold | 0.80 | 0.80 |
| Maximum RMSD severe warning | 0.250 | 0.250 |
| Maximum RMSD warning | 0.050 | 0.050 |

Tray Information

Title
Location B5