

Sample name: **M05**
Assay name: **UV-metric psKa**
Assay ID: **17I-19008**
Filename: **C:\Sirius_T3\Mehtap\20170918_exp04_uv_M01-M14\17I-19008_M05_UV-metric psKa.t3r**

Experiment start time: **9/19/2017 6:09:04 AM**
Analyst: **Dorothy Levorse**
Instrument ID: **T311053**




Yasuda-Shedlovsky result

Extrapolation type	pKa 0%	SD	Intercept	Slope	R ²	Ionic strength	Temperature
Yasuda-Shedlovsky	4.60	±0.11	11.50	-404.5620	0.9952	0.166 M	24.9°C

Component assay results

Titration	Methanol weight%	Direction	Result type	Dielectric constant	[H2O]	Ionic strength	Temperature	psKa 1
17I-19008 Points 4 to 42	49.21 %	Up	UV-metric pKa	56.8	24.9 M	0.157 M	24.9°C	✓ 2.99
17I-19008 Points 44 to 86	39.44 %	Up	UV-metric pKa	61.3	30.3 M	0.167 M	24.9°C	✓ 3.37
17I-19008 Points 88 to 130	29.71 %	Up	UV-metric pKa	65.7	36.0 M	0.173 M	24.9°C	✓ 3.80

Graphs

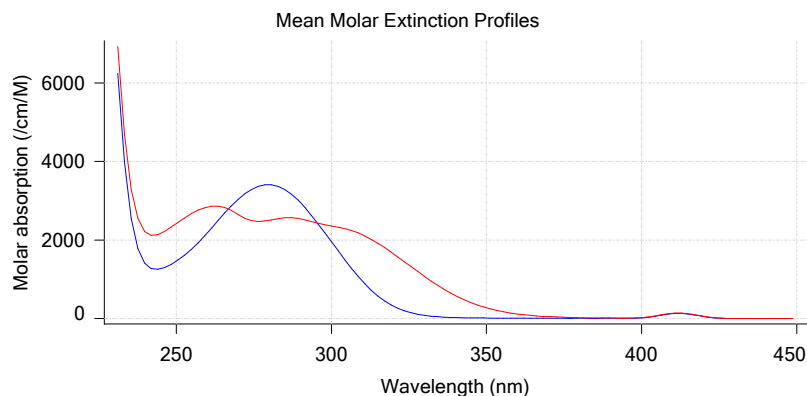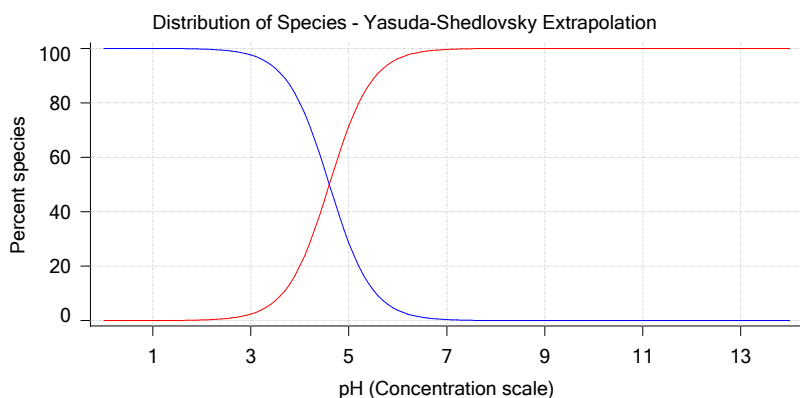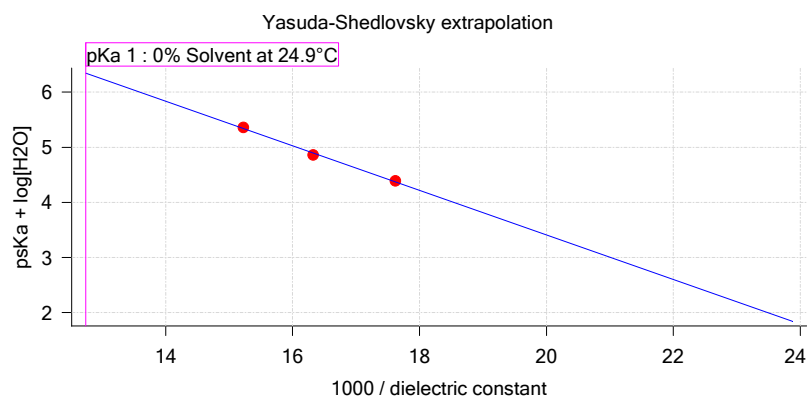

UV-metric psKa Titration 1 of 3 17I-19008 Points 4 to 42

Results

pKa 1	2.99
RMSD	0.001 0.003
Chi squared	0.0046
PCA calculated number of pKas	1
Average ionic strength	0.157 M
Average temperature	24.9°C
Analyte concentration range	88.5 µM to 83.0 µM
Methanol weight %	49.2 %
Dielectric constant	56.8
Water concentration	24.9 M
Number of pKas source	Manual (1)
Wavelength clipping	230.0 nm to 450.0 nm

Sample name: **M05**
 Assay name: **UV-metric psKa**
 Assay ID: **171-19008**
 Filename: **C:\Sirius_T3\Mehtap\20170918_exp04_uv_M01-M14\171-19008_M05_UV-metric psKa.t3r**

Experiment start time: **9/19/2017 6:09:04 AM**
 Analyst: **Dorothy Leverse**
 Instrument ID: **T311053**

Results (continued)

pH clipping 1.478 to 12.516


Warnings and errors

Errors None
 Warnings PCA calculation disagrees with predicted number of pKas

Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
Assay Medium				
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			


Graphs

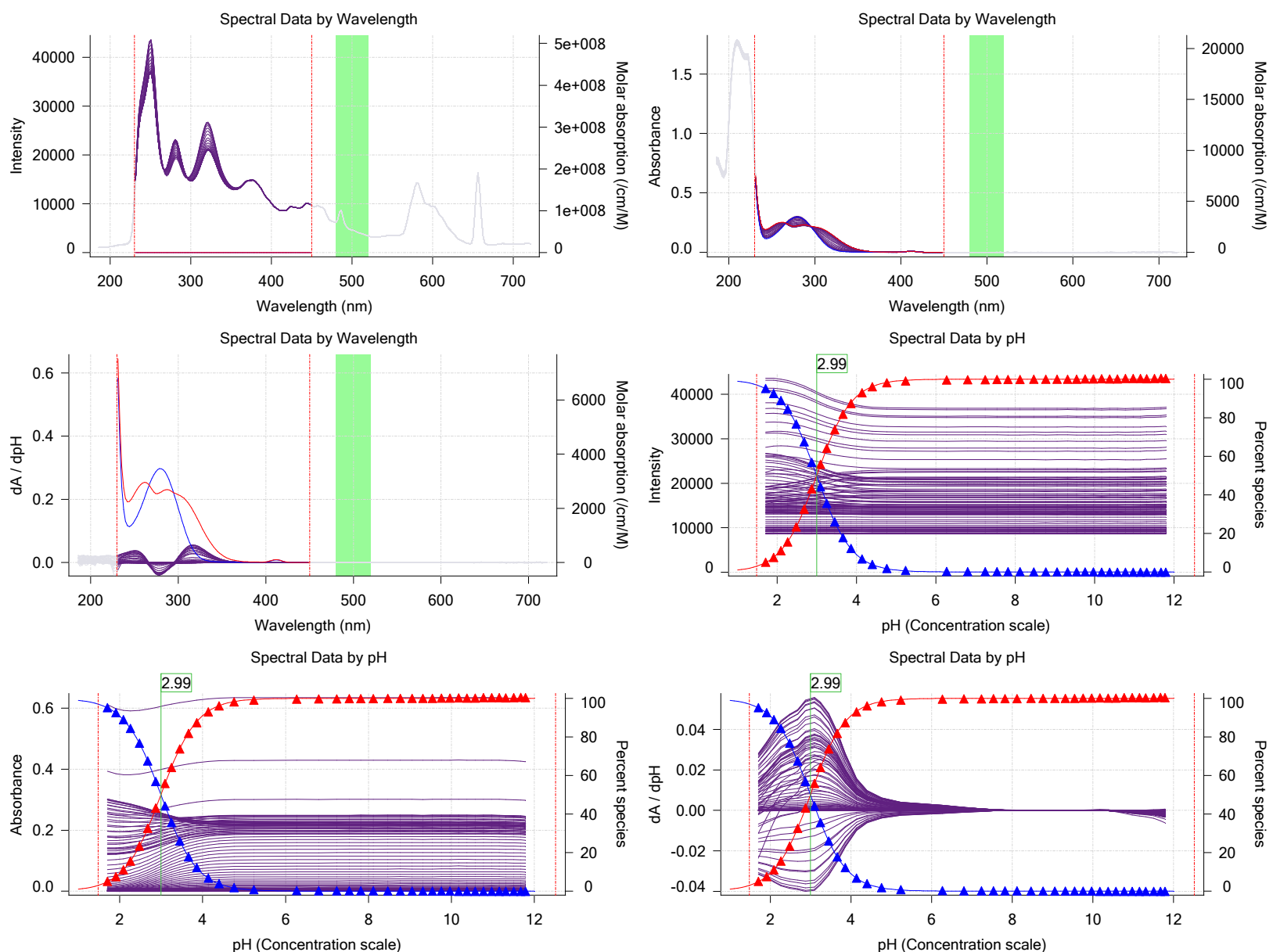

Sample name: **M05**
Assay name: **UV-metric psKa**
Assay ID: **171-19008**
Filename: **C:\Sirius_T3\Mehtap\20170918_exp04_uv_M01-M14\171-19008_M05_UV-metric psKa.t3r**

Experiment start time: **9/19/2017 6:09:04 AM**
Analyst: **Dorothy Levorse**
Instrument ID: **T311053**

Graphs (continued)

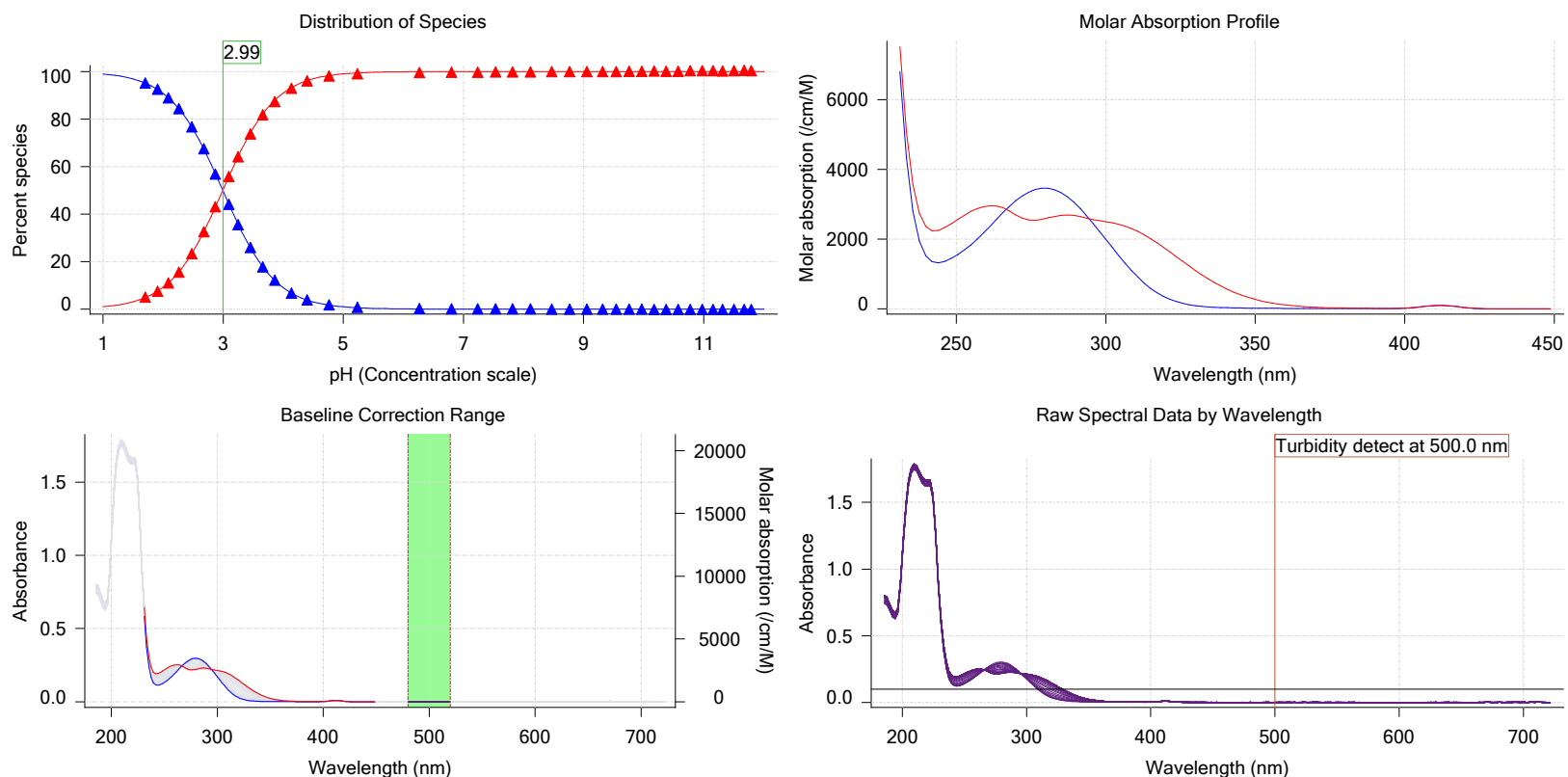

UV-metric psKa Titration 2 of 3 171-19008 Points 44 to 86

Results

pKa 1	3.37
RMSD	0.001 0.002
Chi squared	0.0021
PCA calculated number of pKas	2
Average ionic strength	0.167 M
Average temperature	24.9°C
Analyte concentration range	72.1 µM to 67.9 µM
Methanol weight %	39.4 %
Dielectric constant	61.3
Water concentration	30.3 M
Number of pKas source	Manual (1)
Wavelength clipping	230.0 nm to 450.0 nm
pH clipping	1.438 to 12.520

Warnings and errors

Errors None
Warnings None

Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			

Assay Medium


Sample name: **M05**
 Assay name: **UV-metric psKa**
 Assay ID: **171-19008**
 Filename: **C:\Sirius_T3\Mehtap\20170918_exp04_uv_M01-M14\171-19008_M05_UV-metric psKa.t3r**

Experiment start time: **9/19/2017 6:09:04 AM**
 Analyst: **Dorothy Leverse**
 Instrument ID: **T311053**

Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			


Graphs

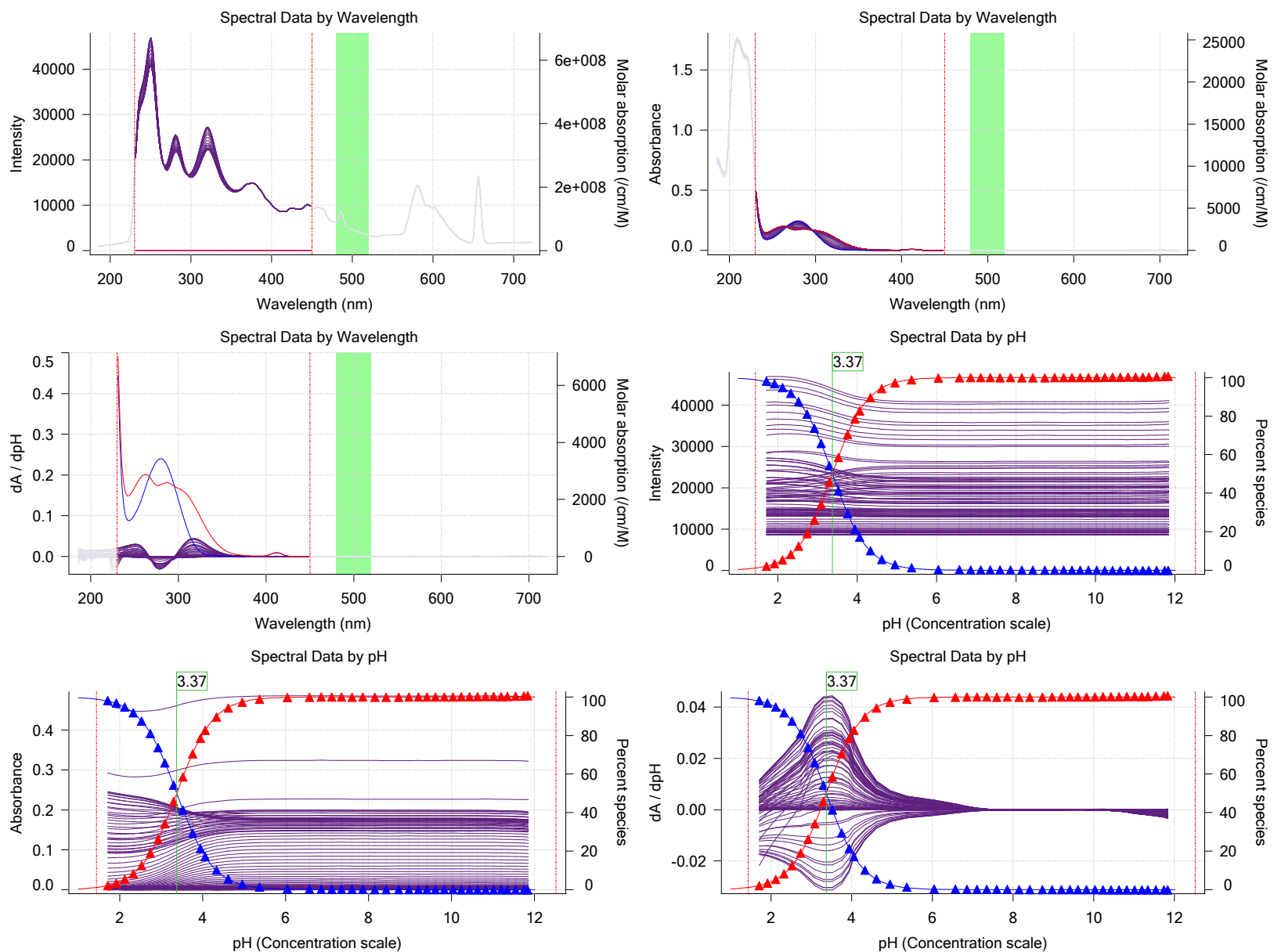

Sample name: **M05**
Assay name: **UV-metric psKa**
Assay ID: **171-19008**
Filename: **C:\Sirius_T3\Mehtap\20170918_exp04_uv_M01-M14\171-19008_M05_UV-metric psKa.t3r**

Experiment start time: **9/19/2017 6:09:04 AM**
Analyst: **Dorothy Levorse**
Instrument ID: **T311053**

Graphs (continued)

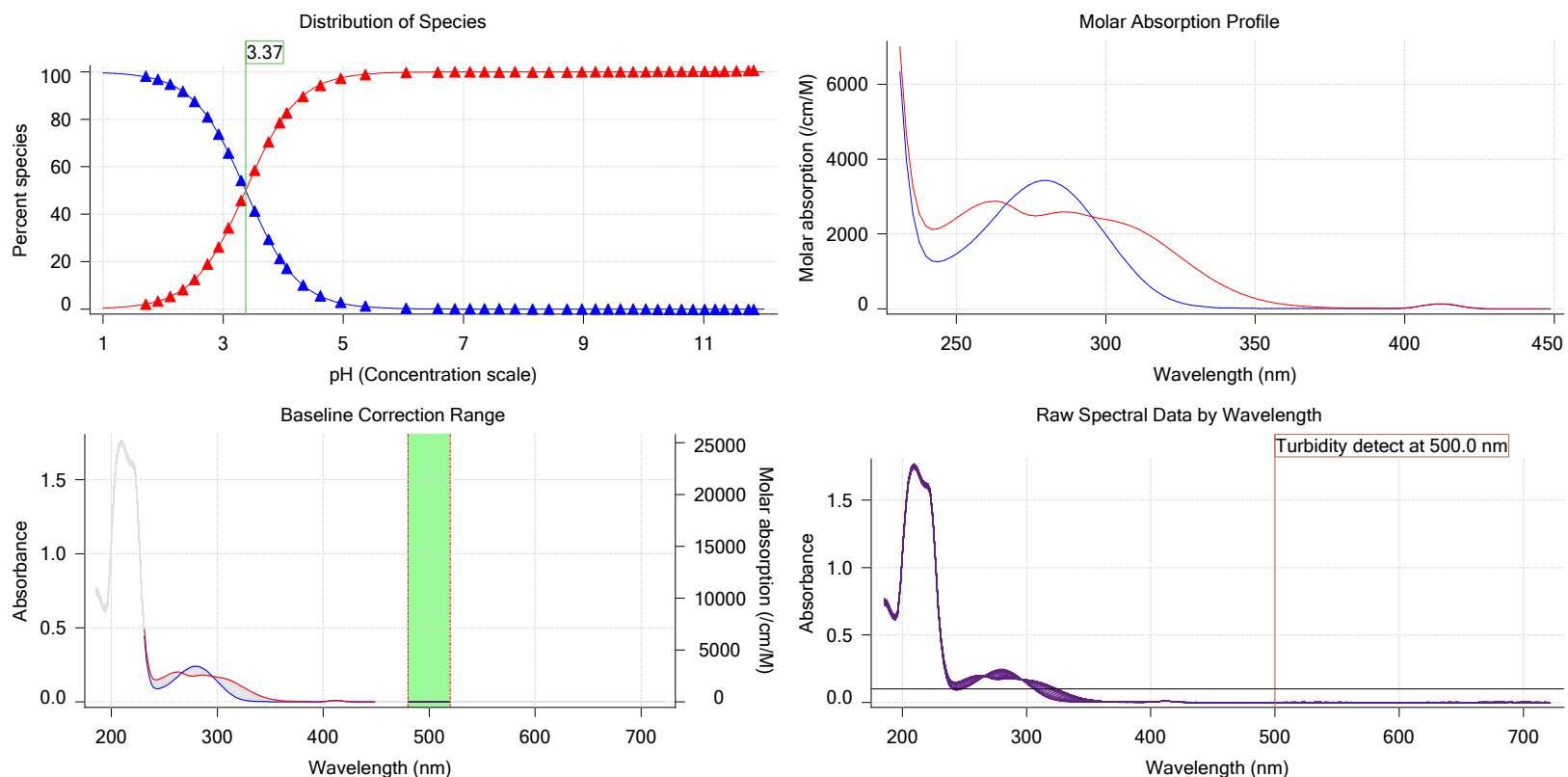

UV-metric psKa Titration 3 of 3 171-19008 Points 88 to 130

Results

pKa 1	3.80
RMSD	0.002 0.003
Chi squared	0.0026
PCA calculated number of pKas	2
Average ionic strength	0.173 M
Average temperature	24.9°C
Analyte concentration range	55.4 µM to 52.3 µM
Methanol weight %	29.7 %
Dielectric constant	65.7
Water concentration	36.0 M
Number of pKas source	Manual (1)
Wavelength clipping	230.0 nm to 450.0 nm
pH clipping	1.453 to 12.514

Warnings and errors

Errors None
Warnings None

Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
Buffer in use	Yes			
Buffer type	Phosphate Buffer			

Assay Medium


Sample name: **M05**
 Assay name: **UV-metric psKa**
 Assay ID: **171-19008**
 Filename: **C:\Sirius_T3\Mehtap\20170918_exp04_uv_M01-M14\171-19008_M05_UV-metric psKa.t3r**

Experiment start time: **9/19/2017 6:09:04 AM**
 Analyst: **Dorothy Leverse**
 Instrument ID: **T311053**

Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			


Graphs

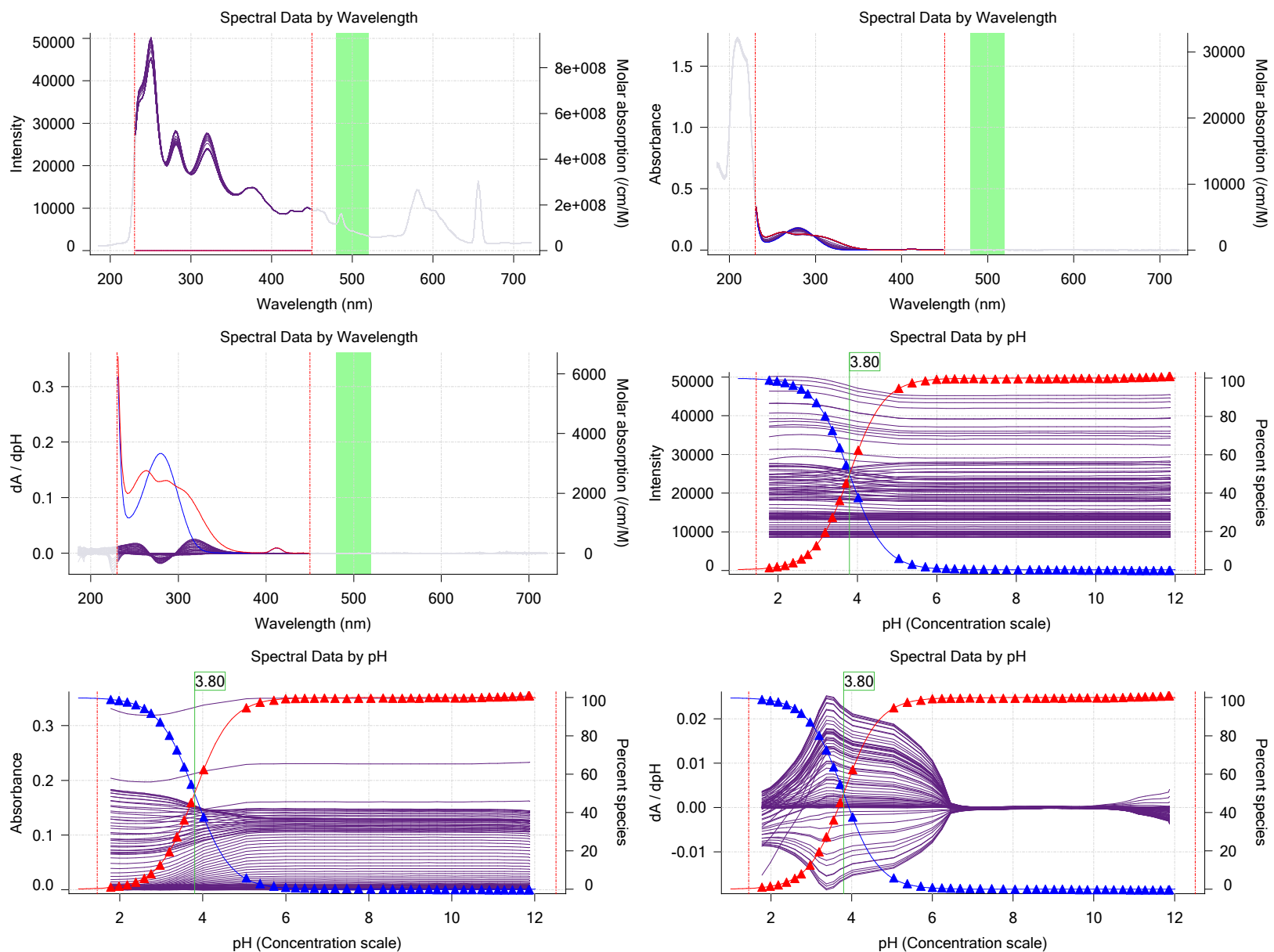

Sample name: **M05**
Assay name: **UV-metric psKa**
Assay ID: **171-19008**
Filename: **C:\Sirius_T3\Mehtap\20170918_exp04_uv_M01-M14\171-19008_M05_UV-metric psKa.t3r**

Experiment start time: **9/19/2017 6:09:04 AM**
Analyst: **Dorothy Levorse**
Instrument ID: **T311053**

Graphs (continued)

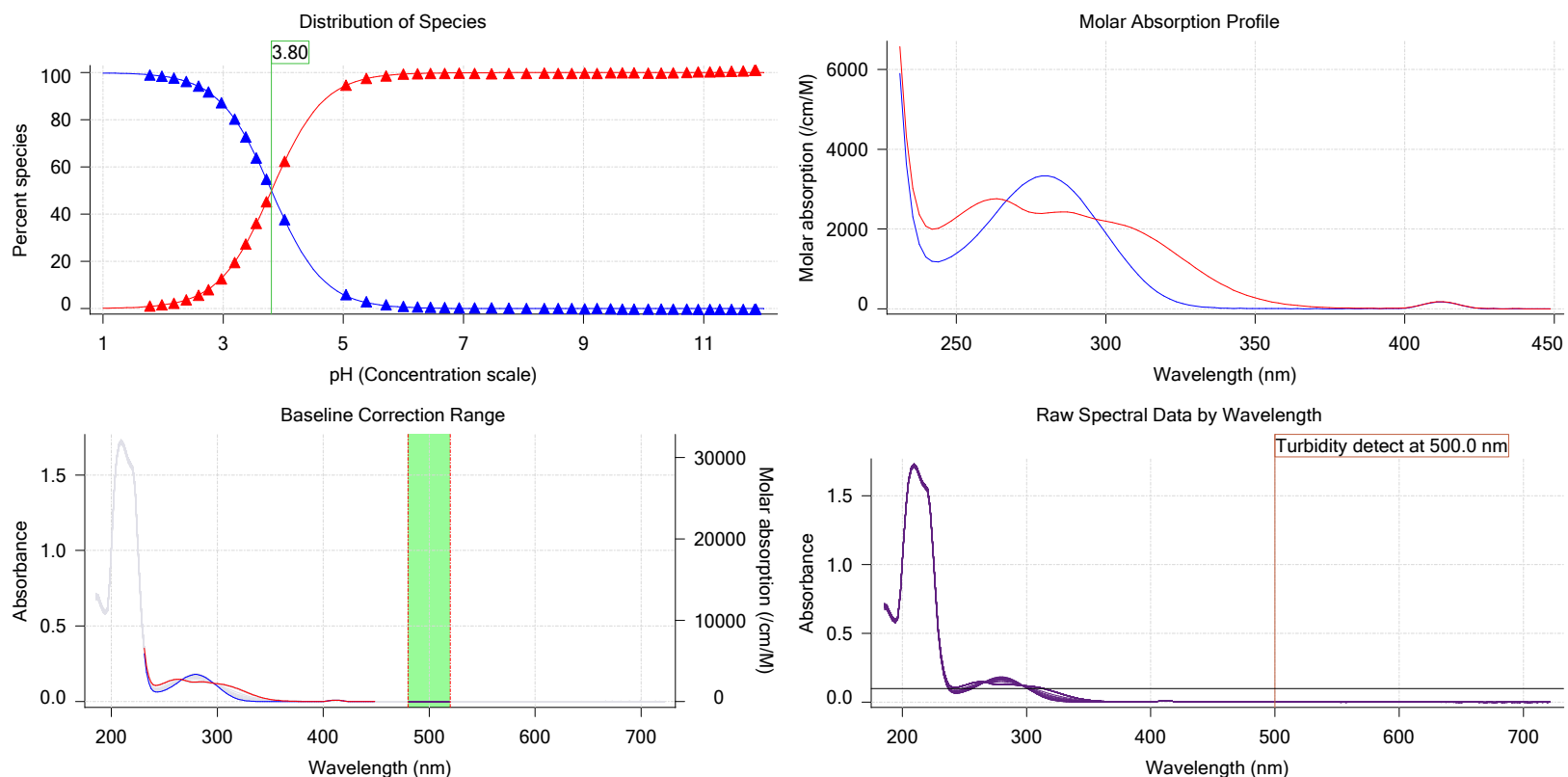

Assay Model

Settings	Value	Date/Time changed	Imported from
Sample name	M05	9/18/2017 4:10:29 PM	User entered value
Sample by	Volume		Default value
Sample volume	0.0030 mL	9/18/2017 4:10:29 PM	User entered value
Solvent	DMSO		Default value
Sample concentration	0.047100 M	9/18/2017 4:10:29 PM	User entered value
Solubility	Unknown		Default value
Molecular weight	304.77	9/18/2017 4:10:53 PM	User entered value
Individual pKa ionic environments	No		Default value
Number of pKas	2	9/18/2017 4:10:29 PM	User entered value
Sample is a	Ampholyte	9/18/2017 4:10:29 PM	User entered value
pKa 1	5.77	9/18/2017 4:10:29 PM	User entered value
Type	Base	9/18/2017 4:10:29 PM	User entered value
pKa 2	12.00	9/18/2017 4:10:29 PM	User entered value
Type	Acid	9/18/2017 4:10:29 PM	User entered value
logp (XH2 +)	-10.00		Default value
logP (neutral XH)	-10.00	9/18/2017 4:10:29 PM	User entered value
logP (X -)	-10.00		Default value

Assay Settings

Setting	Value	Original Value	Date/Time changed	Imported from
General Settings				
Analyst name	Dorothy Levorse			
Separate reference vial	Yes			
Standard Experiment Settings				
Number of titrations	3			
Minimum pH	2.000			
Maximum pH	12.000			

Sample name: **M05**
 Assay name: **UV-metric psKa**
 Assay ID: **171-19008**
 Filename: **C:\Sirius_T3\Mehtap\20170918_exp04_uv_M01-M14\171-19008_M05_UV-metric psKa.t3r**

Experiment start time: **9/19/2017 6:09:04 AM**
 Analyst: **Dorothy Levorse**
 Instrument ID: **T311053**

Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
pH step between points of	0.200			
Minimum titrant addition	0.00002 mL			
Maximum titrant addition	0.10000 mL			
Argon flow rate	100%			
Start titration using	Cautious pH adjust			
Advanced General Settings				
Detect turbidity using	Spectrometer			
Monitor at a wavelength of	500.0 nm			
Absorbance threshold of	0.100			
Collect turbidity sensor data	No			
Stir after titrant addition for	5 seconds			
For titrant addition, stir at	15%			
Titrant Pre-Dose				
Titrant pre-dose	None			
Assay Medium				
Cosolvent in use	Yes			
Cosolvent type	Methanol			
Cosolvent volume	1.15 mL			
Cosolvent added	Automatic			
ISA water volume	0.35 mL			
Water added	Automatic			
After water addition, stir for	5 seconds			
At a speed of	15%			
Buffer in use	Yes			
Buffer type	Phosphate Buffer			
Volume of buffer introduced	0.025000 mL			
Add buffer manually	Manual			
After medium addition, stir for	5 seconds			
Sample Sonication				
Sonicate	No			
Sample Dissolution				
Perform a dissolution stage	No			
Carbonate purge				
Perform a carbonate purge	No			
Temperature Control				
Wait for temperature	Yes			
Required start temperature	25.0°C			
Acceptable deviation	0.5°C			
Time to wait	60 seconds			
Stir speed of	15%			
Titration 1				
Titrate from	Low to high pH			
Adjust to start pH	Yes			
After pH adjust stir for	10 seconds			
Titration 2				
Titrate from	Low to high pH			
Additional cosolvent volume	0.00 mL			
Add additional water	0.15 mL			
Additional water added	Automatic			
After pH adjust stir for	10 seconds			
Titration 3				
Titrate from	Low to high pH			
Additional cosolvent volume	0.00 mL			
Add additional water	0.34 mL			
Additional water added	Automatic			
After pH adjust stir for	10 seconds			
Data Point Stability				
Stir during data point collection	Yes			

Sample name: **M05**
Assay name: **UV-metric psKa**
Assay ID: **171-19008**
Filename: **C:\Sirius_T3\Mehtap\20170918_exp04_uv_M01-M14\171-19008_M05_UV-metric psKa.t3r**

Experiment start time: **9/19/2017 6:09:04 AM**
Analyst: **Dorothy Levorse**
Instrument ID: **T311053**

Assay Settings (continued)

Setting	Value	Original Value	Date/Time changed	Imported from
For point collection, stir at	15%			
Delay before data point collection	0 seconds			
Number of points to average	20 points			
Time interval between points	0.50 seconds			
Required maximum standard deviation	0.00500 dpH/dt			
Stability timeout after	60 seconds			
Experiment cleanup				
Adjust pH to cleanup	To start pH			
And then stir for	60 seconds			
For cleaning, stir at	20%			
Then add water volume	0.25 mL			
And then stir for	30 seconds			

Calibration Settings

Setting	Value	Date/Time changed	Imported from
Four-Plus alpha	0.094	9/19/2017 6:09:04 AM	C:\Sirius_T3\171-18009_Blank standardisation.t3r
Four-Plus S	1.0023	9/19/2017 6:09:04 AM	C:\Sirius_T3\171-18009_Blank standardisation.t3r
Four-Plus jH	0.8	9/19/2017 6:09:04 AM	C:\Sirius_T3\171-18009_Blank standardisation.t3r
Four-Plus jOH	-0.5	9/19/2017 6:09:04 AM	C:\Sirius_T3\171-18009_Blank standardisation.t3r
Base concentration factor	1.015	9/19/2017 6:09:04 AM	C:\Sirius_T3\KOH17111.t3r
Acid concentration factor	1.006	9/19/2017 6:09:04 AM	C:\Sirius_T3\171-18009_Blank standardisation.t3r

Instrument Settings

Setting	Value	Batch Id	Install date
Instrument owner	Merck		
Instrument ID	T311053		
Instrument type	T3 Simulator		
Software version	1.1.3.0		
Dispenser module		T3DM1100253	3/31/2009 6:24:52 AM
Dispenser 0	Water		3/31/2009 6:25:05 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Water (0.15 M KCl)	8-18-17	9/18/2017 9:13:04 AM
Dispenser 2	Acid		3/31/2009 6:25:11 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Acid (0.5 M HCl)	166940	9/8/2017 9:21:27 AM
Dispenser 1	Base		3/31/2009 6:25:21 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Base (0.5 M KOH)	01/06/17	9/8/2017 9:20:03 AM
Dispenser 5	Cosolvent		3/31/2009 6:26:24 AM
Syringe volume	2.5 mL		
Firmware version	1.2.1(r2)		
Distribution valve 5	Distribution Valve		3/31/2009 6:28:19 AM
Firmware version	1.1.3		
Port A	Methanol (80%, 0.15 M KCl)	8-15-17	9/13/2017 12:23:11 PM
Dispenser 3	Buffer		8/3/2010 6:05:16 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Phosphate Buffer		9/12/2017 12:32:29 PM
Dispenser 6	Octanol		10/22/2010 11:52:43 AM
Syringe volume	0.5 mL		
Firmware version	1.2.1(r2)		
Titrant	Octanol	9-14-17	9/14/2017 10:30:38 AM
Titration		T3TM1100153	3/31/2009 6:24:17 AM

Sample name: **M05**
 Assay name: **UV-metric psKa**
 Assay ID: **171-19008**
 Filename: **C:\Sirius_T3\Mehtap\20170918_exp04_uv_M01-M14\171-19008_M05_UV-metric psKa.t3r**

Experiment start time: **9/19/2017 6:09:04 AM**
 Analyst: **Dorothy Levorse**
 Instrument ID: **T311053**

Instrument Settings (continued)

Setting	Value	Batch Id	Install date
Horizontal axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Probe I/O firmware version	1.1.1		
Electrode	T3 Electrode	T3E0769	8/15/2017 10:21:54 AM
E0 calibration	-7.47 mV		9/19/2017 6:09:28 AM
Filling solution	3M KCl	KCL095	9/18/2017 9:17:15 AM
Liquids			
Wash 1	50% IPA:50% Water		9/18/2017 9:09:36 AM
Wash 2	0.5% Triton X-100 in H2O		9/18/2017 9:09:39 AM
Buffer position 1	pH7 Wash		9/18/2017 9:09:41 AM
Buffer position 2	pH 7		9/18/2017 9:09:44 AM
Storage position			9/18/2017 9:10:43 AM
Wash water	7.4e+003 mL	9-18-17	9/18/2017 8:54:32 AM
Waste	2.6e+003 mL		9/18/2017 8:54:39 AM
Temperature controller			8/5/2010 7:35:13 AM
Turbidity detector			3/31/2009 6:24:45 AM
Spectrometer		072390	11/23/2010 12:22:28 PM
Dip probe		11086	
Wavelength coefficient A0	185.563		
Wavelength coefficient A1	2.17439		
Wavelength coefficient A2	-0.000285622		
Total lamp lit time	143:53:30		11/23/2010 12:22:28 PM
Calibrated on	9/18/2017 9:35:14 AM		
Integration time	11		
Scans averaged	10		
Autoloader		T3AL1100237	11/10/2015 10:34:13 AM
Left-right axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Front-back axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Vertical axis firmware version	1.17 AI1DI2DO2 Stepper 2		
Chassis I/O firmware version	1.11 AI1DI0DO4 Norgren I/O		
Configuration			
Alternate titration position	Titration position		
Alternate reference position	Reference position		
Maximum standard vial volume	3.50 mL		
Maximum alternate vial volume	25.00 mL		
Automatic action idle period	5 minute(s)		
Titrant tube volume	1.3 mL		
Syringe flush count	3.50		
Flowing wash pump volume	20.0 mL		
Flowing wash stir duration	5 s		
Flowing wash stir speed	30%		
Solvent wash stir duration	5 s		
Solvent wash stir speed	30%		
Surfactant wash stir duration	5 s		
Surfactant wash stir speed	30%		
E0 calibration minimum number of points	10		
E0 calibration maximum standard deviation	0.01500		
E0 calibration timeout period	60 s		
E0 calibration stir duration	5 s		
E0 calibration preparation stir speed	30%		
E0 calibration buffer wash stir duration	5 s		
E0 calibration buffer wash stir speed	30%		
E0 calibration reading stir speed	0%		
Spectrometer calibration stir duration	5 s		
Spectrometer calibration stir speed	30%		
Spectrometer calibration wash pump volume	20.0 mL		
Spectrometer calibration wash stir duration	5 s		



Assay Settings

Sample name: **M05** Experiment start time: **9/19/2017 6:09:04 AM**
Assay name: **UV-metric psKa** Analyst: **Dorothy Levorse**
Assay ID: **17I-19008** Instrument ID: **T311053**
Filename: **C:\Sirius_T3\Mehtap\20170918_exp04_uv_M01-M14\17I-19008_M05_UV-metric psKa.t3r**

Instrument Settings (continued)

Setting	Value	Batch Id	Install date
Spectrometer calibration wash stir speed	30%		
Overhead dispense height	10000		

Refinement Settings

Setting	Value	Default value
Turbidity detection method	Spectrometer	Spectrometer
Turbidity wavelength to assess	500.0 nm	500.0 nm
Turbidity maximum absorbance	0.100	0.100
Turbidity probe threshold	50.00	50.00
Exclude turbid points	Yes	Yes
Low intensity warning threshold	100	100
Minimum absorbance change threshold	0.100	0.100
Eigenvector autocorrelation threshold	0.80	0.80
Maximum RMSD severe warning	0.250	0.250
Maximum RMSD warning	0.050	0.050

Tray Information

Title
Location E3